



US005690692A

United States Patent [19]

Fleming

[11] Patent Number: 5,690,692
[45] Date of Patent: Nov. 25, 1997

[54] BIO-ACTIVE FREQUENCY GENERATOR AND METHOD

[75] Inventor: Janet E. Fleming, San Diego, Calif.

[73] Assignee: Regeneration Technology, Tijuana, Mexico

[21] Appl. No.: 787,158

[22] Filed: Jan. 22, 1997

Related U.S. Application Data

[63] Continuation of Ser. No. 541,182, Oct. 11, 1995.

[51] Int. Cl.⁶ A61N 1/32

[52] U.S. Cl. 607/50; 607/66

[58] Field of Search 607/50, 66, 67, 607/142

[56] References Cited

U.S. PATENT DOCUMENTS

4,255,790 3/1981 Hondeghe .
4,738,250 4/1988 Fulkerson et al. .
4,895,154 1/1990 Bartelt et al. .
4,919,139 4/1990 Brodard .
5,048,523 9/1991 Yamasawa et al. .
5,285,781 2/1994 Brodard .
5,447,526 9/1995 Karsdon .

Primary Examiner—William E. Kamm


Assistant Examiner—Carl H. Layno

Attorney, Agent, or Firm—Gilliam, Duncan, Harms

[57] ABSTRACT

A system and method for generating bio-active frequencies. Such generated frequencies are very precise and are accurate within 0.001 Hz. The frequencies are used to advantage in health science applications such as killing microorganisms and viruses and enhancing tissue regeneration and may be used in various commercial and industrial applications, such as food processing. A specific precise frequency synthesizer is controlled by a programmable control that instructs the synthesizer to generate a specific precise frequency or a series of precise frequencies. A keyboard is actuable by a user to select a specific frequency or sequence of frequencies or a preprogrammed series of frequencies. Depending upon the specific frequency being generated, a circuit gates the generated signal ON or OFF with a predetermined determinable periods being selected for such ON/OFF periods. Immediately prior to initiation of a specific precise frequency signal a warning auditory signal is made to alert a user of the invention and the frequency signal is gradually increased at a predetermined rate from zero to a predetermined level controllable by a user. A display indicates information representative of the specific frequencies selected and the period of running of each. Application electrodes for the device are disclosed as well as an indicator as to whether the electrodes are in proper contact with a subject.

7 Claims, 23 Drawing Sheets

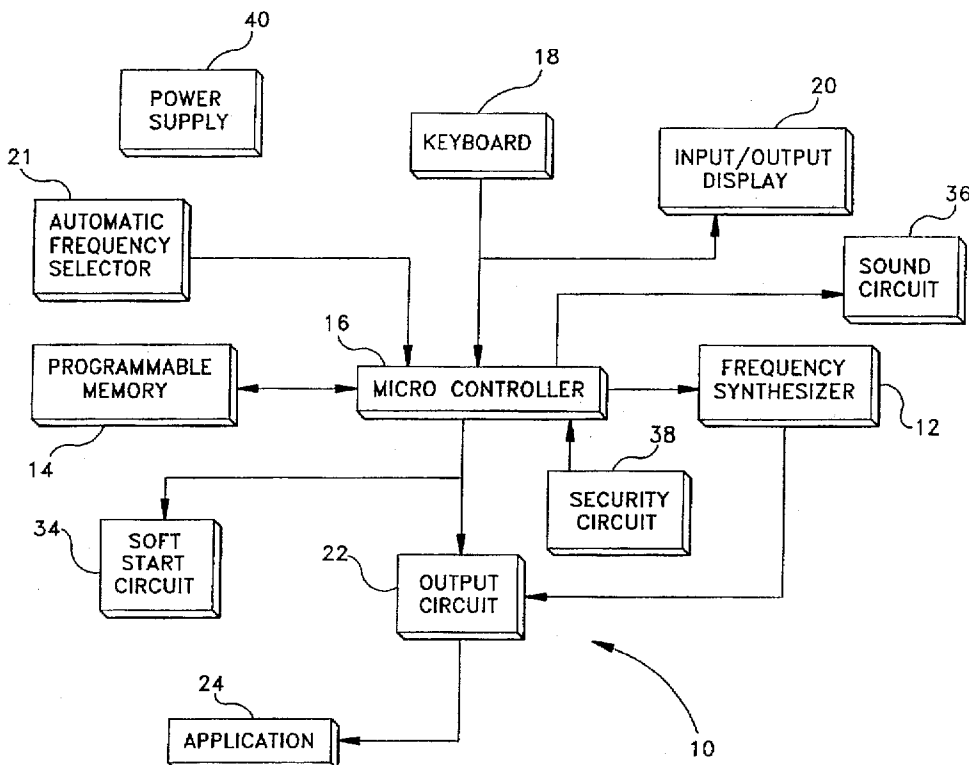

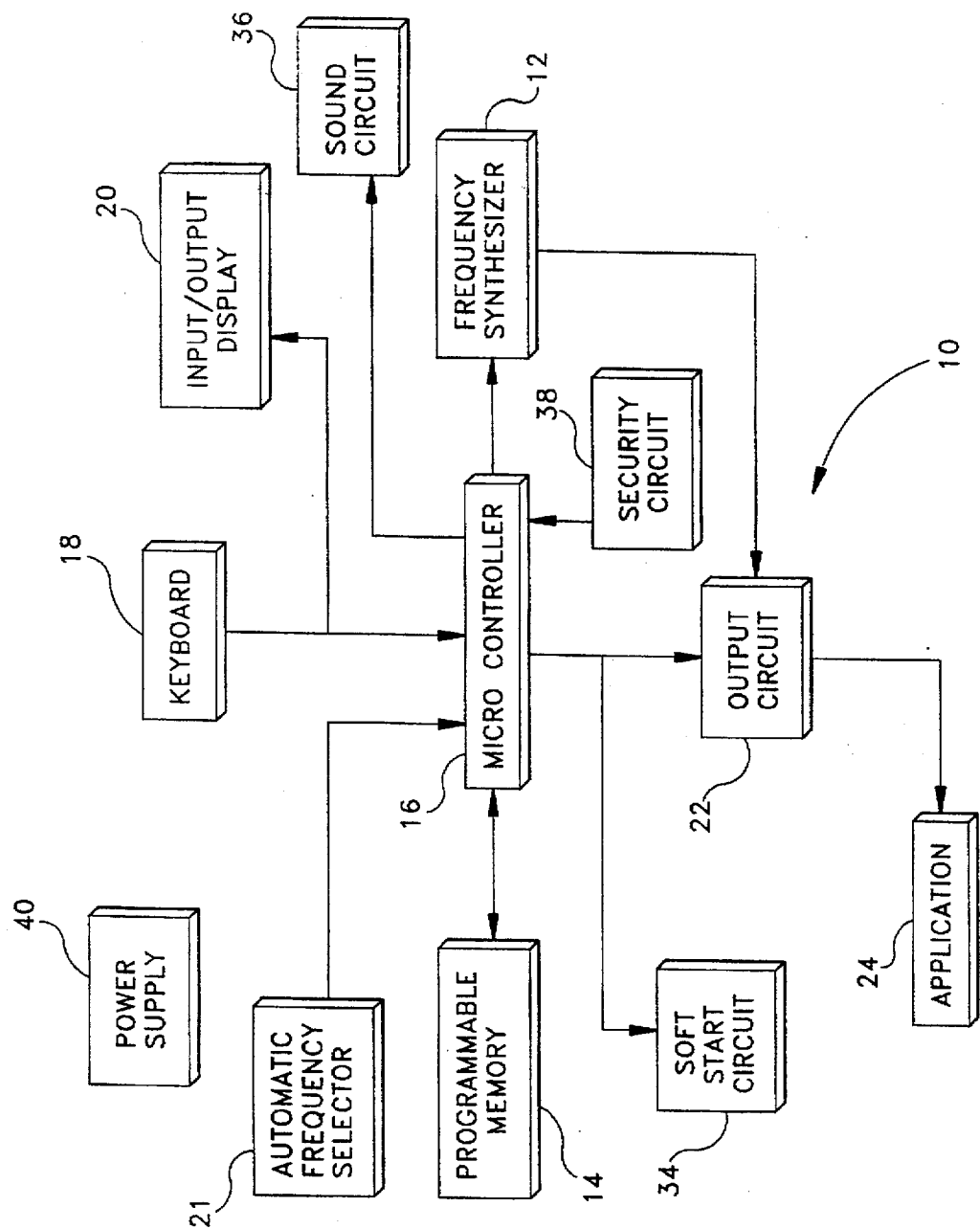


FIG. 1

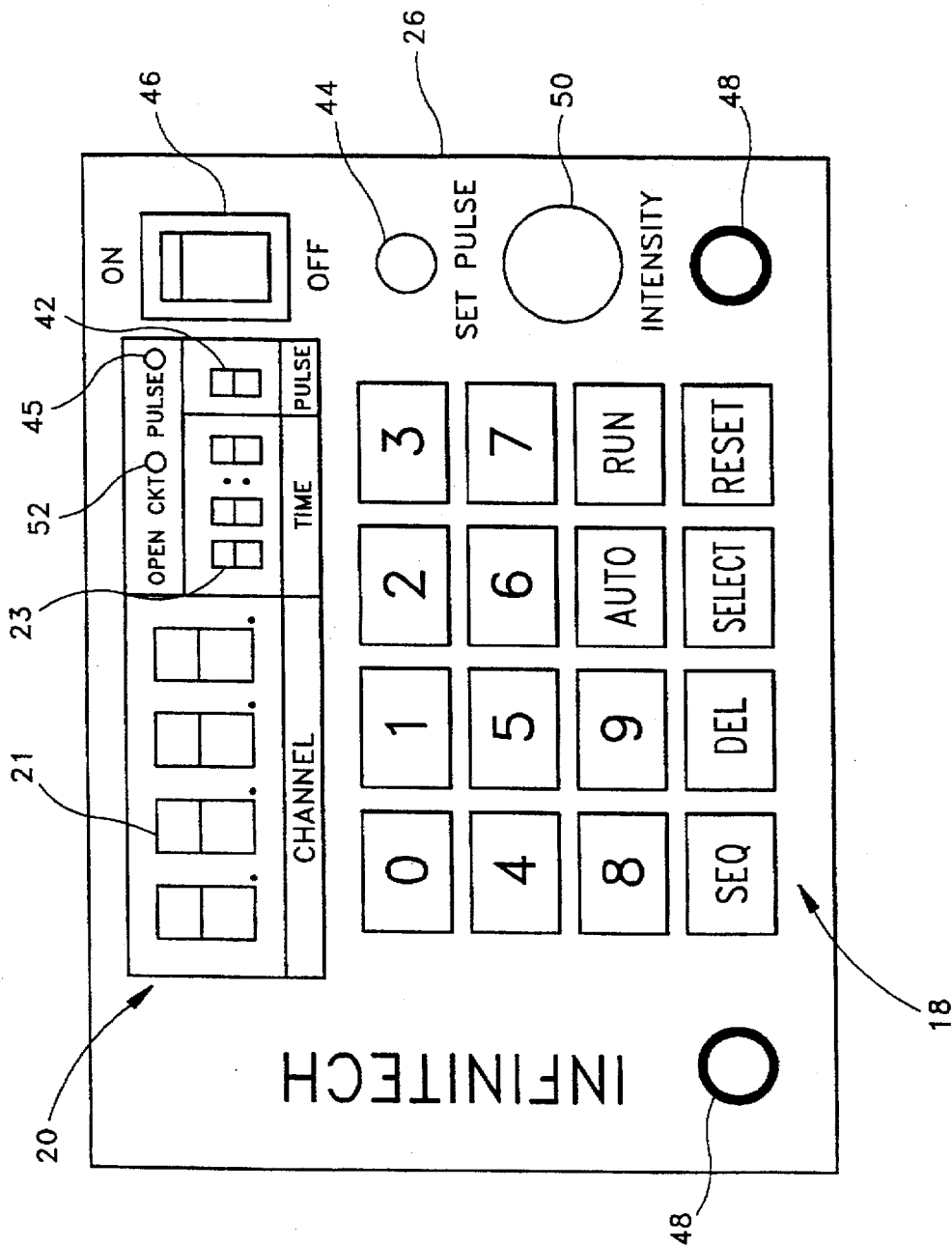


FIG. 2

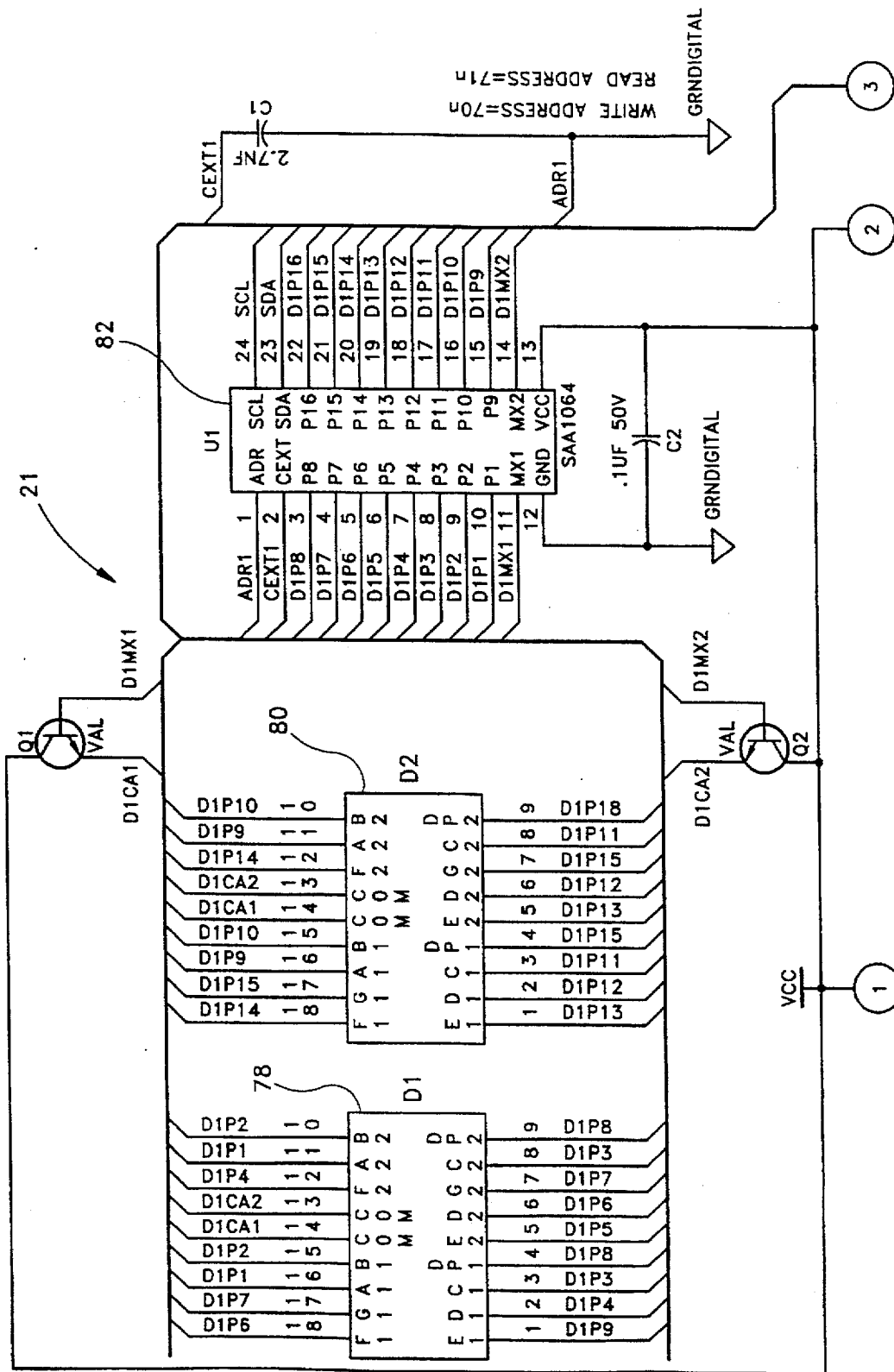


FIG. 3a

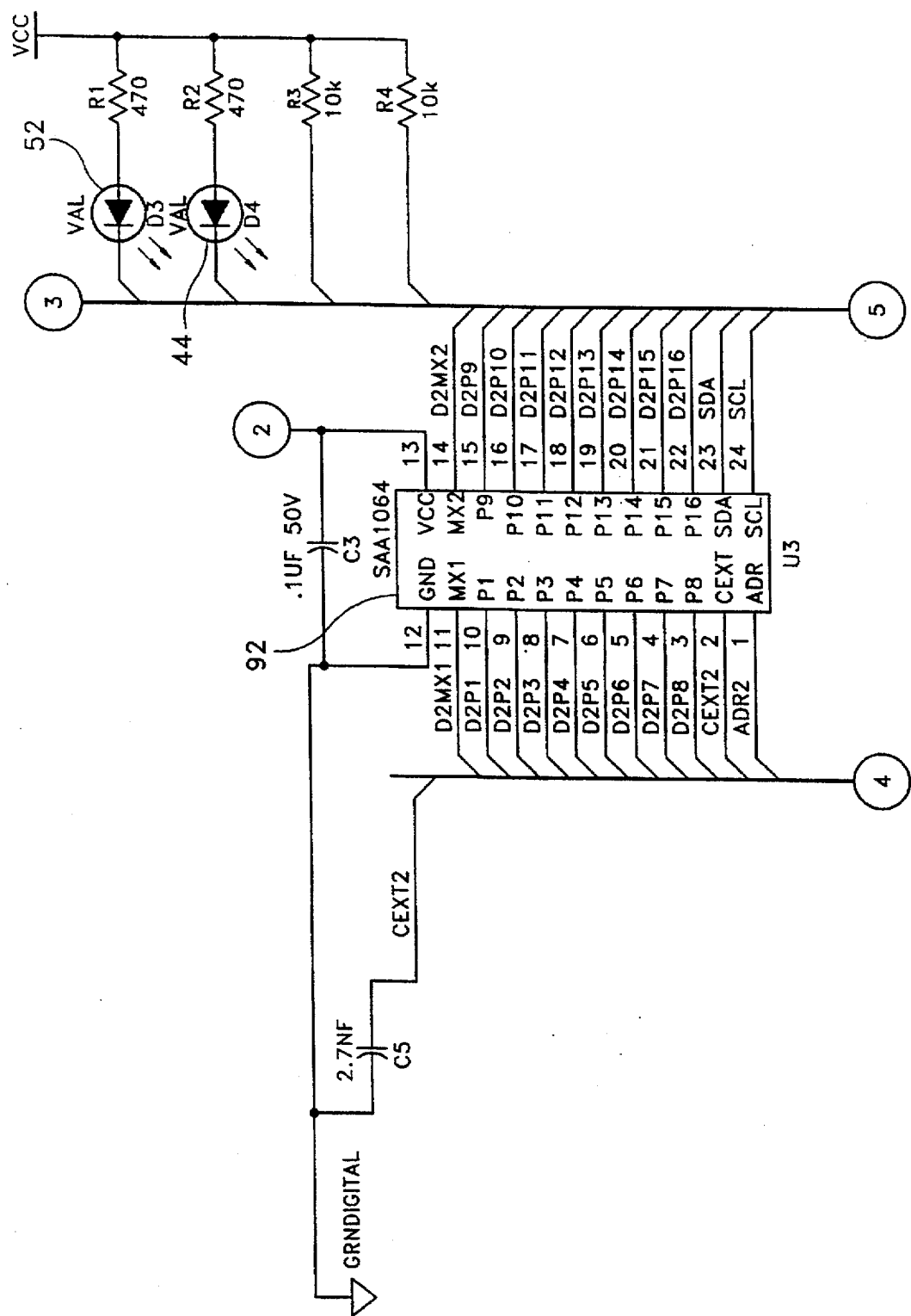


FIG. 3b

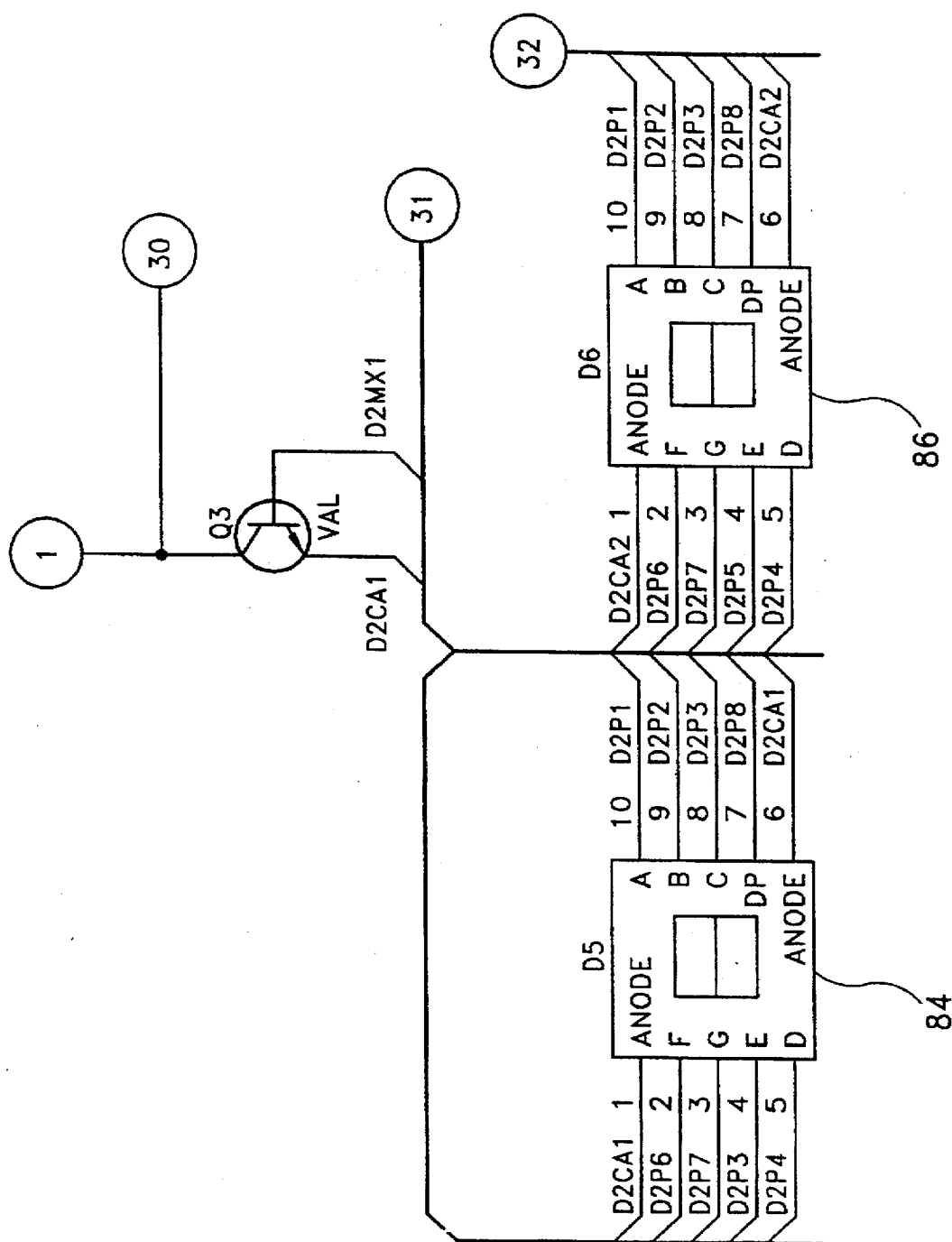


FIG. 3c

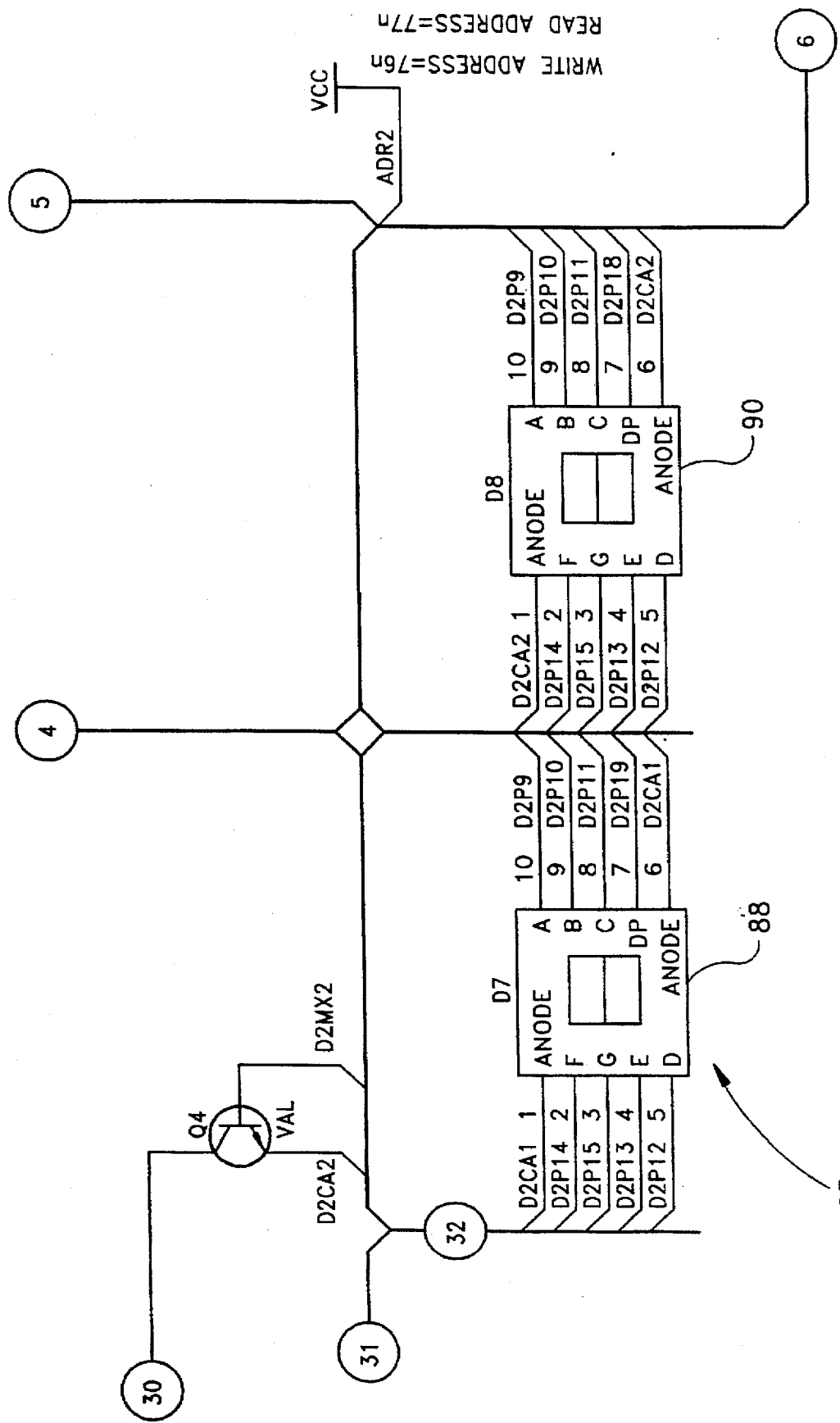
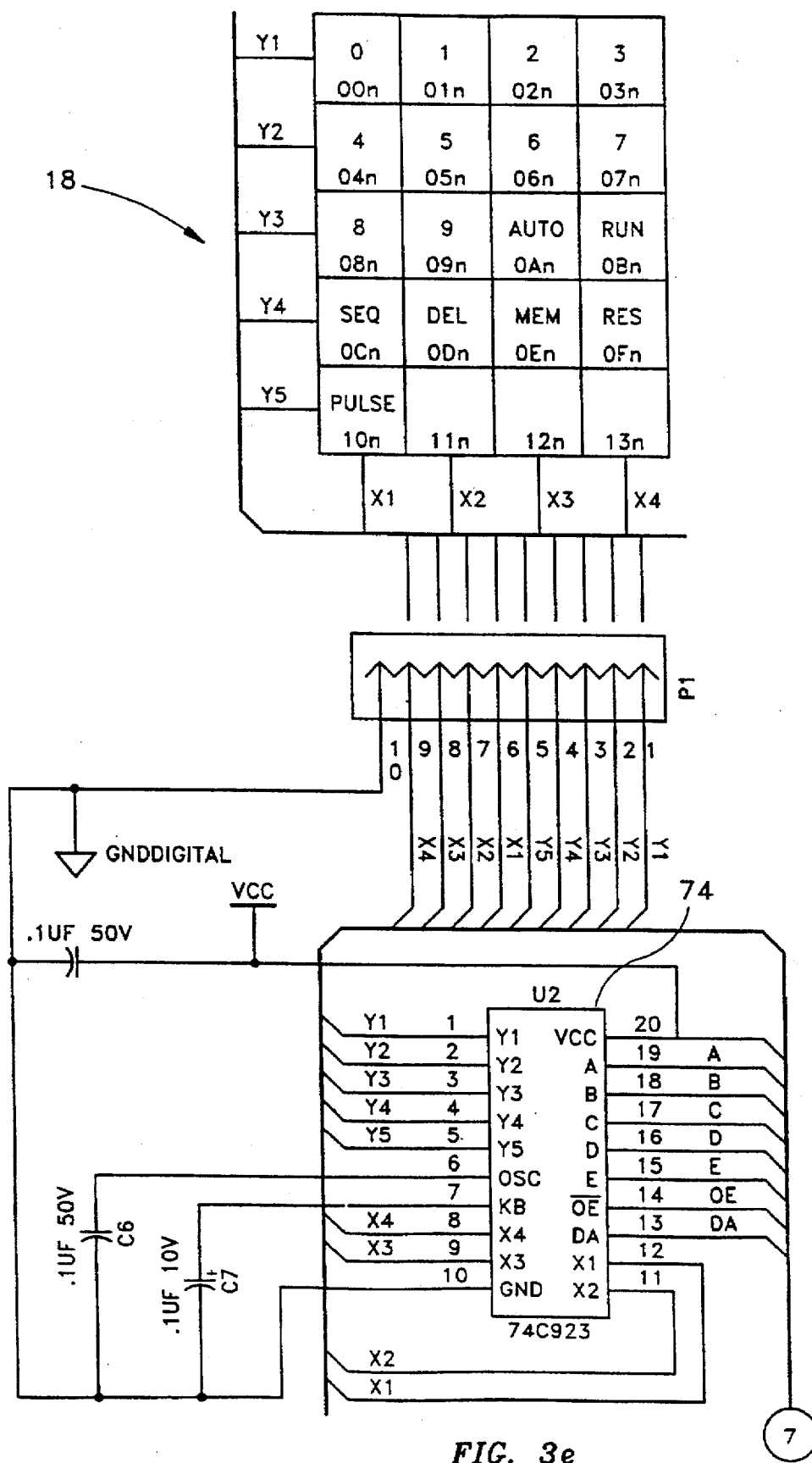


FIG. 3d



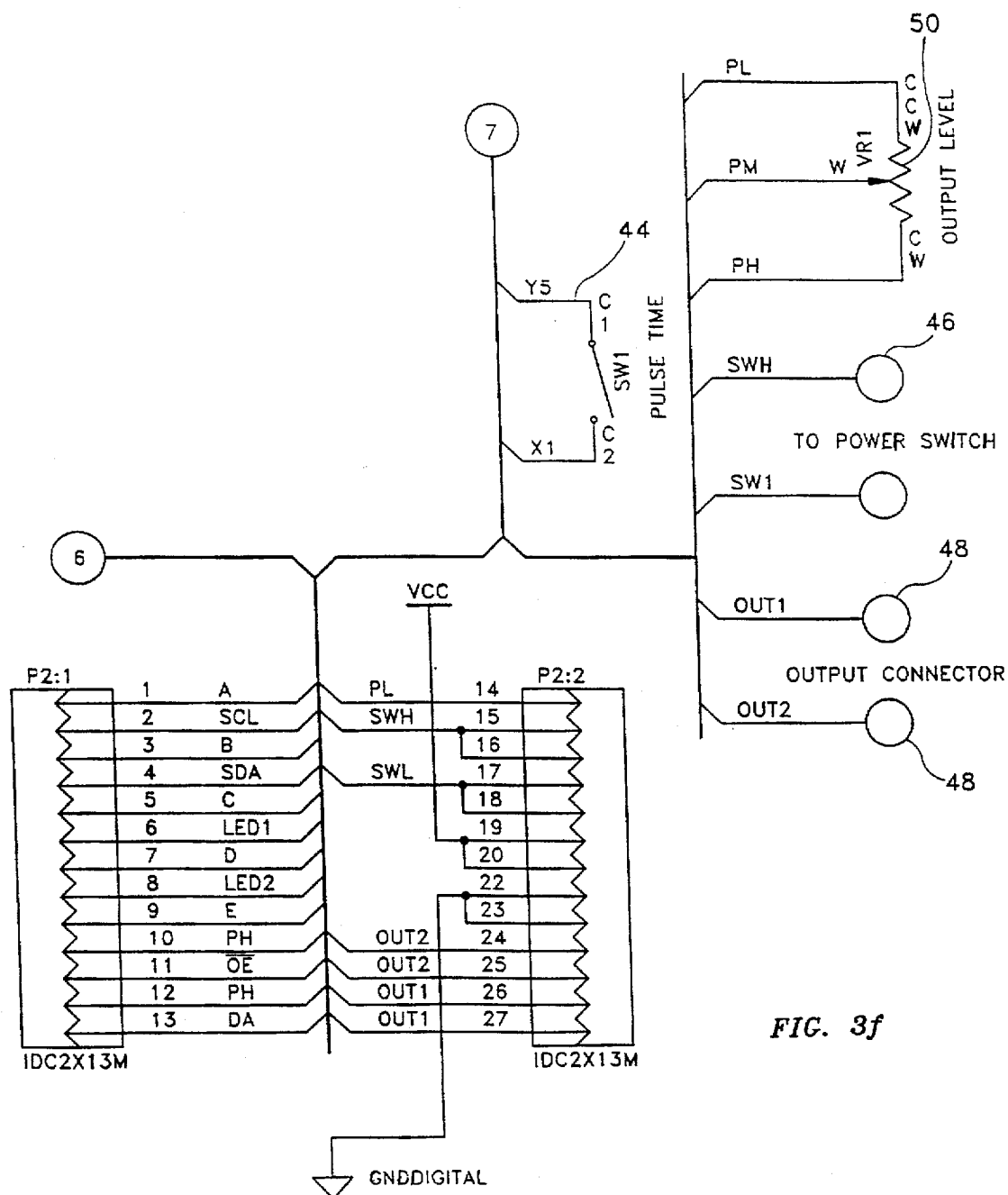


FIG. 3f

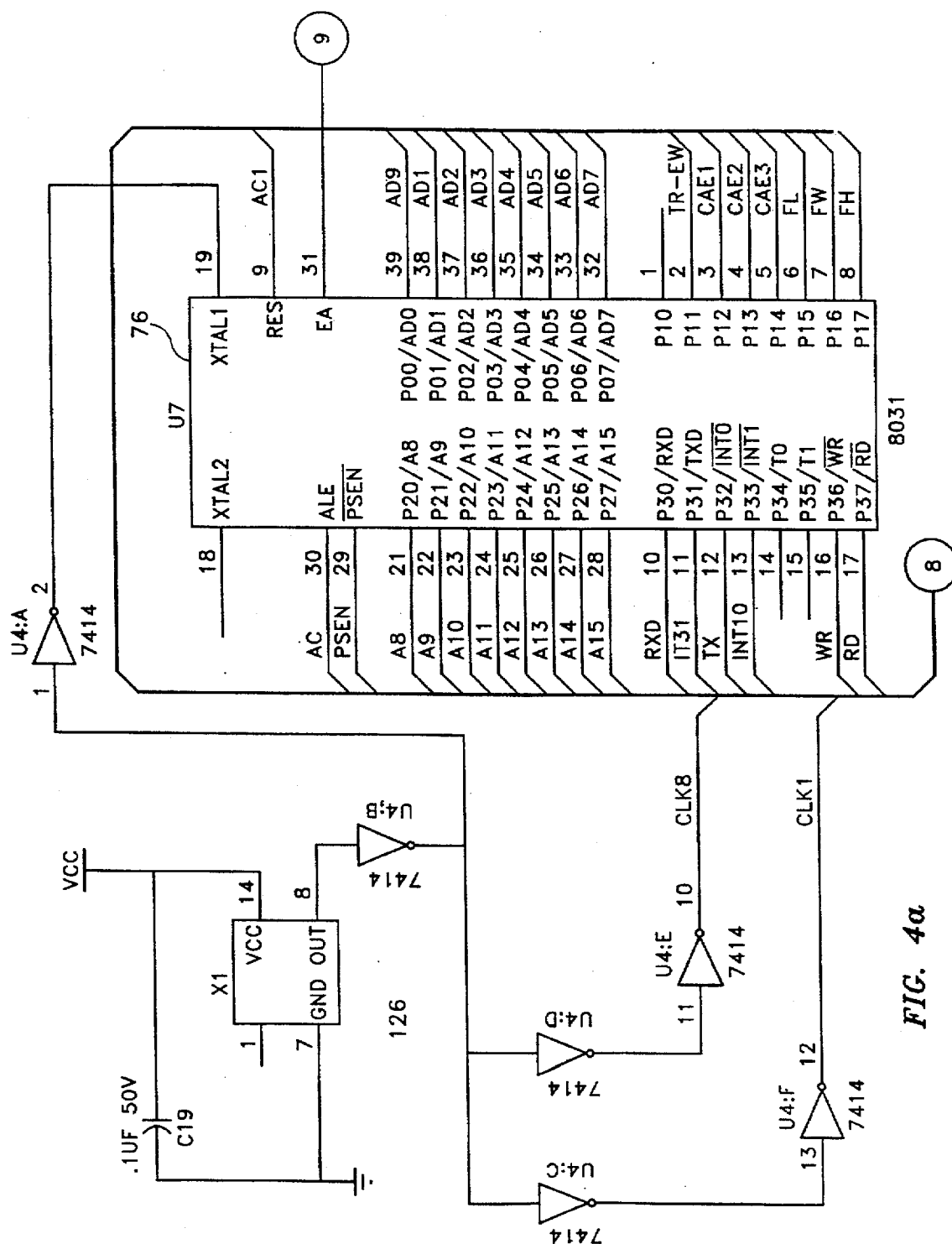


FIG. 4a

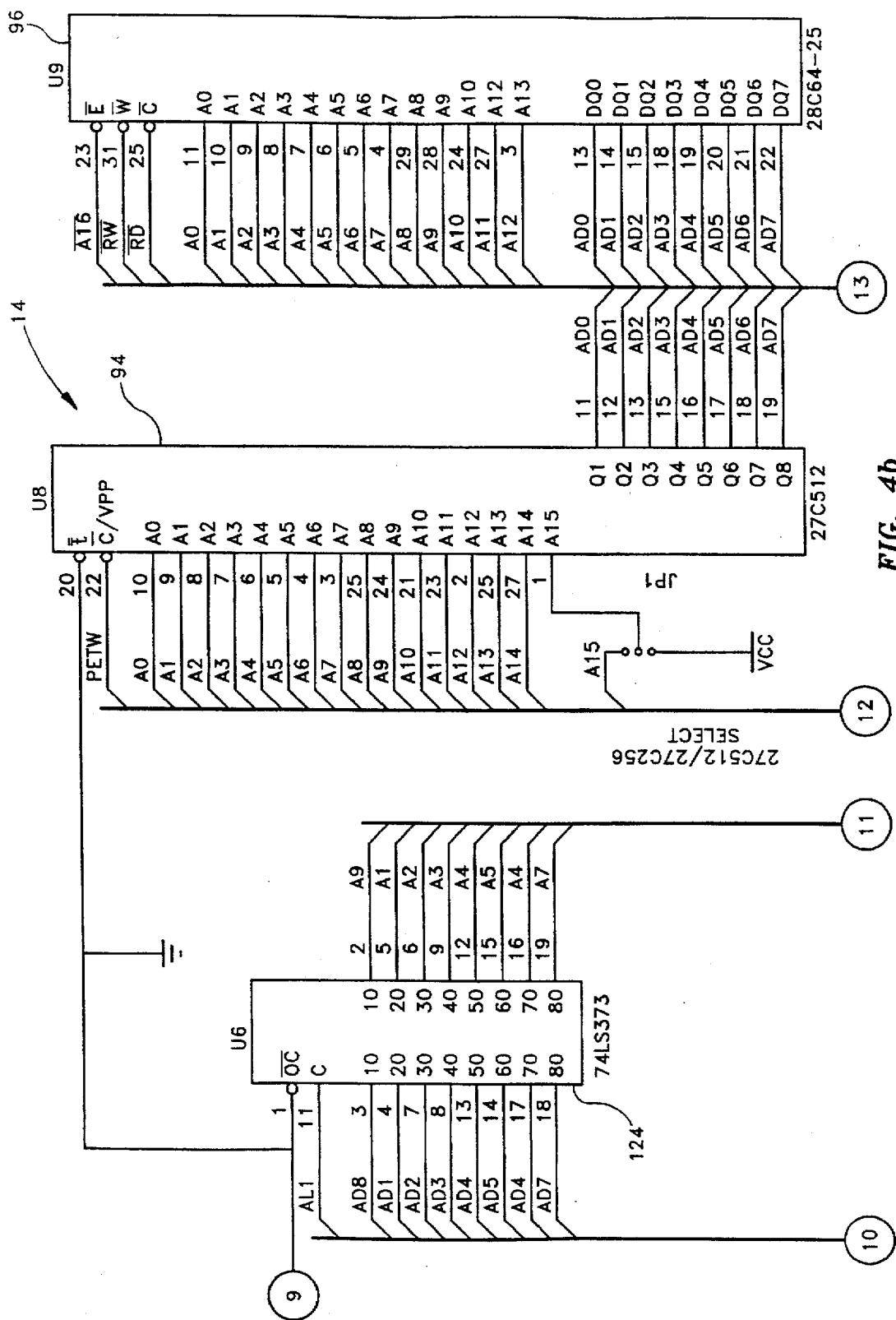


FIG. 4b

ROM (U8) ADDRESSES 0000h TO FFFFh
RAM (U9) ADDRESSES 0000n TO 1FFFh
U12 ADDRESS 6000n
U16 ADDRESSES 2000n TO 2001h
U17 ADDRESSES 4000n TO 4003n

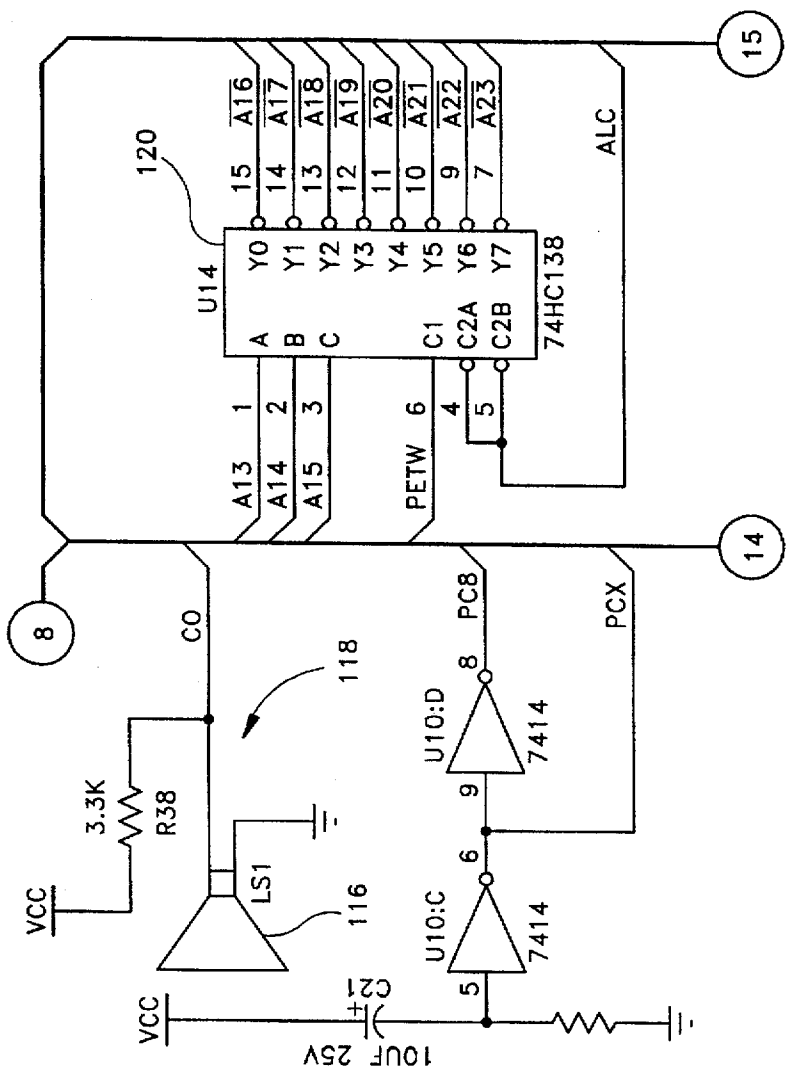


FIG. 4c

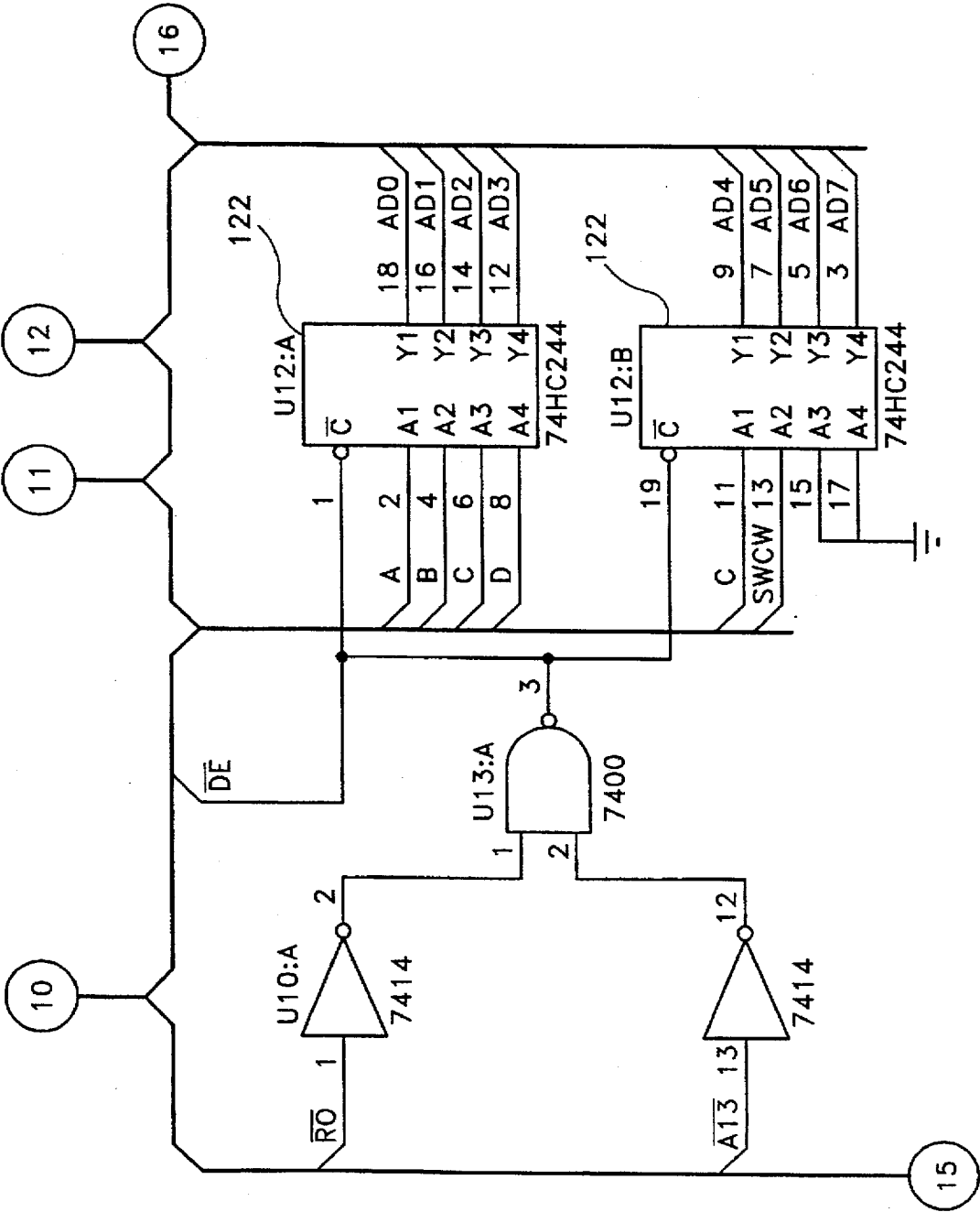


FIG. 4d

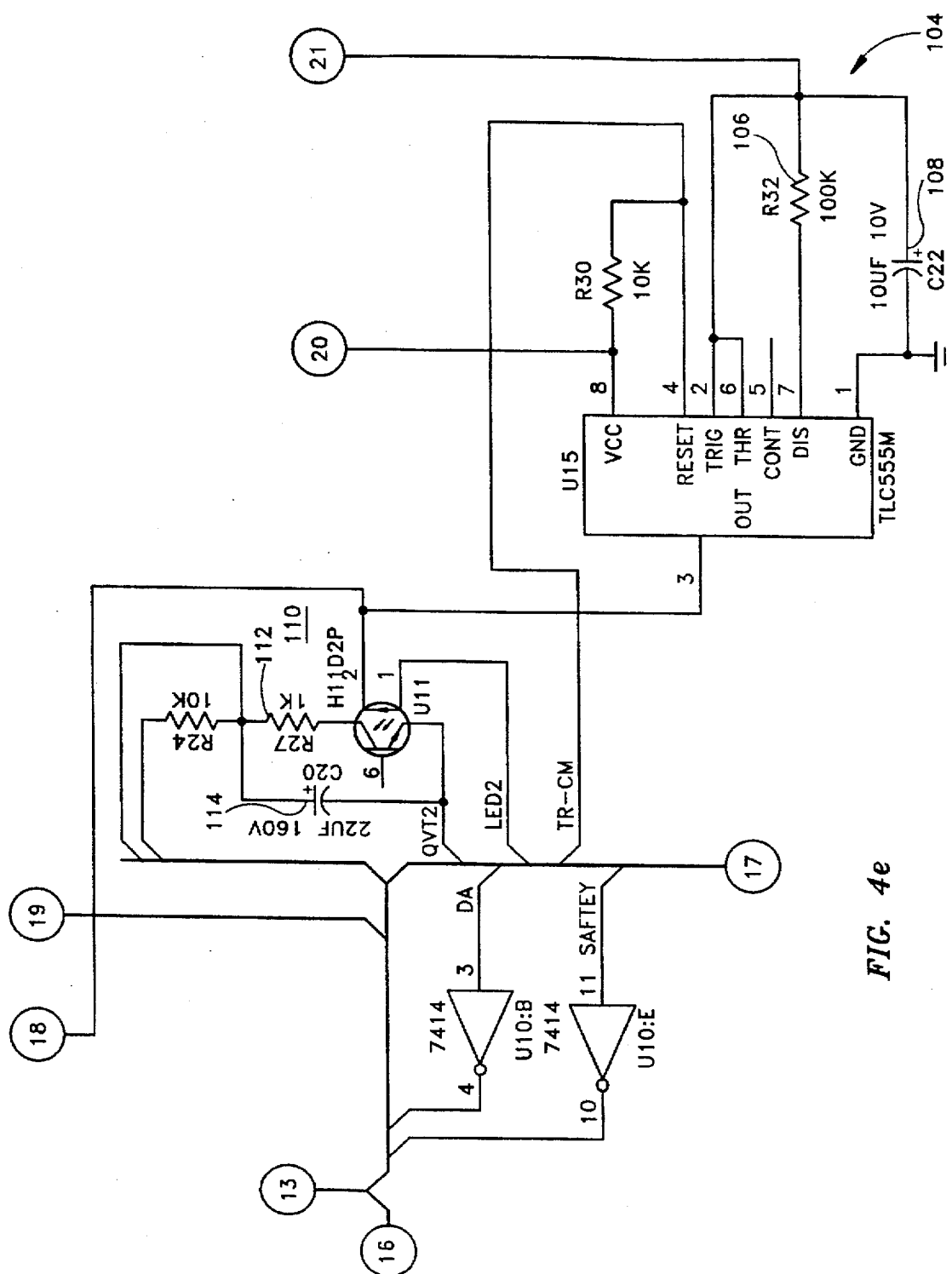


FIG. 4e

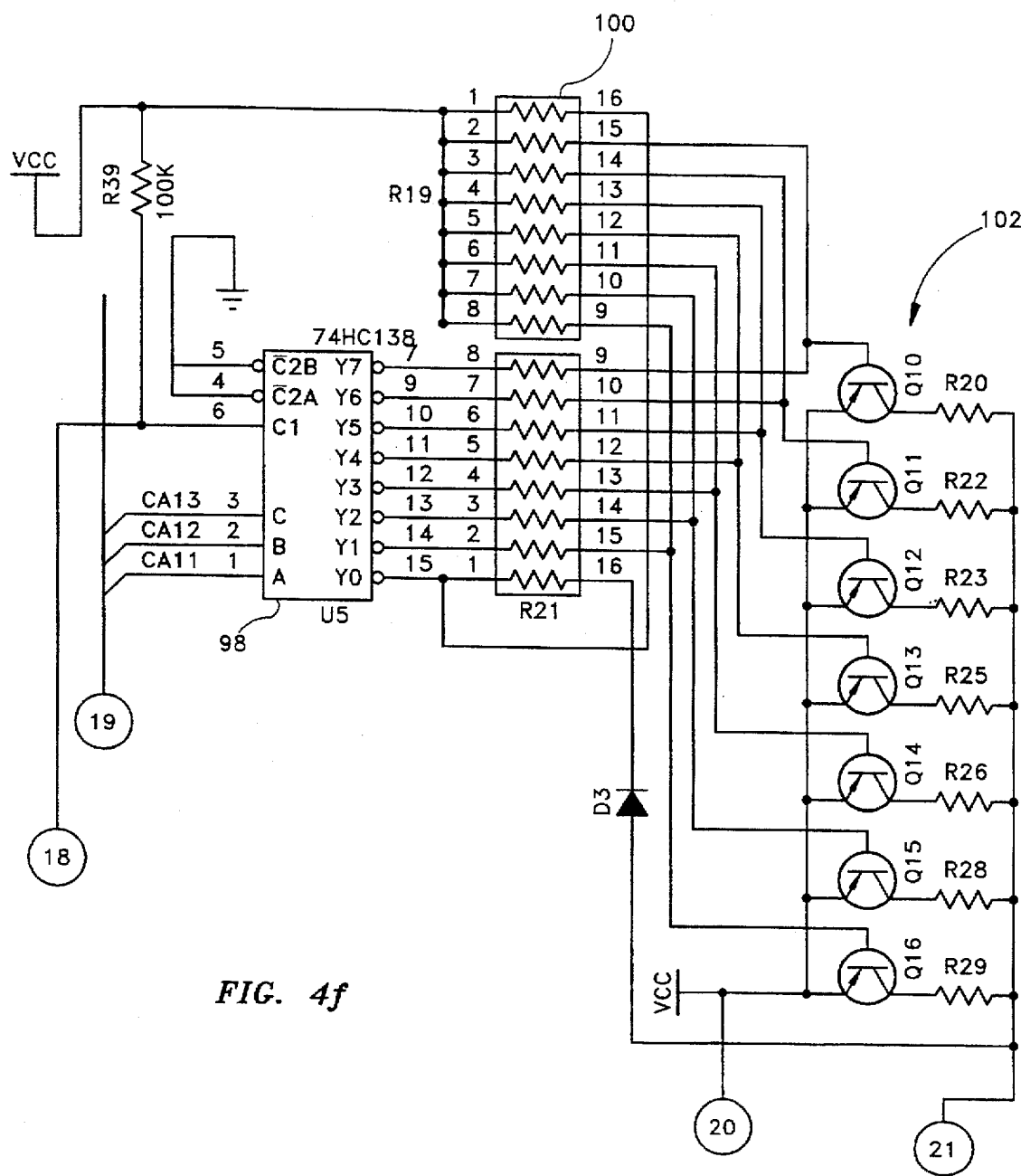


FIG. 4f

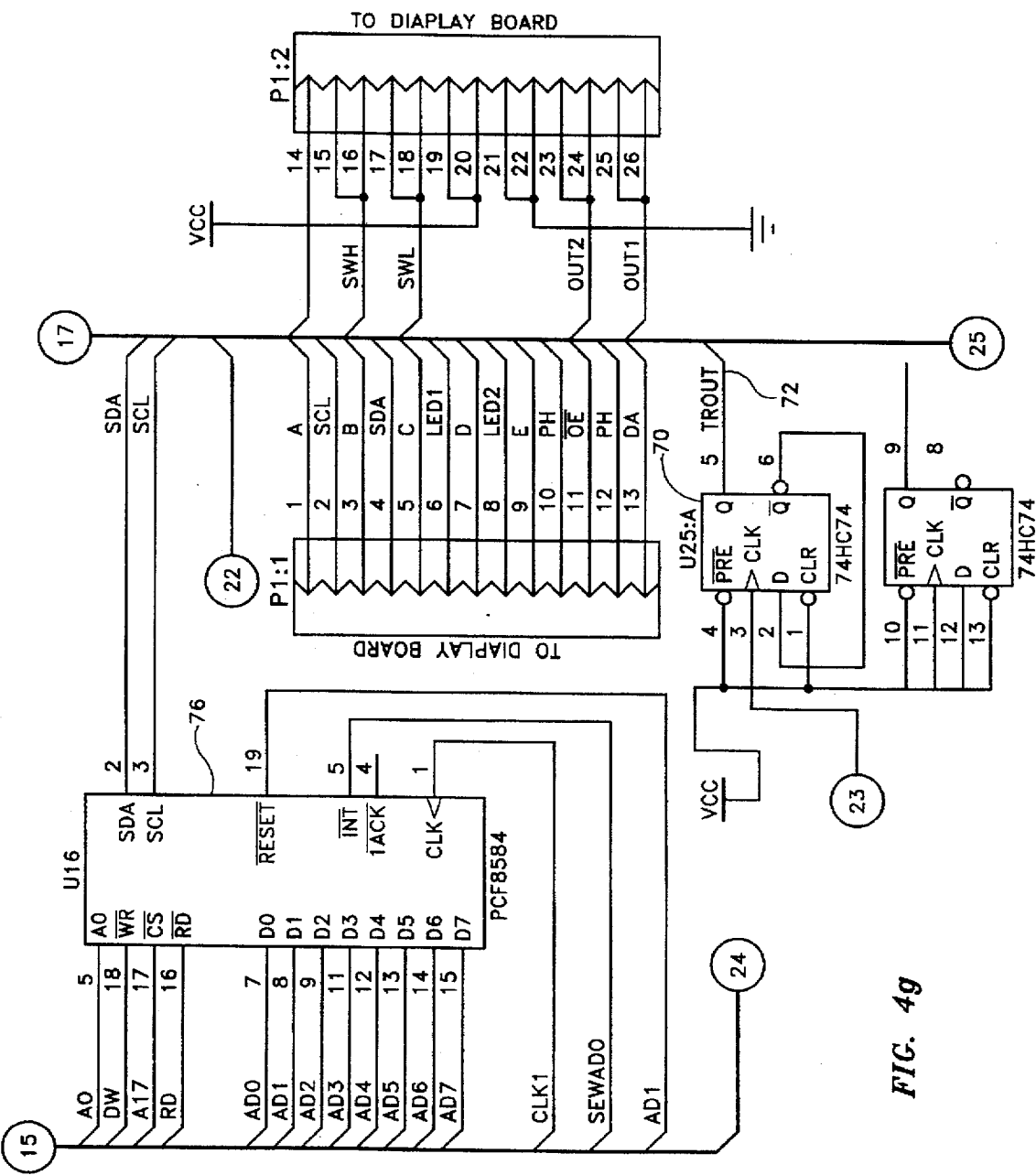


FIG. 4g

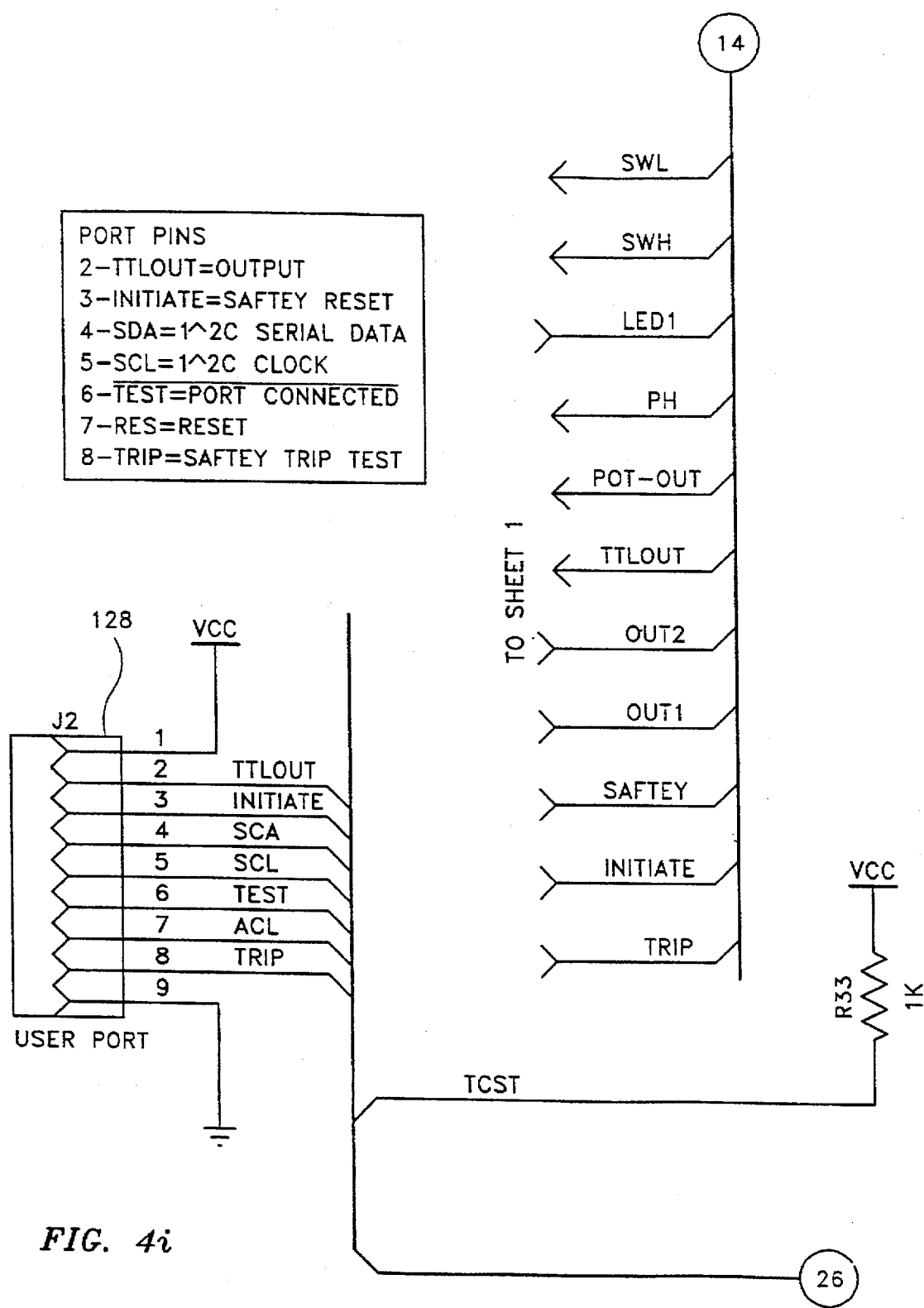
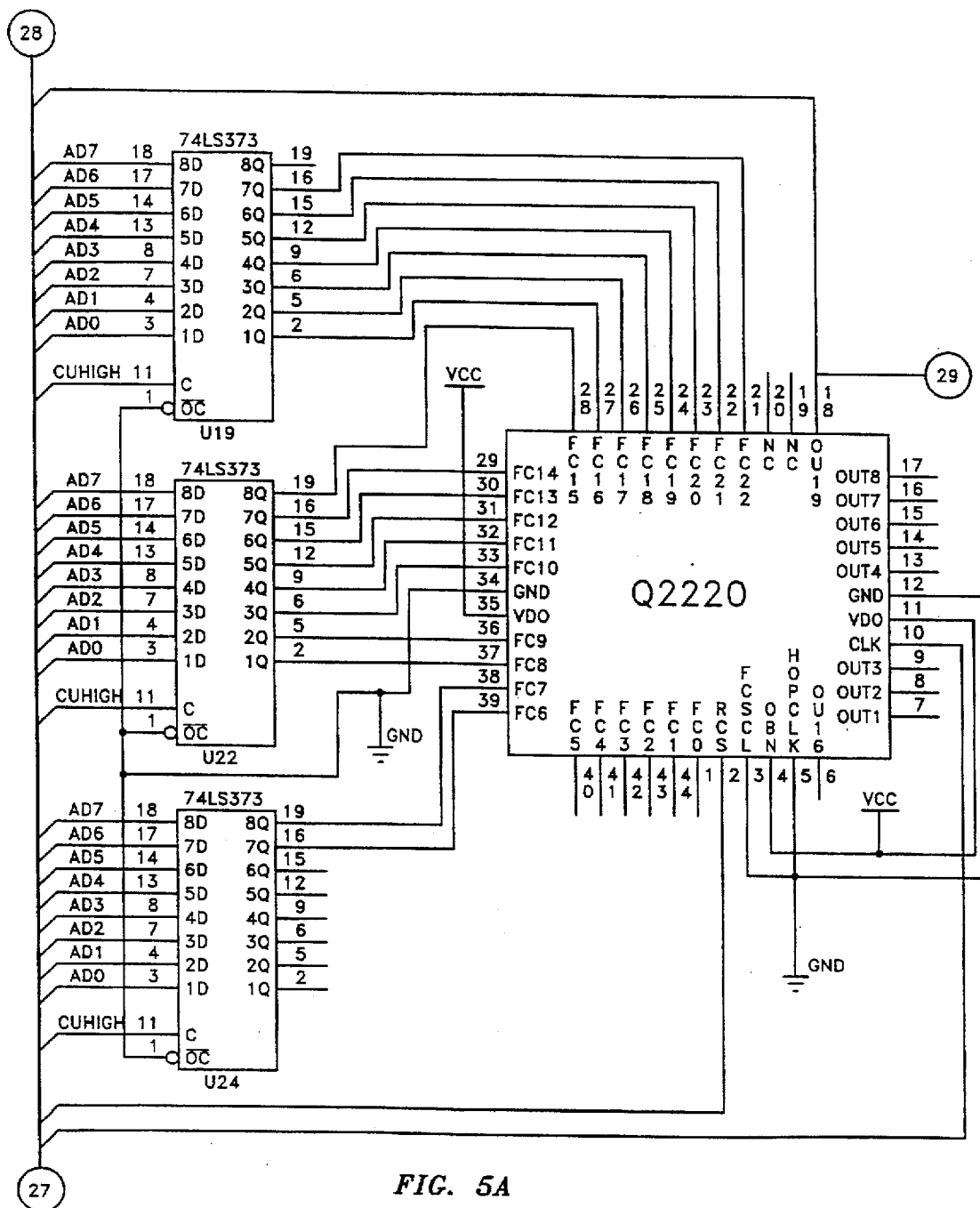


FIG. 4i



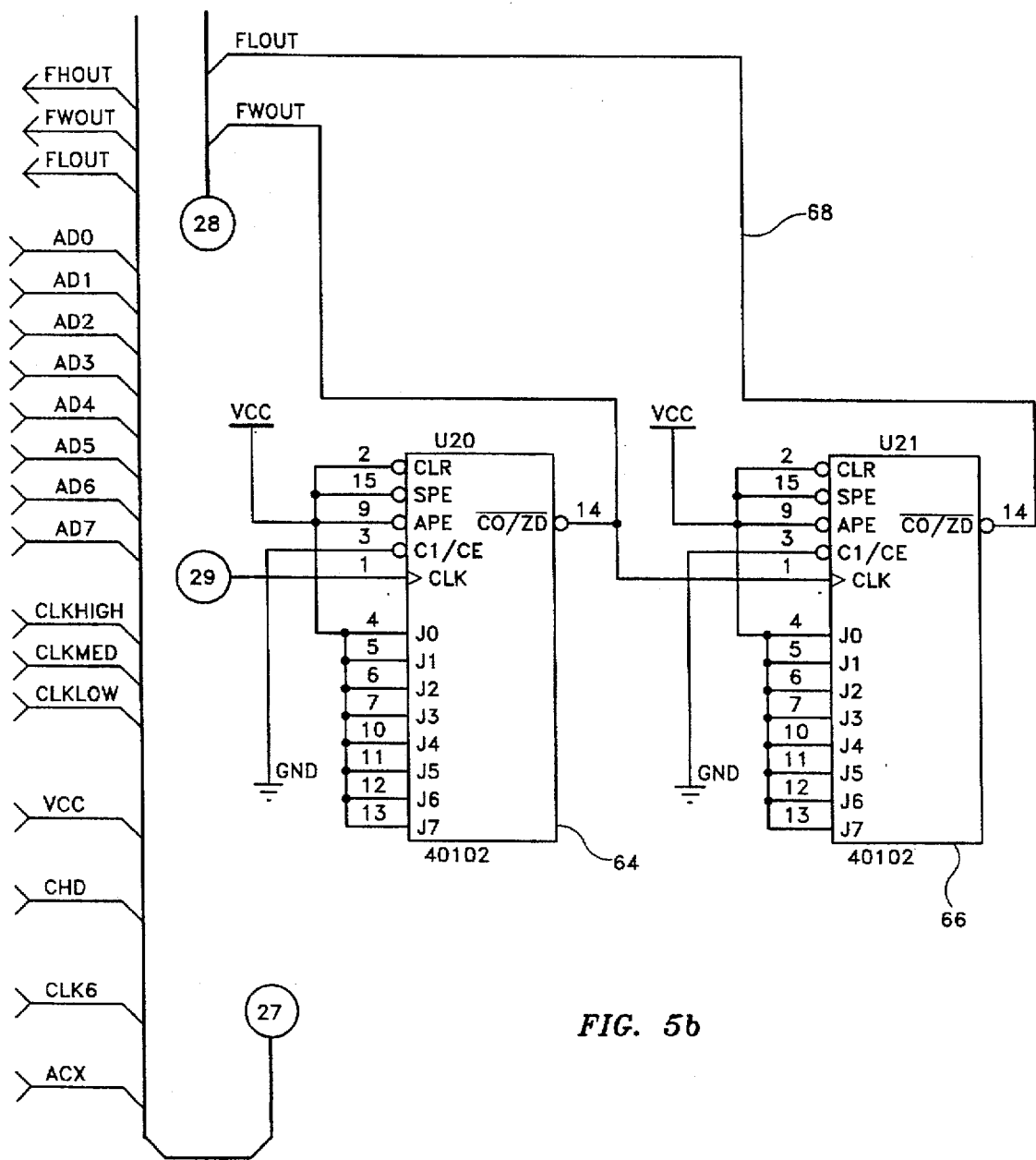


FIG. 5b

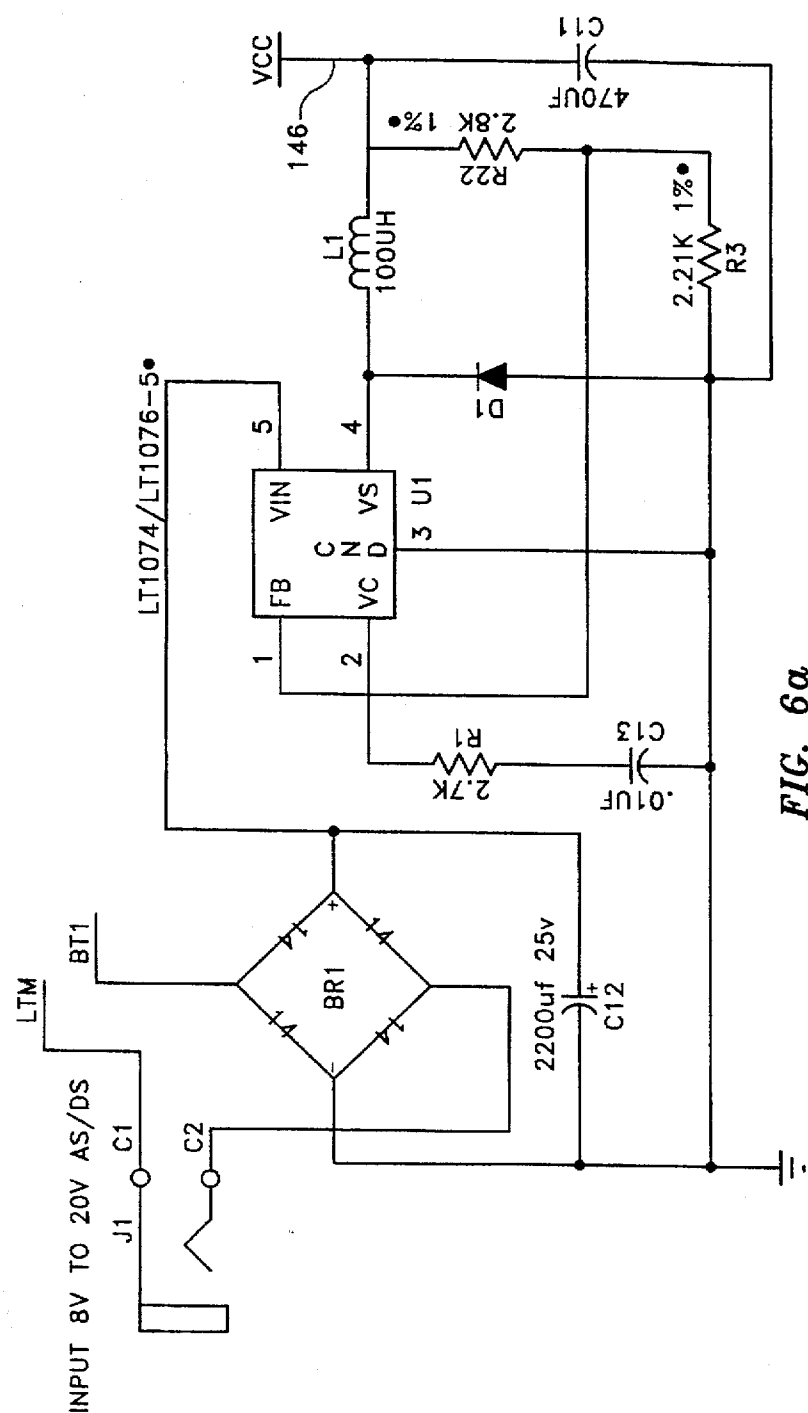


FIG. 6a

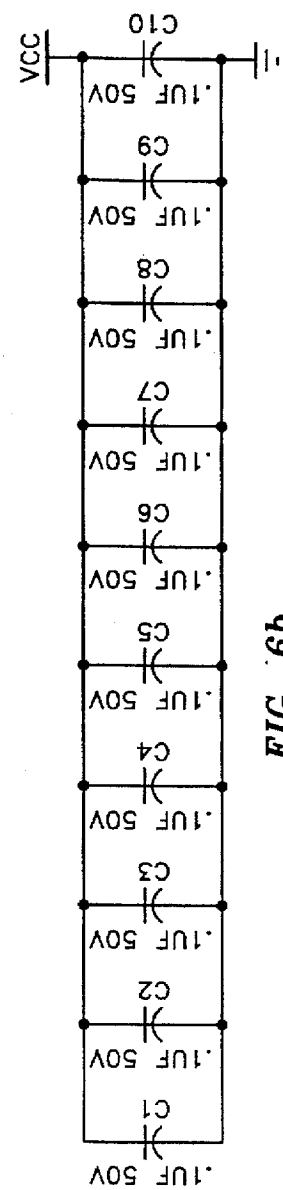


FIG. 6b

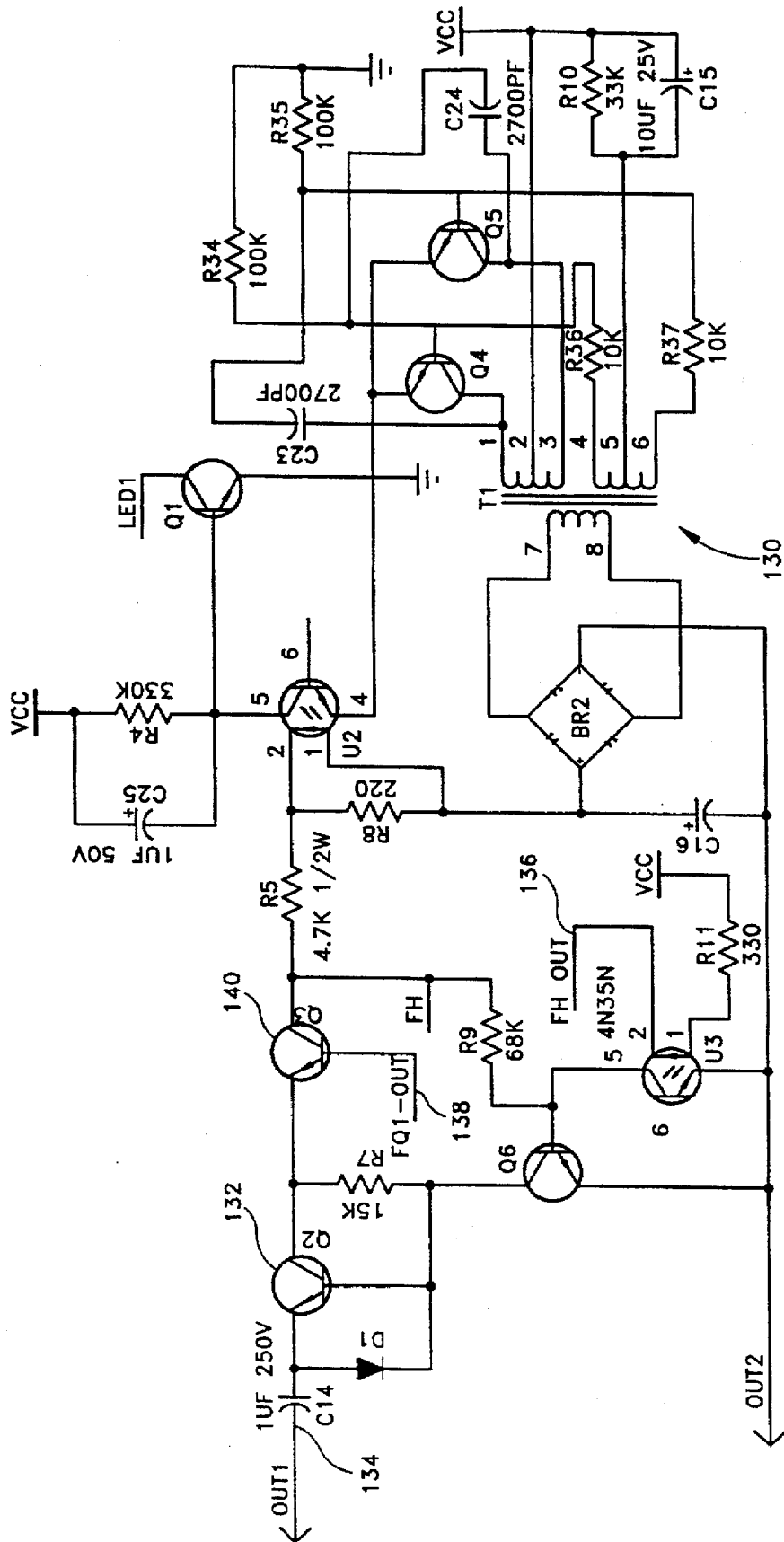
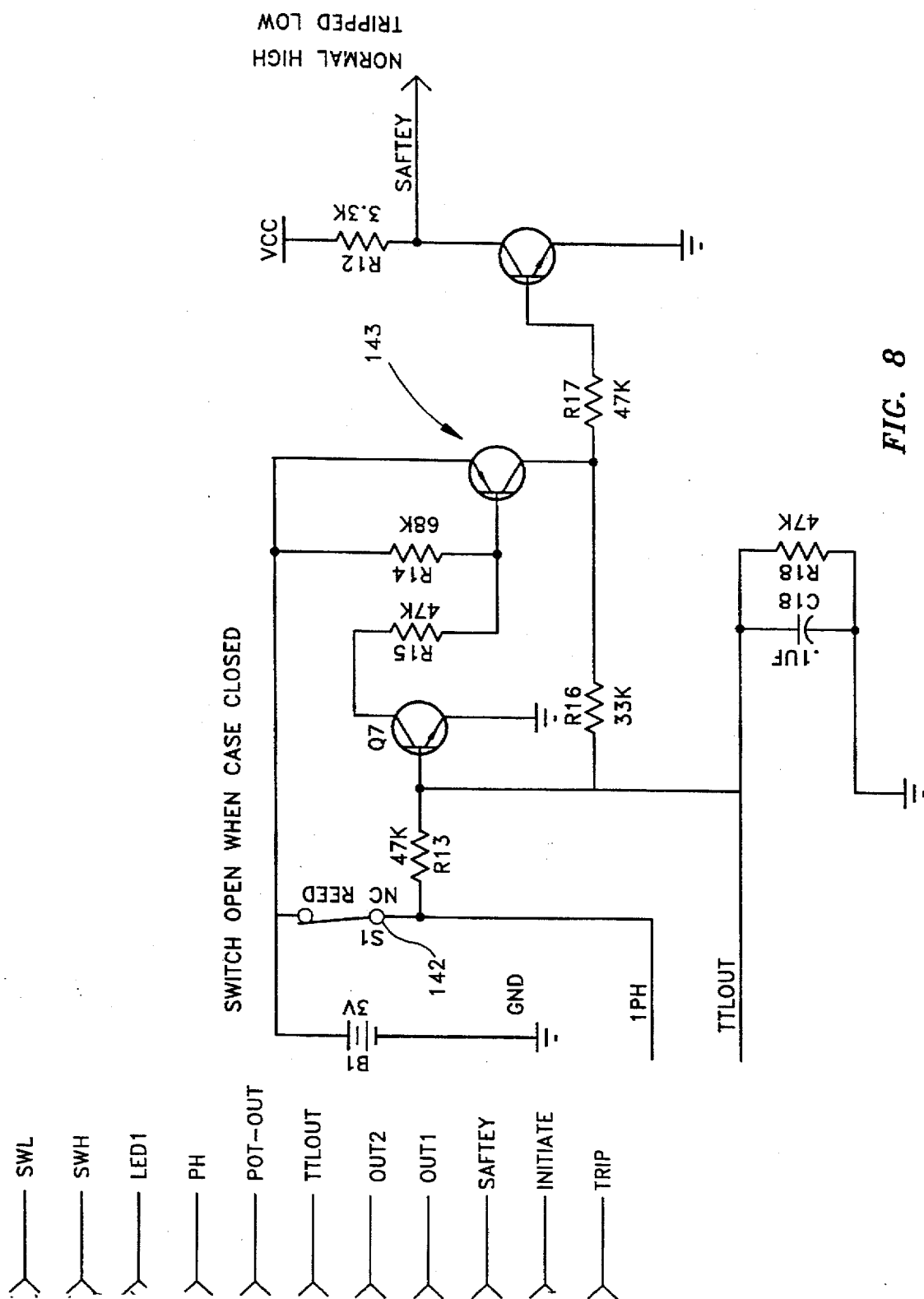


FIG. 7



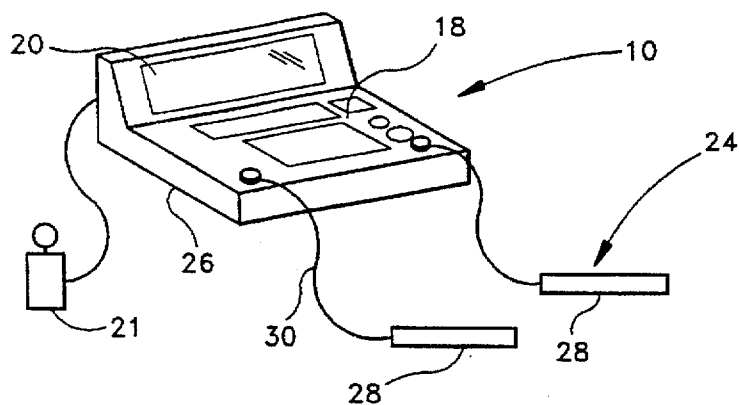


FIG. 9

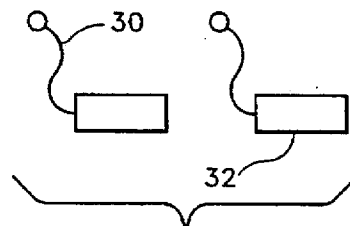


FIG. 10

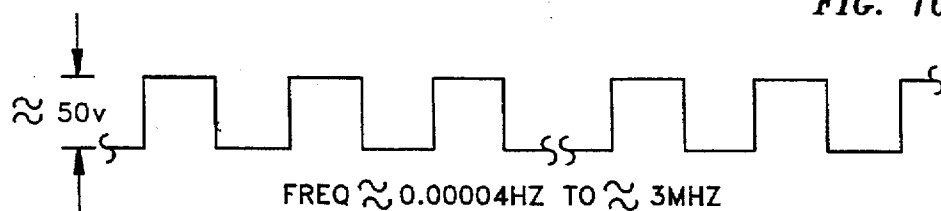


FIG. 11

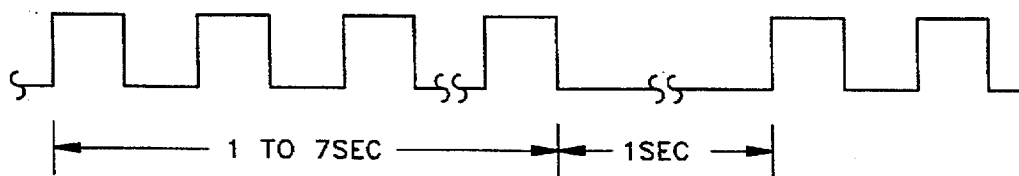


FIG. 12



FIG. 13

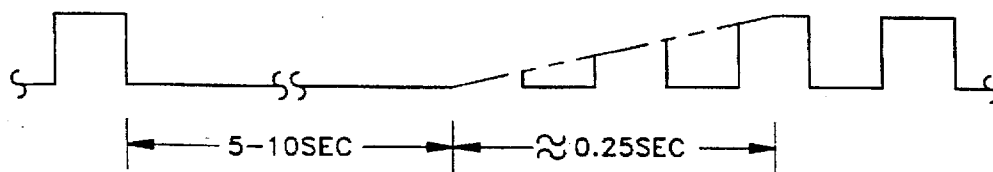


FIG. 14

BIO-ACTIVE FREQUENCY GENERATOR AND METHOD

This is a continuation of copending application Ser. No. 08/541,182 filed on Oct. 11, 1995.

BACKGROUND OF THE INVENTION

1. Field of the Invention

This invention relates in general to an improvement in the field of devices for health science applications and other applications that include industrial control of micro-organisms, and more particularly but not by way of limitation, to a device and method for generating with extremely high accuracy specific frequencies and sequences of specific frequencies that are particularly adapted for killing micro-organisms and viruses and for enhancing metabolic tissue functions in mammals.

2. Description of the Prior Art

The use of electrical stimulation in health science applications involving humans has long been known and it has been reported that one of the earliest attempts to suppress organic pain and other neurophysical effects utilizing electrical stimulation may have occurred as early as 2,000 years ago, when it was reported that gout could apparently be successfully treated by placing a diseased extremity in a tub of water containing electric eels. Apparently, attempts were also made to treat headaches with a similar treatment.

Since that early time many devices and methods have been devised to use electricity and electromagnetic fields for medical purposes. For example, U.S. Pat. No. 4,598,713 was issued on Jul. 8, 1986 to Achim Hansjurgens et al for ELECTROSTIMULATION THERAPY DEVICE AND METHOD. This device relates to electrostimulation therapy of the electrical interference therapy type. In this therapy medium frequency voltages with frequencies in the range of a few thousand Hz are applied over at least two circuits, each with one pair of electrodes, which differ from each other by a low frequency range of from a few to a few hundred Hz, so that stimulations, and thereby therapeutically effective oscillations called beats, occur with correspondingly low frequency in the body. This device stores data representative of signals having predetermined instantaneous value variations with time. These values are used to construct electrical signals having curves that represent both the medium frequency output signals and the interference beats with sufficient stability of frequency that makes possible small spacing of adjacent frequencies.

U.S. Pat. No. 5,304,486 issued to Donald C. Chang for a METHOD OF AND APPARATUS FOR CELL PORATION AND CELL FUSION USING RADIO FREQUENCY ELECTRICAL PULSE relates to the poration and fusion of cells using high-power radio frequency electrical pulses. A power function generator applies a continuous AC electrical field and/or a high-power pulsed radio frequency electrical field to induce cell congregation by the process of dielectrophoresis which porates or fuses the cells. The method can be used with a variety of cells including animal cells, human cells, plant cells, protoplasts, liposomes, bacteria and yeasts. It also can introduce into such cells a large variety of chemical agents including DNA, RNA, antibodies, proteins, drugs and inorganic chemicals.

U.S. Pat. No. 4,917,092 was issued on Apr. 17, 1990 to G. A. Todd et al for TRANSCUTANEOUS NERVE STIMULATOR FOR TREATMENT OF SYMPATHETIC NERVE DYSFUNCTION. This nerve stimulator is to be used in T.E.N.S. (Transcutaneous Electrical Nerve Stimulation)

therapy and provides for generating pulses with the pulse width, rate, and amplitude being modulated over a modulation cycle to stimulate the autonomic and central nervous system of a patient. A sympathetic nerve dysfunction is treated with a low frequency of pulses, while through the strength-duration modulation at higher pulse rates the patient's nociceptors are addressed so that the autonomic nervous system is stimulated for more than 5% of the amount of time as said central nervous system.

U.S. Pat. No. 4,509,520 was issued on Apr. 9, 1985 for an ELECTRICAL STIMULATING APPARATUS and is directed toward stimulating osteogenesis in a living body. The device is a light weight, low drain, battery operated stimulating instrument to be carried by a patient during treatment for applying a high frequency alternating current treatment signal to the skin of the patient according to a treatment program and creating a record of such treatment. A stimulating signal of 60 Hz is generated and has a square or sine waveform.

U.S. Pat. No. 4,305,390 issued to M. R. Swartz for a METHOD FOR GENERATING OXYGEN IN AN EXCITED ELECTRONIC STATE AND INACTIVATION OF MICROORGANISMS and contemplates a system for inactivation of such microorganisms and their products, for example, Herpes simplex viruses by superposition of a component such as methylene blue plus light, oxygen and electricity. This system generates the superoxide radical anion and consequently hydrogen peroxide and the hydroxyl radical, both of which enter into the inactivation process.

U.S. Pat. No. 4,428,050 discloses a tanning device for monitoring time integrated exposure to UV radiation and indicating the achievement of preselected dosages of such exposure for obtaining a desired tan.

U.S. Pat. No. 4,535,775 discloses a method healing bone fractures non-invasively by applying to electrodes coupled to the skin of a living body in the vicinity of a bone fracture an alternating voltage having a wave form that is symmetrical with respect to the axis, a frequency in the range of 20-100 KHZ and a value in the range from about 2 to 10 volts peak to peak.

It has also been found that each species of life has its own unique electronic signature. Accordingly, every micro-organism has its own specific molecular oscillation pattern. It has also been found that if such micro-organism is subjected to a specific precise electrical frequency signal at a predetermined amplitude it is possible to inactivate or kill such micro-organism while not effecting any other micro-organism or tissue. While the process by which such inactivation of the micro-organism takes place when subjected to such precise frequency signal is not fully understood, it is apparent that such procedure is efficacious. Also, it has been found that subjecting the tissue of a mammal to a specific precise frequency signal does enhance tissue regeneration. For present purposes, such precise frequency signals are termed bio-active frequencies and it is apparent that a need exists for a device to generate upon command such precise bio-active frequency electrical signals for desired applications which may include health science applications and other industrial and commercial applications. Such a device would be able to generate specific bio-active frequency signals selected by a user of the device as well as to generate specific bio-active frequency signals that were preprogrammed into the device for selection by a user for a particular health science, industrial or commercial application. It is believed that such need is fulfilled by the present invention and the particular method employed by it to generate such bio-active frequencies as desired.

SUMMARY OF THE INVENTION

Briefly stated, the present invention contemplates a system and method to generate an electrical output signal having a precise frequency and particularly adapted for use in applications for which it may be particularly suited such as health science applications and other applicable industrial and commercial applications. The particular precise frequency signal generated may be a single frequency signal that is generated for a predetermined period of time or may be a first frequency signal that is generated for a predetermined period of time and which is then followed by a second or further frequency signals that are also generated for predetermined periods of time or may be a beat frequency created by heterodyne action of two carrier frequencies as in an optional antenna array designed for use with this device. The precise frequency output signal may also be a continuous signal or may be interrupted for predetermined periods of time to provide a signal that is ON for a controllable period of time and then OFF for a constant recurring interval. The continuous precise frequency signal is believed to be particularly useful in the inactivation of microorganisms while the frequency signal that is periodically (pulsed) interrupted is believed to be useful for tissue regeneration purposes and metabolic tissue functions for mammals.

The system of the present invention for generating such precise frequency output signals includes a frequency synthesizer means that is a response to a control signal to generate a precise frequency output signal. A programmable control means is coupled to the frequency synthesizer means and is operated to generate a determinable control signal selected by a user of the system for application to the frequency synthesizer means. The control means may be programmed to generate single precise frequency signal or it may be programmed to generate a sequence of control signals for generation of a series of frequency output signals, such sequence of precise frequency output signals may be termed the SEQUENCE programs. The control means may also be preprogrammed so that a particular command from a user of the system is operative to generate a predetermined sequence of precise frequency output signals, such series of precise frequency output signals may be termed the AUTO programs. The programmable control means also stores instructions associated with a particular specific frequency which indicate whether the precise frequency output signal is to be generated as a continuous signal or whether periodic interruption intervals are to be interposed into the signal. If an interruption interval is to be interposed into the precise frequency output signal, that interval is determined to be generally around one second in duration. The period that the precise frequency output signal is ON between intervals is controllable by a user of the system and will be from about one to seven seconds.

An output circuit is coupled to the frequency synthesizer means for generation of an amplified output signal that is coupled to application means for use in a particular desired application such as health science applications and other selected industrial and commercial applications. In the instance of the application of the amplified output signal to mammals, the amplitude of the output signal may be adjusted for maximum efficiency while avoiding undue discomfort to the subject of such signal. Further, the present invention provides circuitry associated with the output signal that provides for the amplified output signal, at the start of each specific frequency, to gradually have its amplitude increased from zero to a predetermined amplification level at

a predetermined rate to avoid any inadvertent momentary discomfort to a recipient of the amplified output signal in a "soft" start of each precise frequency output signal. In addition, a no output signal interval is interposed at the start of each frequency signal and between each precise frequency signal in a sequence of frequencies so that an audio circuit may be activated to provide an auditory alert to the recipient of the amplified output signal of the ensuing change of the frequency signal.

The system is provided with a keyboard means for actuation by a user of the system to select a specific precise frequency output signal, or a sequence of precise frequency output signals or a preprogrammed sequence of frequency signals, the AUTO mode. Two separate controls are also actuable to determine amplitude level of the output signal and the ON pulse time of a precise frequency output signal that is subject to intervals of no signal. A display means is associated with the keyboard means to display information relating to actuation of the keyboard and to status of a particular output signal then being generated.

The amplified output signal generated by the output means is applied to a subject such as a mammal by spaced electrode means that are placed in direct contact with the subject in a predetermined manner. The electrode means may take several forms according to that judged most applicable for the particular application. For example, the object to be treated may also be immersed in a vat or container of a conductive fluid medium such as water, which may or may not contain electrolytes, to improve conduction of the electrical precise frequency output signal therethrough. In an application such as food processing, the use of a conductive fluid medium is particularly well suited since it permits amplification of the precise frequency signal to a high levels that achieves greater penetration than is safe or easily tolerated by living subjects such as mammals.

The precise frequency output signal is generated in a manner to ensure maximum stability and is stable to at least 0.001 Hz. Such generated precise frequency output signal is also accurate to 0.000001 Hz at the lower service range of frequencies as will be set forth hereinafter.

The system of the present invention, other than the application means, is housed in a separable enclosure which is not intended to be opened by unauthorized parties. A security circuit is provided which is operative upon opening of the enclosure to disable the system.

Other features and attendant advantages of the present invention will become apparent to those skilled in the art from a reading of the following detailed description constructed in accordance with the accompanying drawings and wherein:

BRIEF DESCRIPTION OF THE DRAWINGS

FIG. 1 is a block diagram a system constructed in accordance with a preferred embodiment of the present invention.

FIG. 2 a perspective of the display and keyboard of the block diagram of FIG. 1.

FIGS. 3a-3f is a schematic diagram of primarily the keyboard and display aspects of the system.

FIGS. 4a-4i is a schematic diagram of primarily the micro controller and programmable control portions of FIG. 1.

FIGS. 5A-B is a schematic diagram of primarily the frequency synthesizer of FIG. 1.

FIG. 6a is a schematic diagram of primarily the input power supply of the system.

FIG. 6b is a schematic diagram of another aspect of the power supply shown in FIG. 1.

5

FIG. 7 is a schematic diagram of primarily the output circuit of the system.

FIG. 8 is a schematic diagram of primarily the security circuit of the system.

FIG. 9 is a simplified perspective of the invention and an exemplary application electrode means.

FIG. 10 is a simplified perspective of an alternative application electrode means.

FIG. 11 is a waveform of one specific precise frequency output signal.

FIG. 12 is a waveform illustrating the gated ON/OFF intervals of a frequency output signal.

FIG. 13 is a waveform illustrating a soft start for the initiation of a frequency output signal.

FIG. 14 is a waveform illustrating the interval between two successive frequency output signals and the soft start of the succeeding frequency signal.

DETAILED DESCRIPTION OF THE DRAWINGS

Referring now to the drawings in detail and in particular to FIG. 1, the reference character 10 generally designates a system constructed in accordance with a preferred embodiment of the invention. The system 10 includes a frequency synthesizer means 12 that is responsive to a control signal to generate a precise frequency output signal having a frequency selected from a range of about 0.00004 Hz. to 3 Mhz. The output waveform of the frequency synthesizer preferably has a square wave output with a 50% duty cycle. The frequency generated by the frequency synthesizer is very precise since for some health science applications a variance of 1 Hz. is acceptable while for other health science applications a variance of no more than 0.001 Hz. would be desirable. Accordingly, it is desirable for the frequency range of 0.00004 Hz to 3 Mhz. contemplated to be generated by the synthesizer 12 such that a variance of less than 0.001 Hz be achieved, which requirement is met by the frequency synthesizer means shown in greater detail hereinafter. A representation of such precise frequency output signal is seen in FIG. 11.

A programmable memory control means 14 is coupled to the frequency synthesizer circuit 12 through a suitable micro controller means 16 to instruct the frequency synthesizer 12 as to the specific frequency output signal the frequency synthesizer 12 is to generate. The control means 14, which will be seen in greater detail hereinafter, includes EPROM and EEPROM integrated circuits, and is operable to store control signals representative of a single specific precise frequency output signal and a sequence of precise frequency output signals. The control signals for the programmable storage means 14 for the generation of a particular frequency or frequencies may be generated by the actuation of a keyboard means 18 which is coupled through the micro controller 16 to the programmable memory control means 14.

Actuation of the keyboard means, which is seen most clearly in FIG. 2, by a user permits the selection of a single specific frequency and the selection of a sequence of predetermined frequencies. In addition, a user of the system 10 may also choose a predetermined sequence of specific frequencies by choosing only one number which is representative of a number of such predetermined frequencies that are arranged in a predetermined sequence.

The keyboard means 18 is coupled to a suitable input/output display means 20 through the micro controller 16. The display means 20, as seen in FIG. 2, displays informa-

6

tion representative of the specific frequency output signal or sequence of output signals and a number, if chosen, that is representative of a predetermined number of specific frequency output signals that are automatically selected by selection of one representative number by a user. In addition, the display means 20 displays information relating to the time selected for a particular specific frequency signal to run and to the time remaining for a particular specific frequency signal to be generated by the system 10, either singly or as part of a chosen sequence, during the running of the system 10. The system 10 also includes an automatic frequency or channel selector 21 which is coupled to the micro controller 16. The term channel may also be used at times herein to designate a precise frequency output signal. The automatic frequency selector 21 is arranged to be placed in contact with a subject for purposes of scanning such subject and selecting automatically a particular sequence of precise frequency output signals for application to such subject. Such selection is based on sensing the frequency pattern of a subject and then determining the appropriate sequence of precise frequency output signals to be generated automatically in response to the condition sensed.

The frequency synthesizer 12 is coupled to a suitable output circuit 22 for amplification to a desired value. The output circuit 22 is then connected to a suitable application means 24 for applying the amplified specific precise frequency output signal to a subject for a desired application, which may be health science, industrial, commercial or other as desired. As seen in FIGS. 9 and 10, the system 10, except for the application means 24, is contained within a suitable enclosure 26. The application means may take the form of electrically conductive electrodes 28 in the form of elongated stainless steel cylinders that are connected to the system by suitable leads 30. The conductive cylinders are adapted to be held in the hands of a user of the system 10. The electrodes 28 may also take any suitable form such as the electrically conductive metal foil strips 32 shown in FIG. 10. Preferably, the foil strips 32 are covered in cloth coverings that are wet by a suitable electrically conductive fluid such as water to ensure that the amplified output signal properly flows through the electrode means. In such applications as may require the use of the metal foil strips as seen in FIG. 10, the strips 32 are placed in contact with the body or an extremity of a subject and are arranged in a directly opposing manner so as to position the afflicted extremity or afflicted portion of an extremity or the like directly between such electrode means. The term application means 24 should also be construed to include means to apply the specific frequency signals for such applications as industrial and commercial uses. In such instances a suitable transmitter of such frequency signals could be wire screens immersed in containers or vats of a conductive fluid, such as water, for the transmission of the frequency signals to a subject such as food which is to be sanitized or disinfected by these specific frequency signals.

The system 10 has additional affirmative features. Upon application of a control signal to the frequency synthesizer means 12 to generate a specific precise frequency output signal for application to the output circuit 22, a signal is directed to a circuit 34 which is termed a soft start circuit and which is coupled to the output circuit 22. As seen most clearly in FIG. 13, the soft start circuit 34 is operative at the initiation of a specific precise frequency signal to gradually increase the value of the amplified frequency output signal from the output circuit 22 at a predetermined rate from zero to a predetermined ultimate determinable output value. As shown in FIG. 13, this ramp up of the amplified output signal

takes place in a predetermined time period which for present purposes has been chosen to be a time period around 0.25 sec. Such a ramp up or soft start of the amplified output signal permits the subject to become quickly acclimated to the application of such signal through the output means without undue surprise to the subject. The soft start application of the amplified output signal, as seen also in FIG. 14, is interposed at the start of each specific precise frequency signal, including all frequency signals which appear in a predetermined sequence.

Another feature of the system 10 is an audio warning circuit 36 which coupled to the micro controller 16 that directs the sound circuit 36 to sound an audible sound or alarm prior to the initiation of the application of a specific precise frequency output signal to the output means 22. The purpose of the audible signal or alarm by the circuit 36 is to alert a subject to which the application means 24 have been applied of the initiation of application of the amplified output signal. As seen in FIG. 14, a time period of around 5 to 10 seconds is interposed between specific frequency signals appearing in a sequence and at the commencement of a single specific precise frequency signal to permit the auditory signal to be sounded for a desired predetermined period of time to permit a subject to which the application means 24 is applied to anticipate the commencement of the amplified output signal without surprise.

A security circuit 36 is connected to the micro controller 16. The enclosure 26 for the system 10 is separable for access to the system as may be required. To ensure integrity of the system 10, the security circuit 38 is operative upon opening of the enclosure 10 to disable the system 10 from further operation until the circuit has been reset by a an authorized person.

A suitable power supply 40 for the system 10 is included to supply requisite electrical power to all components of the system. The power supply 40 preferably receives input electrical power at 8 to 20 V DC or AC and with a Hz of 50 to 400, with at least around 1 ampere.

Referring again to FIG. 10, the actuable portion of the keyboard means 18 and the visual display portion of the display means 20 will be seen on the upper portion of the enclosure 26. The visual portion, illustrated by four single digit numbers 21, of the display means 20 displays information representative of a specific bio-active precise frequency signal either being selected by a user of the system 10 prior to starting the system to generate bio-active frequency signals or is representative of a specific precise frequency signal then being generated by the system 10. The first two digits of the three digits 23 of the display 20 also display the time selected by a user by actuation of the keyboard means 18 for a particular frequency or channel to be run. The full three digits 23 of the display 20 also display the time remaining for a specific precise frequency output signal to be generated to suitably inform the user of the system.

As previously discussed with respect of FIG. 12, an OFF signal may be generated by the system 10 to provide a desired ON/OFF cycling of certain frequency output signals with the OFF period remaining constant at around 1 sec. and the ON period of the frequency output signal being variable between about 1 to 7 sec. according to the control of a user. The period for the frequency signal to be generated, as noted above, is determined by actuation of the keyboard means 18 and would generally be considered to be around three minutes in duration.

In the illustrated embodiment of the invention, the ON pulse period for a particular frequency signal is displayed by

the single digit display 42 and such ON pulse period may be determined by a user by selective actuation of a switch 44 to increment the pulse period displayed to the desired period of time. In addition, a LED display 45 provides a visual indication to a user of the system 10 of the existence of an ON period for the precise frequency output signal.

The enclosure 26 is provided with a suitable ON/OFF switch 46 for the system 10. In the present instance, output receptacles 48 are provided for the electrode means 24 at the upper portion of the enclosure 26.

As was previously described, the amplitude of the amplified output signal is determinable by a user of the invention. A suitable amplitude control 50, such as for, example, a variable resistor, permits a user to adjust the amplitude of the amplified frequency output signal to a comfortable level commensurate with sufficient amplitude to accomplish the desired health science and/or commercial application. Generally, this amplitude is not expected to exceed a level of around 50 volts. Another LED 52 displays further information concerning the application of the amplified output signal to a subject and indicates if the electrode means 24 are not in proper contact with the subject by indicating that an open circuit condition has been sensed by the invention so that proper corrective action may be taken to ensure proper application of the electrode means 24 and correct flow of the output signal.

The illustrated portion of the keyboard means shows various keys AUTO, RUN, SEQ, DEL, SELECT, and RESET that may be actuated by a user of the invention to choose a specific frequency signal or a sequence of such signals or to choose a preprogrammed sequence of signals or to store for future access a predetermined sequence of frequency signals.

As an example of an exemplary application of the present invention, the following table sets forth a number of conditions that may be encountered in the health science field. The table also indicates a particular code or channel which is representative of a specific frequency output signal that may be applied to advantage in alleviating such condition or inactivating a particular micro-organism which is responsible for such condition. As will be noted, in some instances a plurality of frequencies or sequences of frequencies may be used to advantage.

CONDITION	REPRESENTATIVE FREQUENCY NUMBER
Adynamia, geriatric (fatigue of age)-- (applicators placed on solar plexus and behind animal's head)-- 49, 56	
Abdominal inflammation--31, 82, 4, 5, 6, 9, 10, 12 + 16, 14, 15 17, 18, 20, 28, 29 + 86, 88, 42, 45, 48, 58, 103, 104, 105, 106, 107, 112, 118, 119, 121, 125	
Abdominal pain (applied to solar plexus and belly)--78, 104, 45, 40, 57, 87, 88, 108, 109, 114, 1, 3. Also for parasites + 86, 42, 45, 48	
Abscesses--4, 6, 14, 17, 20, 39, 101	
Acidosis (urine too acid)--1, 14, 16, 17, 18, 20, 58, 88	
Acne--4, 6, 9, 10, 12, 14, 17, 20	
Actinomyces--28, 1, 17, 20, 58	
Acupuncture disturbance field (scar focus)--59	
Acute Pain--40, 57, 3, 45, 87, 88, 108, 109, 114, 1, 14, 17, 20 + 22 + 24	
Adenoids--11 + 16, 14, 17, 18, 20, 29, 58, 105, 106, 107, 109, 110, 112, 116, 117, 121, 123, 4, 6	
Adhesions--4, 6, 11 + 16, 14, 17, 18, 19, 20, 39	
Adrenal stimulant--58	
Allergy--1, 78 and 32, 87, 88, 108, 109, 114, 42, 45, 48, 58, 29 + 86	

-continued

CONDITION	REPRESENTATIVE FREQUENCY NUMBER
Alopecia (loss of hair)--1, 14, 17, 20, 28, 88	
Amenorrhea (absence of menstruation)--29, 87, 88, 108, 109, 114, 42, 45, 48, 14, 1, 3, 17, 19, 20, 58	
Anal itching (Pruritus)--28, 29, 42, 45, 48, 29 + 86, 58, 112, 87, 88, 108, 109, 114 (Secondary: 1, 14, 17, 19, 20)	
Angina (quinsy--in throat)--17, 18, 20, 28, 105, 106, 107, 109, 110, 112, 117, 118, 119, 121	
Ankylosing spondylitis--3, 45, 11 + 16, 14, 17, 18, 20, 25 + 26 + 27, 55, 82, 61, 54, 55, 69, 82, 43, 44, 49, 123, 105, 106, 109, 110, 112	
Antiseptic effect--11 + 16, 14, 17, 19, 20, 28, 29, 103, 105, 106, 107, 108, 109, 110, 112, 118, 119	
Apoplexy (stroke paralysis)--often caused by netsitter geopathic zone 53, 58, 9, 25 + 26 + 27, 42, 45, 48, 58, 29 + 86, 87, 88, 108, 109, 114, 14, 17, 20	
Appetite, lack of--28, 29 + 86, 42, 45, 48, 58, 1, 14, 17, 20	
Arteriosclerosis (hardening of arteries: regeneration takes time but is accelerated by full-body sunbaths; stopped before reddening)--1, 4, 6, 9, 10, 12, 14, 17, 18, 20, 58, 123	
Arthritis (since many different organisms can cause arthritis, many frequencies must be tested. See other types of arthritis below)--1, 11, 14, 16, 17, 20, 55, 58, 69, 78, 82	
Arthritis, rheumatoid of the muscles and tendons--34, 82, 25 + 26 + 27, 17, 20, 123	
Arthritis, arthrosis, and parathyroid disturbances affecting calcium metabolism--62	
Arthritis, arthralgia due to gout--63	
Arthritis: focal origin (gastrogenic, tonsillogenic, and syphilis or paresis)--63	
Asthma--4, 6, 9, 10, 12, 14, 17, 20, 123	
Astrocytoma--65, 67, 69, 6, 7, 14, 22, 24	
Ataxia (incoordination of muscles--slow results in some cases)--4, 6, 9, 10, 12, 14, 17, 18, 20, 25, 26, 27, 29 + 86, 42, 45, 48, 58, 121, 125 spastic ataxia--65, 67, 69	
Athlete's Foot--Try 28 first, then 11, 14, 16, 20, 58	
Autoimmune Disorders--1, 11 + 16, 14, 17, 18, 20, 25, 26, 27, 28, 34, 55, 58, 62, 63, 69, 78, 87, 88, 102	
Autointoxication--85, 87, 88, 1, 14, 17, 20, 58, 108, 109, 114	
Back pain--Best is 57 for 15 min. + 40, 3, 87, 88, 108, 109, 114. Others: 1, 19, 11 + 16, 14, 17, 20	
Bad Breath (halitosis)--11 + 16, 14, 17, 20, 58	
Bedsore--14, 11 + 16, 17, 20, 28, 58; then 82 and 47	
Biliousness--11 + 16, 1, 14, 15, 17, 20, 28	
Bladder and prostate complaints--14, 11 + 16, 17, 20, 28, 58, 63	
Boils--14, 11 + 16, 17, 20, 28, 106, 107, 109	
Bone trauma (cuts, fractures)--31, 11 + 16, 1, 14, 17, 20	
Bone/periodontal disease (also osteomalacia, etc., greatly aided by full-body sun exposure) 51, 9, 10, 25 + 26 + 27, 14, 17, 18, 20	
Bone protuberance (spurs)--82 and 34	
Brachial Neuralgia (electrodes placed on elbow and head)--85	
Breast: fibroid cysts--14, 11 + 16, 17, 18, 20 + 22 + 24	
Bright's Syndrome (Nephritis)--12, 14, 17, 20	
Bronchial asthma--85, 87, 88, 42, 45, 48, 29 + 87, 58	
Bronchitis--14, 20, 63, 64	
Bronchial Pneumonia--11 + 16, 14, 17, 18, 20	
Bubonic Plague--101, plus secondary infections: 14, 17, 20, 58	
Burns--39, 1, 14, 17, 20, 28, 38	
Bursitis (may be caused by many viruses; also experiment with arthritis frequencies)--14, 17, 20	
Cancer: see Carcinoma, Leukemia, and Sarcoma	
Candida--28	
Carcinoma--7, 14, 17, 20 + 24 + 22, 107	
Cardiac Edema (congestive heart failure)--65	
Carpal Tunnel Syndrome--11 thru 28, 71, 41, 80, 20, 87, 88, 29	
Cataract--1, 14, 17, 20	
Catarh--11 # 16, 14, 17, 20, 29, 58	
Cerebral Palsy--14, 17, 20, 87, 88	
Chicken Pox (Varicella)--11 + 16	
(and for secondary infections)--9, 10, 12, 14, 17, 20, 58	
Cholera--14, 16, 103, 14, 17, 20	
Chronic Fatigue Syndrome (generic disorder)--1, 4, 9, 14, 16, 17, 18, 20, 28, 29 + 86, 32, 41, 42, 45, 47, 48, 54, 56, 58, 71, 74, 105, 106,	
Circulatory stasis--53	
Circulation disturbance/problems--53, 63	

-continued

CONDITION	REPRESENTATIVE FREQUENCY NUMBER
5 Concentration (to increase mental)--68	
Cold in head, chest (Common cold/Adenoviruses mutate constantly; there are too many new strains to include a complete list for every cold)--1, 11, 14, 16, 17, 18, 20, 29, 58	
Colic--11, 15, 16, 17, 20, 58	
10 Colitis (irritation of colon)--11, 14, 15, 16, 104	
Constipation--11, 14, 16, 17, 18, 20, 58	
Convulsions--1, 14, 17, 20 (most commonly used)	
Contusions--66 and 43 over the focus	
Costalgia (rib pain)--40, 57, 3, 14, 17, 20, 1, 11 + 16	
Cramps--1, 14, 17, 20, 57	
15 Cystitis (of urinary bladder)--11, 14, 16, 17, 20, 28, 58	
Deafness (partial to complete)--1, 11, 14, 16, 17, 20, 58	
Dental foci (neglecting this can prevent recovery from ANY illness)--3, 45, 39, 51, 10, 11 + 16, 12, 14, 15, 17, 18, 19, 20, 24 + 25 + 27, 28, 4, 5, 9, 87, 88, 108, 109, 114	
Depression, anxiety, trembling, weakness--77	
Depression (due to exogenous circumstances--54 and 17	
20 Depression (due to drugs or toxins)--83 and 47	
Detoxification--71, 41, 80, 20, 87, 88, 29, 24, 105, 14, 17, 108, 109, 114	
Diabetes--54 + 111 + 17, and 61 + 34, 25 + 26 + 27, 42, 45, 48 + 58 + 86, 87, 105 and 53 + 63 for circulation problems. for secondary infections: 1, 4, 6, 9, 11, 14, 16, 20, 28, 58	
25 Diabetic loading--54 and 21	
Diarrhea/Dysentery--11, 14, 16, 17, 20, 28	
Distorsion (Twisting of muscles, spine)--66 and 43	
Dizziness--73	
Dupuytren's Contracture--82 and 34	
Dysmenorrhea--For rapid, lasting relief of painful menstruation: after a pure water douche, use #57 over and under uterus; leave running after pain disappears. When cause is uterine infections: 74, 11 + 16, 14, 17, 20, 28	
30 Dyspepsia (indigestion)--11, 14, 16, 17, 20	
Bars--various maladies--discharges, vertigo, ringing, hearing loss--65, 14, 17, 20, 58	
35 Eczema--65, 11, 16, 17, 20	
Eczema in vascular and lung functional disturbances--63	
Enuresis--1, 14, 17, 20, 28	
Edema--87, 88, 71 and 41, 29, 104, 14, 17, 20, 28	
Epicondylitis (pain in elbow)--82 and 34, 3, 40, 57	
Epididymitis (inflammation of testicle area/ducts)--12, 14, 17, 20, 58	
40 Epilepsy--1, 14, 16, 17, 20, 21, 25, 27, 42, 58	
Epstein Barr--105, 106, 18, 28, 14, 17, 20	
Erysipelas (skin inflammation) caused by strep pyogenes, etc.--106, 112, 1, 14, 17, 20, 28, 58	
Erythema nosodum--63	
Esophagus (congestion)--14, 17, 20	
45 Eustachian tube inflammation--11, 14, 16, 17, 18, 20, 28	
Eye inflammation (on left and right temples)--82 and 46	
Eye Disorders (blurred, cataracts, crossed, diplopia, infected, etc.)--10, 1, 14, 17, 20, 58	
Eyes (Glaucoma)--10, 14, 17, 20	
Facial paralysis--1, 14, 17, 20, etc.	
50 Fatigue--4, 9, 105, 106, 28, 42, 45, 48, 56, 58, 29 + 86	
Fever (various causes)--14, 17, 20, 58	
Fibroma . . . 7, 8, 20, 22, 24 (secondary: 11, 16, 28)	
Fibrosis of lung--56 and 37 and 30 on chest	
Fistula, Ulcer--14, 15, 17, 20	
Flashes; hot (complications)--1, 14, 17, 20	
55 Flatulence (intestinal gas)--11, 14, 16, 17, 19, 20, 28	
Flu--11, 14, 16, 17, 20, 58	
Foot--blisters--1, 14, 17, 20, 28	
Fractures--37 and 36, 1, 14, 17, 20	
Frostbite--14, 17, 20	
Frozen shoulder--1, 14, 16, 17, 20	
Fungal infection--28, 11, 14, 16, 20, 58	
60 Furunculosis herpes, skin diseases--38 and 13, 11, 16, 112; secondary 17, 20	
Gall bladder dystonia with osteitis--79, 3, 14, 17, 20, 58	
Gallstones--79, 3, 14, 17, 20, 58	
Gangrene--14, 17, 20, 58, 47	
Gas (intestinal)--11, 14, 16, 17, 19, 20, 28	
65 Gastritis and flatus--14, 15, 17, 20, 58	
Gout--63, 3, 14, 17, 20, 58	

CONDITION	REPRESENTATIVE FREQUENCY NUMBER
Gravel in urine--79, 3, 14, 17, 20, 58	
Gums (inflammation, gingivitis, pyorrhea)--11, 14, 16, 17, 20, 28, 58	
Hair loss (alopecia)--1, 14, 17, 20, 28, 88	
Hallucinations--1, 14, 17, 20, 58	
Hangover--1, 87, 88	
Hay Fever (only on some types)--14, 17, 20, 58	
Headaches--Best freq. for very rapid relief of symptoms--#40 applied near hairline. Continue for 10 min. or until pain disappears; then massage scalp in CW direction. Nevertheless, if headaches persist despite elimination of all possible causes, check the possibility of tumor. For headaches of unknown cause or toxins--61, 71, 87, 88, 74, 3, 20, 41, 42, 45, 48, 58, 61, 73, 80, 29, 23, 14, 17, 25, 26, 27, 3, 105, 106, 108, 109, 114	
For headaches caused by parasites--40, 42, 45, 47, 58, 20, 3	
Headaches: urogenitaliy caused--40, 63, 3, 108, 109, 114	
Headaches caused by vertebral misalignment-(not a substitute for chiropractic adjustment)--40, 62, 3	
Head Injuries (seek immediate medical attention)--40, 62, 3, 14, 17, 20, 87, 48, 73, 74	
Hearing problems--1, 11 + 16, 14, 17, 19, 20, 58	
Heart (lab animals only) (see Pericarditis)--46, 40, 58, 47, 75, 1, 14, 17, 20, 28, 42, 45, 58	
Heartburn--15, 4, 6, 7, 9, 10, 11 + 16, 12, 14, 17, 20, 23, 28, 29 + 86, 42, 45, 48, 58	
Hemorrhage (uncontrolled bleeding)--11, 16	
Hemorrhoids--11, 14, 16, 20	
Hepatitis (generic--placement is over and behind liver)--82 and 55, 11, 14, 16, 20	
Hernia--1, 17, 20	
Herpes (zoster)--4, 6, 9, 10, 12, 11, 16, 86 and (secondary) 14, 17, 20, 58	
High blood pressure, hypertension--1, 14, 17, 20, 65	
Hip Pain (as in coxarthritits)--3, 14, 17, 20, 40, 57	
Hives (Urticaria)--9, 14, 17, 20, 87, 88, 74	
Hoarseness--14, 19, 20	
Hyperacidity of stomach--68 on solar plexus and 36 on stomach	
Hydrocele (fluid in testicle, etc)--14, 17, 20	
Hyperosmia (taste--morbid)--58, 1, 87, 88, 124	
Hyperthyroid--78 and 85	
Hypoacidity of stomach--111, 58 over & behind pancreas	
Hypothyroid (Thyroid deficiency)--60 behind head, and 54 on thyroid	
Hypertension--58, 45, 65, 72, 1, 14, 17, 20,	
Some other varieties of hypertension: (kidney-induced, red-high, diastolic high pressure)--65 and 72	
Hypertension, spastic--45	
Impotence (many causes)--63, 7, 8, 28, 14, 16, 17, 20 + 22 + 24, 42, 45, 47, 48, 58, 25 + 26 + 27	
Indigestion--1, 14, 17, 20, 28, 29, 58, 42, 45, 48, 74	
Infantile Paralysis--12, 14, 17, 20	
Infections (many classes)--10, 11, 12, 14, 15, 16, 17, 18, 19, 20, 21, 22, 24, 25, 26, 27, 28, 29, 42, 45, 48, 58, 86, 101, 103, 104, 105, 106	
Infertility--7, 8, 28, 14, 16, 17, 20 + 22 + 24, 25 + 26 + 27, 63	
Inflammation--81	
Influenza (flu viruses mutate frequently; these codes are only a partial list)--11, 16, 12, 14, 17, 20, 58	
Insomnia--76, 78, 11, 12, 14 and 16	
Intercostal neuralgia (pain in rib musculature)--3, 40, 57, 11, 14, 16, 17, 18, 20, 42, 58, 86, 29	
Intermittent claudication (behind the head)--52 and 50	
Itching (Pruritis)--14, 17, 20, 28, 29, 42, 45, 48, 58, 29 + 86, 112	
Jaundice (liver)--10, 11, 12, 14, 16, 25, 26, 27, 29 + 86, 88, 34, 42, 45, 48, 58	
Kidney insufficiency (many viruses and toxins cause this)--61 and 53, 104, 10, 11, 12, 14, 16, 17, 20, 25, 26, 27, 28, 47, 53, 61, 29 + 86, 87, 88, 34, 41, 42, 45, 48, 58, 108, 109, 114	
Knee/Joint pain--3, 40, 57, 11, 14, 16, 17, 20, 55, 58, 69, 78, 82, 34, 82, 62, 63	
Larynx--104, 28, 29, 11, 14, 16, 17, 20, 55, 69, 78, 82, 34, 82, 62, 63	
Leprosy (more below)--27	
for infection of lesions: 10, 11, 12, 14, 15, 16, 17, 18, 19, 20, 21, 22, 24, 25, 26, 27, 28, 29, 58, 101, 103, 104 to 106	

CONDITION	REPRESENTATIVE FREQUENCY NUMBER
5 Leukemia--7, 8, 14, 17, 20, 22, 24, 107	
Leukodermia--14, 17, 20, 29, 58, 107, 112	
Locomotor Dysfunction (slow results if nerve damage exists)	
incoordination--1, 14, 17, 18, 20, 25, 26, 27, 29 + 86, 42, 45, 48, 58	
convulsions, spasticity--65, 67, 69	
10 Low Blood Pressure, hypotension--14, 17, 20, 58	
Lumbago--1, 14, 17, 20, 42, 45, 48, 29 + 86, 65, 67, 69	
Lungs--11, 16, 14, 17, 18, 20, 42, 45, 48, 58, 29 + 86, 103, 107	
Lupus Erythematosus--14, 17, 20, 18, 102	
Luxation (dislocation of organs or joints)--66 and 43	
Lymph stasis--71 and 41, 80, 20, 87, 88, 29, 105, 14, 17	
15 Malaria--108	
Measles (all)--109	
Mental Disorders (toxins usually the cause)--29, 41, 58, 71, 80, 87, 88, 1, 42, 45, 48, 74, 20, 105, 29, 11 to 28	
Menieres Disease--11 + 16, 14, 17, 20, 28, 105, 107, 109	
Meningitis--11 + 16, 14, 15, 17, 20, 25, 26, 27, 28, 29 + 86, 42, 45, 48, 58, 105, 106	
20 Menstrual problems--plain water douche first, then #57 applied above and behind uterus until pain disappears.	
Other frequencies--3, 14, 17, 20, 28, 40, 58	
Migraine--40 and 61, usually applied near hairline.	
Motion Sickness--25, 26, 27, 28, 29 # 86, 87, 88, 42, 45, 48, 58	
Mouth--eruptions; white patches (precancerous, leukoplakia)--28, 25, 7, 8, 20, 22, 24	
herpes sores . . . 5, 9, 28, 11, 12, 14, 17, 20, 102, 105	
Mucous membrane inflammation--31	
Multiple Sclerosis--use "Autoimmune" frequencies plus 14, 17, 20	
Muscles, to relax (myospasms)--19 and 70	
Muscular Dystrophy (suggested)--1, 11 + 16, 14, 17, 18, 20, 25, 26, 27, 28, 34, 55, 58, 62, 63, 69, 78, 87, 88, 102, 87, 88, 108, 109, 114	
Muscular pain (reported as applied to the focus of injury or pain for short periods, from 30 seconds to 3 minutes, depending on need)--3, 57, 40, 33, 34, 35, 46, 53, 61, 73, 80, 81, 82, 84, 85	
Muscle pain from parasites--29 + 86, 42, 445, 48, 58	
35 Mumps--1, 20, 4, 5, 7, 8, 110, 105, 14, 17, 20, 58	
Nausea--14, 15, 17, 20, 58, 74	
Nephritis nephrosis--14, 17, 20, 61, 53, 47 and 28	
Nervousness, from geopathic netzgitter zones, drugs, Prozac agitation (akathisia) . . . 78	
Nerve disorders--1, 4, 5, 6, 9, 10, 11 + 16, 14, 17, 41, 71, 80, 87, 88, 20, 25, 26, 27, 42, 45, 48, 58, 105, 106, 109	
40 Neuralgia--75, 108, 109, 114. Sometimes frequencies for nerves above	
Neurosis--55	
Nose--infection, congestion--11 + 16, 14, 17, 18, 20, 29, 104, 28, 58	
Numbness--1, 4, 5, 6, 9, 10, 11 + 16, 14, 17, 20, 25, 26, 27, 104, 106, 109	
45 Oral Lesions--4, 5, 8, 9, 10, 11 + 16, 14, 17, 18, 19, 20, 28, 29, 87, 88, 108, 109, 114	
Orchitis (Inflammation of testes due to TB, mumps, gonorrhea, cancer, etc.)--4, 5, 6, 7, 8, 9, 10, 11 + 16, 12, 14, 15, 17, 18, 20 + 22 + 24, 25, 26, 27, 42, 45, 48, 58, and 110	
Osteomyelitis (bone infection)--79, 4, 5, 6, 7, 8, 9, 10, 11 + 16, 12, 14, 15, 17, 18, 20 + 22 + 24	
50 Osteomalacia--full body sunbathing (stop before skin reddens)	
Otosclerosis (type of deafness)--65	
Ovarian disorders--25 + 26 + 27, 28, 29, 57, 4, 5, 6, 7, 8, 9, 10, 11 + 16, 12, 14, 15, 17, 18, 20 + 22 + 24, 58	
Pain (look under appropriate affliction for more)--3, 45, 40, 57	
Pain of infection--3, 40, 57, 45, 14, 11 + 16, 17, 18, 20, 55, 74	
Pain of cancer--3, 40, 57, 45, 7, 8, 20 + 22 + 24	
Pancreatic insufficiency--111, 61 over pancreas and 34 under (behind?) pancreas,--(secondary)--25 + 26 + 27, 28, 29, 57, 4, 5, 6, 7, 8, 9, 10, 11 + 16, 12, 14, 15, 17, 18, 20 + 22 + 24, 58	
Paralysis, spastic (slow results in some cases)--1, 14, 17, 18, 20, 25, 26, 27, 29 # 86, 42, 45, 48, 58, 65, 67	
60 Paralysis, nonspastic (langorous) (slow results in some cases)--1, 14, 17, 18, 20, 25, 26, 27, 29 # 86, 42, 45, 48, 58, 65, 67	
Parasites--42, 45, 48, 58, 29 + 86	
Pelvic Inflammatory Disease (PID)--4, 5, 6, 7, 8, 9, 10, 11 + 16, 17, 18, 20 + 22 + 24, 25, 26, 27, 28, 29 + 87, 45, 48, 103, 105, 110, 107, 118	
65 Pericarditis--4, 6, 10, 14, 11 + 16, 17, 19, 20, 26, 28, 40, 42,	

-continued

CONDITION	REPRESENTATIVE FREQUENCY NUMBER
47, 57, 58, 75, 45, 46, 47, 48, 29 + 86	
Pharyngitis--104, 107, 31, 10, 11 + 16, 18, 17, 14, 20, 58, 87, 88,	
Pleurisy--11, 16, 14, 17, 18, 20, 42, 45, 48, 58, 29 + 86, 103, 107	
Pneumonia--11, 16, 14, 17, 18, 20, 58, 103, 107	
Poliomyelitis--121,	
secondary complications--11 + 16, 105, 12, 14, 17, 20	
Polyps--4, 5, 6, 7, 8, 9, 10, 20 + 22 + 24, 25 + 26 + 27, 28, 29, 58, 87, 88	
Pre-op and post operation--prevention and control of nosocomial (hospital-acquired) and idiopathic infection--6, 9, 10, 11 + 16, 12, 14, 15, 16, 17, 18, 20, 28, 29, 87, 88, 105, 107, 109	
Prostatitis, benign prostrate tumor--44 and 30, 87, 88, 4, 5, 6, 7, 8, 11 + 16, 17, 18, 20 + 22 + 24, 28, 42, 45, 48, 58, 29 + 87	
Prostate complaints--63, 7, 8, 20, 22, 24, 28, 14, 16, 17, 20, 42, 45, 47, 48, 58	
Pruritis (itching)--14, 17, 20, 28, 29, 42, 45, 48, 58, 29 + 86, 112	
Prostate tumor (malignant) one electrode against navel, other against rectum--7, 8, 20 + 22 + 24	
Psychosomatic pain (viral, energetic meridian blocks, toxicity, nutritional)--20, 40, 57, 41, 44 + 32, 87, 88, 108, 109, 114	
Pyorrhea (periodontal disease: jawbone infection; receding gums)--4, 5, 8, 9, 10, 11 + 16, 14, 17, 18, 20, 28, 29, 87, 88	
The electrodes are placed so the infection lies on a straight line between them.	
Rabies--122	
Raynaud's Disease (cold hands & feet)--20, 58	
Rheumatoid arthritis of the muscles and tendons--82 and 34, 62, 63, 25 + 26 + 27, 1, 11 + 16, 14, 17, 20, 55, 58, 69, 78, 82	
Rhinitis or Runny nose--1, 11 + 16, 12, 14, 17, 18, 20, 29, 58, 87, 88, 28	
Sarcoma--7, 8, 14, 17, 20, 22, 24, 107	
Scars--59 + 11 to 28	
Scarlet Fever--14, 17, 20, 22 + 24	
Sciatica or Ischia	
(severe cases require direct application of electrodes over and behind the afflicted area at a high intensity)--11 + 16, 14, 17, 20, 22, 24, 61	
Sedative Effect (reported use on bleeding, bruises, insomnia, sinusitis (also reported use on lymph stasis/edema, auricular treatment)--80	
Sexual Dysfunction--63, 7, 8, 28, 14, 16, 17, 20 + 22 + 24, 42, 45, 47, 48, 58, 25 + 26 + 27	
Shingles #Herpes Zoster--11 + 16, 9, 86 (secondary) 14, 17, 20, 57	
Sinusitis--11 + 16, 14, 17, 20, 42, 48, 87	
Sleeping Sickness--124	
Slipped discs (includes spasms from microbial toxins)--57 (BEST), 42, 14, 17, 109, 17, 20, 45, 48, 87, 88, 108, 109, 114	
Smallpox--125 plus 11, 16, and (secondary) 14, 17, 20	
Sneezing--14, 17, 20, 28, 88	
Sore Throat (pharyngitis) (consider also food allergies)--4, 5, 9, 10, 11 + 16, 14, 17, 18, 20, 28, 107, 123	
Spasms, muscle--70	
Spastic paresis--25, 26, 95 and 50	
Spleen, enlarged--54 and 17	
plus secondary--1, 4, 6, 9, 11, 14, 16, 20, 28, 58	
Spondylitis, acute--82 and 61, or 69 and 55 on the focus	
Staphylococci infection--20	
Stiff neck--74	
Spastic stiff neck--65 and 72	
Stiff muscles in general--33, 34, 35, 40, 42, 46, 53, 58, 61, 73, 80, 81, 82, 84, 85, 9, 11 + 16, 14, 17, 18, 20	
Stomach disorders--7, 8, 14, 17, 20, 20 + 22 + 24, 42, 45, 48, 58, 75, 103, 118, 119	
Stones--77	
Streptococci infection--14	
Streptothrix infection--17	
Sty--1, 14, 17, 20, 58	
Subluxation induced disorders--62	
Sun Allergy (consider also that prescription drugs, such as Psoralen, often produce a "sun allergy" due to phototoxic effects on the liver)--78 and 32	
Sunstroke--29, 105, 39, 3, 45, 87, 88, 14, 58	
Surgery--Prevention and control of nosocomial (hospital-acquired) and idiopathic infection--6, 9, 10, 11 + 16, 12, 14, 15, 16, 17,	

-continued

CONDITION	REPRESENTATIVE FREQUENCY NUMBER
5 18, 20, 28, 29, 87, 88, 105, 107, 109 surgical pain, post-op recovery--3, 45, 40, 57, 44 + 32	
Swelling (Edema)--87, 88, 71 and 41, 29, 105, 14, 17, 20	
Swollen glands--14, 17, 20, 58, etc.	
Syphilis--25 + 26 + 27	
Tachycardia (also used for pain of arthritis, headache, Tui-na, and facial toning)--82	
Teeth (pain) Also see Pyorrhea.	
This can prevent recovery from other illnesses: 3, 45, 39, 51, 10, 11 + 16, 12, 14, 1, 5, 17, 18, 19, 20, 24 + 25 + 27, 28, 4, 5, 9	
Tendomyopathy (applied to the focus of injury or pain for short periods, from 12 seconds to 3 minutes, depending on need)--3, 33, 34, 40, 46, 53, 57, 58, 61, 73, 80, 81, 82, 84, 85, 87, 88	
Tetanus--125, 27	
(for secondary complications)--14, 17, 20	
Thrush (aphtha, sprue, stomatitis)--28	
Tonsillitis--82 and 47, 11 + 16, 12, 14, 15, 17, 18, 20, 25 + 26 + 27, 28, 123, 107	
20 Toothache (hidden dental and jaw infection or foci will prevent recovery from ANY illness)--3, 45, 39, 51, 4, 5, 9, 10, 11 + 16, 12, 14, 15, 17, 18, 20, 24, 25, 27, 28, 87, 88, 108, 109, 114	
Tooth extraction, afterward--3, 45, 51, 68	
Toxins--29, 41, 87, 88, 108, 109, 114. 1, 14, 17, 20, 42, 45, 48, 58	
25 Trauma--3, 45, 39, 19, 14, 17, 20, 28	
Trench Mouth--25 + 26 + 27	
secondary complications--14, 17, 18, 20, 25, 26, 28	
Trigeminal neuralgia--4, 5, 6, 9, 10, 11 + 16, 14, 15, 17, 18, 19, 20, 25, 88, 68, 56, 105 to 107, 109, 114, 116, 112, 121, 122, 125, 124	
30 Tuberculosis--11 + 16, 20 + 22 + 24	
secondary complications--18, 7, 8, 28, 107	
Typhoid--118, 119, 11 + 16, 22, 9	
Ulcers--5, 6, 7, 9, 10, 14, 15, 16, 17, 18, 20, then 82 and 47	
Urethritis--4, 6, 7, 9, 10, 11 + 16, 12, 14, 15, 17, 18, 20, 25 + 26 + 27, 28, 29 + 86, 42, 45, 48, 82	
35 Urticaria (Hives)--(often due to parasite toxins: condition becomes worse until all toxins are excreted)--9, 29, 42, 45, 48, 87, 88, 17, 20, 71	
Varicoses--82 and 55	
Vegetative Dystonia (involuntary muscle dysfunction)--53	
40 Vein thrombosis (blood clot)--23	
Vertigo, giddiness of unknown cause--49 and 73	
Warts--4, 5, 6, 7, 8, 9, 10, 12, 20 + 22 + 24	
Worms (see note under parasites)--42, 45, 48, 58, 29 + 86	
Wound healing--14, 17, 20, 37, 39, 58	
Wound healing, Delayed--53 (plus those above)	
45 Yeast infection (<i>Candida albicans</i> , etc.)--28	
Yellow Fever--114	
MISCELLANEOUS NUMBERS	
50 Requiring special training (acupuncture, aesthetology, biophysics, etc.)	
SCHUMANN FREQUENCY,	
used for entraining the brain in psychic healing experiments, and to scan brain for troubled areas--68	
55 for RELAXATION, MEDIATION AND DEEPER SLEEP--74	
To stimulate MENTAL CLARITY--60	
Facial toning--62	
Kidney meridian (balancing/correction)--65	
60 German and Russian studies have discovered a very high statistical correlation between the personalities of individuals and the position of planets at the time of their birth.	
Dr. Sieger believes this is due to the effect of gravitational interference fields. Two frequency patterns which correlate	
65 with the fields in that time frame are:	
Sagittarius--close to #17, and 38	
Capricorn and Aquarius--37	

As noted above, the table sets forth certain codes or channels that are representative of precise frequency signals. The precise frequency signals represented by such channels are set forth in the table below which also sets forth the minimum time in seconds that the precise frequency signal is gated ON between one second periods that the frequency signal is gated OFF.

CODE	FREQUENCY	GATE	
Auto 1		0	
1	10,000	0	
2	1050	0	
3	3040	0	
4	2720	0	
5	2489	0	5
6	2170	0	
7	2127	0	
8	2008	0	
9	1800	0	
10	1600	0	10
11	1550	0	
12	1500	0	
13	1000	0	
14	880	0	
15	832	0	
16	802	0	
17	787	0	
18	776	0	
19	760	0	
20	727	0	
21	700	0	
22	690	0	
23	685	0	
24	666	0	
25	650	0	
26	625	0	
27	600	0	
28	465	0	
29	444	0	
30	410	3	
31	380	3	
32	330	3	
33	320	3	
34	250	3	
35	240	3	
36	230	3	
37	220	3	
38	200	3	
39	190	3	
40	160	3	
41	148	3	
42	125	3	
43	110	3	
44	100	3	
45	95	3	
46	80	3	
47	73	3	
48	72	3	
49	60	3	
50	48	3	
51	47.5	3	
52	45	3	
53	40	3	
54	35	3	
55	28	3	
56	27.5	3	
57	26	3	
58	20	3	
59	18	3	
60	12	3	
61	10	3	
62	9.6	3	
63	9.4	3	
64	9.35	3	
65	9.2	3	
66	9.1	3	
67	8.25	3	
68	7.83	3	

-continued

CODE	FREQUENCY	GATE
69	7.7	3
70	6.8	3
71	6.3	3
72	6.0	3
73	5.8	3
74	4.9	3
75	3.9	3
76	3.6	3
77	3.5	3
78	3.0	3
79	2.65	3
80	2.5	3
81	1.5	3
82	1.2	3
83	1.1	3
84	1.0	3
85	.5	3
86	1865	0
87	522	3
88	146	3
89	610	3
90	2005	0
91	2025	0
92	.6	0
93	.75	3
94	.4	3
95	30.87	3
96	32.7	3
97	36.71	3
98	41.2	3
99	43.65	3
100	49	3
101	500	0
102	1850	0
103	450	0
104	440	0
105	428	0
106	660	0
107	589	0
108	555	0
109	333	0
110	14	0
111	15	3
112	2000	0
114	999	0
115	.67	0
116	770 + wobble	0
117	780 + wobble	0
118	1570	0
119	1770	0
Auto 2	804.1	3
120	804.1	3
121	805.6	3
122	807.23	3
123	808.79	3
124	810.35	3
125	811.91	3
126	813.48	3
126	815.04	3
127	816.60	3
128	818.16	3
129	819.73	3
130	821.288	3
Auto 3	822.75	
131	822.75	
132	824.22	
133	825.69	
134	827.16	
135	828.63	
136	830.10	
137	831.57	
138	833.04	
139	834.51	
140	835.98	
141	837.463	
Auto 4	838.85	
142	838.85	
143	840.24	

-continued

CODE	FREQUENCY	GATE
144	841.63	
145	843.02	5
146	844.41	
147	845.80	
148	847.19	
149	848.58	
150	849.96	
151	851.35	10
152	852.749	
153	854.07	
154	855.38	
155	856.71	
156	858.02	
157	859.34	15
158	860.66	
159	861.98	
160	863.30	
161	864.62	
162	865.94	
163	867.254	
Auto 6	868.51	20
164	868.51	
165	869.77	
166	871.02	
167	872.28	
168	873.53	
169	874.79	25
170	876.04	
171	877.30	
172	878.55	
173	879.81	
174	881.07	
Auto 7	882.265	30
175	882.265	
176	883.465	
177	884.664	
178	885.863	
179	887.063	
180	888.262	35
181	889.461	
182	890.661	
183	891.860	
184	893.060	
185	894.259	
Auto 8	895.407	40
186	895.407	
187	896.555	
188	897.703	
189	898.851	
190	899.999	
191	901.147	
192	902.295	45
193	903.443	
194	904.591	
195	905.739	
196	906.887	
Auto 9	907.989	50
197	907.898	
198	909.091	
199	910.193	
200	911.295	
201	912.397	
202	913.499	
203	914.602	55
204	915.704	
205	916.806	
206	917.908	
207	919.010	
Auto 10	920.071	60
208	920.071	
209	921.132	
210	922.192	
211	923.253	
212	924.314	
213	925.374	
214	926.435	
215	927.496	65
216	928.556	

-continued

CODE	FREQUENCY	GATE
217	929.617	
218	930.677	
Auto 11	960.211	
219	960.211	
220	989.744	
221	1019.27	
222	1048.81	
223	1078.34	
224	1107.88	
225	1137.41	
226	1166.94	
227	1196.48	
228	1226.01	
229	1255.54	
230	1285.08	
231	1314.61	
232	1344.14	
233	1373.68	
234	1403.21	
235	1432.74	
236	1462.28	
237	1491.81	
238	1521.34	
239	1550.88	
240	1580.41	
241	1609.94	
242	1639.48	
243	1669.01	
244	1698.54	
245	1728.08	
246	1757.61	
247	1787.14	
248	1816.68	
Auto 12	646711	
249	646711	
250	647640	
251	648569	
252	649498	
253	650427	
254	651356	
255	652285	
256	653214	
257	654143	
258	655072	
259	656002	
260	656931	
261	657859	
262	658789	
263	659718	
264	660647	
265	661576	
266	662505	
267	663434	
268	664363	
269	665292	
270	666221	
271	667150	
272	668079	
273	669009	
274	669938	
275	670887	
276	671796	
277	672725	
278	673653	
279	674583	
Auto 13	674583	
280	674583	
281	675477	
282	676371	
283	677265	
284	678159	
285	679053	
286	679947	
287	680841	
288	681735	
289	682629	
290	683523	
291	684417	

-continued

CODE	FREQUENCY	GATE
292	685311	5
293	686206	
294	687100	
295	687994	
296	688888	
297	689782	10
298	690676	
299	691570	
300	692464	
301	693358	
302	694252	15
303	695146	
304	696040	
305	696934	
306	697828	
307	698722	20
308	699616	
309	700510	
310	701404	
Auto 14	702265	
311	702265	25
312	703126	
313	703987	
314	704848	
315	705709	
316	706570	30
317	704432	
318	708292	
319	709153	
320	710015	
321	710876	35
322	711737	
323	712598	
324	713459	
325	714320	
326	715181	40
327	716042	
328	716903	
329	717764	
330	718625	
331	719486	45
332	720347	
333	721209	
334	722070	
335	722931	
336	723792	50
337	724653	
338	725514	
339	726375	
340	727236	
Auto 15	728067	55
341	728067	
342	728898	
343	729730	
344	730562	
345	731393	60
346	732225	
347	733056	
348	733888	
349	734719	
350	735551	65
351	736382	
352	737214	
353	738045	
354	738876	
355	739708	60
356	740539	
357	741371	
358	742202	
359	743034	
360	743865	65
361	744697	
362	745528	
363	746360	
364	747191	
365	748022	65
366	748854	
367	749686	

-continued

CODE	FREQUENCY	GATE
368	750517	Auto 16
369	751349	
370	752180	
371	752985	
372	752985	
373	753790	Auto 16
374	754594	
375	755399	
376	756203	
377	757009	
378	757813	Auto 16
379	758618	
380	759423	
381	760227	
382	761032	
383	761836	Auto 16
384	762641	
385	763446	
386	764251	
387	765055	
388	765860	Auto 16
389	766665	
390	767470	
391	768275	
392	769079	
393	769884	Auto 16
394	770688	
395	771494	
396	772298	
397	773103	
398	773907	Auto 16
399	774712	
400	775517	
401	776322	
402	777103	
403	777103	Auto 16
404	777883	
405	778663	
406	779444	
407	780224	
408	781005	Auto 16
409	781785	
410	782565	
411	783346	
412	784126	
413	784906	Auto 16
414	785687	
415	786468	
416	787248	
417	788028	
418	788809	Auto 16
419	789590	
420	790370	
421	791151	
422	791931	
423	792711	Auto 16
424	793492	
425	794272	
426	795053	
427	795833	
428	796613	Auto 16
429	797394	
430	798175	
431	798955	
432	799735	
433	800493	Auto 16
434	800493	
435	801251	
436	802010	
437	802768	
438	803526	Auto 16
439	804284	
440	805042	
441	805801	
442	806559	
443	807317	Auto 16
444	808075	
445	808833	

-continued

CODE	FREQUENCY	GATE
443	809592	
444	810350	
445	811108	
446	811866	
447	812624	
448	813383	
449	814141	
450	814899	
451	815657	
452	816415	
453	817174	
454	817932	
455	818690	
456	819448	
457	820207	
458	820965	
459	821723	
460	822481	
Auto 19	823217	
461	823217	
462	823953	
463	824688	
464	825424	
465	826160	
466	826896	
467	827632	
468	828367	
469	829103	
470	829839	
471	830575	
472	831310	
473	832046	
474	832782	
475	833518	
476	834254	
477	834989	
478	835725	
479	836461	
480	837197	
481	837932	
482	838668	
483	839404	
484	840140	
485	840876	
486	841611	
487	842347	
488	843083	
489	843819	
490	844554	
Auto 20	845276	
491	845276	
492	845997	
493	846718	
494	847439	
495	848160	
496	848881	
497	849603	
498	850324	
499	851045	
500	851776	
501	852487	
502	853208	
503	853929	
504	854651	
505	855372	
506	856093	
507	856814	
508	857535	
509	858256	
510	858978	
511	859699	
512	860420	
513	861141	
514	861862	
515	862583	
516	863304	
517	864026	
518	864747	

-continued

CODE	FREQUENCY	GATE
519	865468	
520	866189	
521		
522		
523		

Referring now to FIG. 5, it will be seen that the frequency synthesizer means 12 generally includes an integrated circuit 54 provided by a Q220 DDS IC made by Qualcomm. The input data for the control word for the frequency synthesizer circuit 54 is input by means of three eight bit latches 56, 58 and 60. The 23 bit frequency control word provided by the latches 56, 58, and 60 instructs the frequency synthesizer circuit to generate a specific frequency having a frequency resolution of 0.715 Mz. on output 62. Resolution of the frequency signal thus generated is increased further by dividing it by a first Divide by 100 counter 64, a 28C64-15/L counter made by Phillips, and then by a second similar Divide by 100 counter 66. The output 68 of counter 66 is then input to a Divide by 2 flip-flop 70 in FIG. 4 to ensure a 50% duty cycle for the resultant precise frequency output signal on output 72. The precise frequency output signal thus has a range, in the illustrated embodiment, of around 0.00004 Hz to 3 Mhz.

Referring now to FIG. 3a, it will be seen that actuation of a key of the keyboard set 18 provides a particular interrupt control signal to the controller 74 where it is stored. The controller is an integrated circuit Model 74C923 made by National Semiconductor Co. The control signal stored in controller 74 is then read by the microcontroller in FIG. 4. The microcontroller 16 is an Intel P80C31BH IC controller. The microcontroller 16 operates through a bus controller 76, a Phillips PCF8584 controller, which has the ability to operate in both master and slave modes to communicate with the 4 digit display 21 in FIG. 3a to display information representative of the particular frequency output signal that has been chosen by a user of the system 10 or which is presently running on the system. The microcontroller 76 also communicates similarly with the 4 digit display 23 in FIG. 3d to display the time remaining for a specific frequency signal to continue to be generated and the particular pulse ON period for that specific frequency signal. The 4 digit display 23 is provided by display circuits 78 and 80 under the control of a suitable display controller 82, a Phillips SAA1064 controller. The display circuits are HDSP-A011 units made by Hewlett Packard.

The display 23 is provided by Hewlett Packard HDSP-A011 units 84, 86, 88, and 90 under the control of a display controller 92, a Phillips SAA1064 controller. The LED 52 and the LED 44 are activated under the control of the microcontroller 16, as seen in FIG. 1. Also, in FIG. 3f it will be seen that the pulse time switch 44 of the keyboard means and the voltage control 50 of the keyboard communicate with the microcontroller 16 in FIG. 1.

Referring again to FIG. 4b, the programmable memory 14 of the system 10 is provided by a first programmable memory 94, an EPROM model NM27C512Q-150 made by National Semiconductor, and a second programmable memory 96, an EEPROM model 28C64A-15/L made by Microchip. The first programmable memory 94 stores the control signals for the preprogrammed sequences of frequency signals that may be accessed by a user actuating a particular control number or channel number that is representative of a particular sequence of frequencies. The programmable memory 96 contains control signals that may be

accessed by a user when selecting a particular specific precise frequency signal or when creating a particular customized sequence of specific frequency signals for application.

A first decoder 98 and a second decoder 100 communicate with the programmable control means 14 and cooperate with a suitable digital to analog converter circuit 102 to provide analog signals to the R-C circuit 104 formed by resistance 106 and capacitor 108 to provide the one second delay between ON periods in certain specific frequency signals as heretofore described.

The circuit 110 includes the resistor 112 and the capacitor 114 and provides the soft start circuit 34 that upon commencement of a specific precise frequency signal gradually ramps up the output signal from zero to a predetermined value at a constant predetermined rate.

A speaker 116 is included in the audio circuit 118 to provide an auditory signal to a user of the system 10 upon the initiation of each specific precise frequency signal whether singly or part of a sequence. The decoder 120 provides address DEMUX and the integrated circuit 122 acts as a receiver to clean up signals received from the keyboard to preclude erroneous signals from key bounce for example. Reference character 124 refers to an integrated circuit that provides multiplexing of the data bus for the particular integrated circuits chosen for implementation of the system 10. The circuit generally designated by 126 controls the slew rate of the pulses. A suitable user port 128 is provided for further applications of the amplified output signal beyond the particular spaced electrode means illustrated in the previously described Figures.

Referring now to FIG. 7, the output circuit 22 is illustrated. A transformer coupled inverter 130 generates around 100 V DC. Isolation of the output is accomplished by means of the transformer 17 and optically coupling signals to the output circuit which includes the emitter follower 132 and the output capacitor 134. The specific precise frequency signal including such OFF intervals as may be interposed appears on lead 136 and is optically coupled to the emitter follower 132 for control thereof. The output level of the amplified frequency signal is determined by the control signal 138 established by the intensity control 50 and for present purposes may be considered not to exceed a level of around 50 volts. The control signal 138 controls a series pass transistor 140.

FIG. 8 illustrates the security circuit 22 and includes a normally open reed relay 142 that is configured in such a manner so as to detect any opening of the enclosure 26. Any momentary contact of the relay 142 causes a flip-flop 143 to be set to Low when tripped. A Low signal for the flip-flop 144 is directed to the microcontroller 16 to disable the system 10 until the system is reset by an authorized party.

FIGS. 6a-6b illustrates the power supply 40 which supplies power to all components of the system 10. The power supply is a conventional buck converter capable of operation from 7 to 20 V DC or AC. The voltage VCC appears on lead 146 and is +5 V. FIG. 6b is a further depiction of charging capacitors to maintain VCC during usage.

From the foregoing, it is apparent that the present invention provides a novel system and method for generating a plurality of specific precise bio-active frequencies and sequences of such frequencies for appropriate applications such as health science, industrial and other commercial uses. The term "health science applications" should be considered broadly and not limiting as the bio-active frequency generator of the present invention may also be used for such applications as destroying micro-organisms such as are found in water supplies and in various oil supplies.

Changes may be made in the combination and arrangement of parts or elements as heretofore set forth in the specification and shown in the drawings, it being understood that changes may be made in the precise embodiment disclosed without departing from the spirit and scope of the invention as defined in the following claims.

What is claimed is:

1. A system to generate an output signal having a predetermined precise frequency and particularly adapted for use in applications including health science and industrial uses, comprising:

frequency generation means responsive to a control signal to generate an output signal having a predetermined precise frequency within a frequency range from at least 0.00004 Hz to 3 Mhz;

said frequency generation means includes a frequency synthesizer means, and a programmable memory means coupled to a keyboard means whereby actuation of the keyboard means to generate a signal representative of a predetermined frequency selected from the memory means for application to the frequency synthesizer means to cause such frequency synthesizer means to generate the selected predetermined precise frequency,

keyboard means coupled to the frequency generation means and actuable to generate a control signal representative of a predetermined frequency; said keyboard means is actuable to store within the memory means instructions representative of a sequence of predetermined precise frequencies whereby a number of frequencies can be generated serially upon command,

output means coupled to the frequency generation means to amplify the voltage of the generated frequency signal to a predetermined level, and

application means coupled to the output means for application of the output signal for applications including health science and industrial uses the shape of said output signal is a square wave having a 50% duty cycle.

2. The invention as defined in claim 1 further comprising gating means coupled to the frequency generation means and responsive to the specific precise frequency being generated to permit selective control of the generated frequency signal.

3. The invention defined in claim 1 further comprising display means coupled to the keyboard means for displaying information representative of actuation of the keyboard means and the output signal.

4. The invention as defined in claim 1 wherein the frequency of the output signal has an accuracy of at least 0.001 Hz.

5. The system as defined in claim 1 wherein said frequency of said generated signal comprises a selected frequency in cycles per second from the group of specific frequencies for inactivating micro-organisms and viruses which may be present and enhancing specific tissue regeneration in a mammal, said group of frequencies and the purpose therefor consisting of Adynamia, geriatric 60, 27.5; Abdominal inflammation 380, 1.2, 2727, 2489, 2170, 1800, 1600, 2302, 880, 832, 787, 776, 727, 465, 444, 1865, 146, 125, 95, 72, 20, 450, 440, 428, 660, 589, 2000, 1570, 1770, 805.6, 811.91; Abdominal pain 3.0, 440, 95, 160, 26, 522, 146, 555, 333, 999, 10000, 3040, 1865, 125, 95, 72; Abscesses 2720, 2170, 880, 787, 727, 190, 500; Acidosis 10000, 880, 802, 787, 776, 727, 20, 146; Acne 2720, 2170, 1800, 1600, 1500, 880, 787, 727; Actinomycosis 465, 10000, 787, 727, 20; Acupuncture disturbance field 18;

Acute Pain 160, 26, 3040, 95, 522, 146, 555, 333, 999, 10000, 880, 787, 727, 690, 666; Adenoids 1550, 802, 880, 787, 776, 727, 444, 20, 428, 660, 589, 333, 14, 2000, 770, 780, 805.6, 808.79, 2720, 2170; Adhesions 2720, 2170, 1550, 802, 880, 787, 776, 760, 727, 190; Adrenal stimulant 20; Allergy 10000, 3.0 and 330, 522, 146, 555, 333, 999, 125, 95, 72, 20, 444, 1865; Alopecia 10000, 880, 787, 727, 465, 146; Amenorrhea 444, 522, 146, 555, 333, 999, 125, 95, 72, 880, 10000, 3040, 787, 760, 727, 20; Anal itching 465, 444, 125, 95, 72, 444, 186, 20, 2000, 522, 146, 555, 333, 999, 10000, 880, 787, 760, 727; Angina 787, 776, 727, 465, 428, 660, 589, 333, 14, 2000, 780, 1570, 1770, 805.6; Ankylosing spondylitis 3040, 95, 3352, 880, 787, 776, 727, 1875, 28, 1.2, 10, 35, 28, 7.7, 1.2, 110, 100, 60, 808.79, 428, 660, 333, 14, 2000; Antiseptic effect 2352, 880, 787, 760, 727, 465, 440, 450, 428, 660, 589, 555, 333, 14, 2000, 1570, 1770; Apoplexy 40, 20, 1800, 1875, 125, 95, 72, 20, 2309, 522, 146, 555, 333, 999, 880, 787, 727; Appetite, lack of 465, 2309, 125, 95, 72, 20, 10000, 880, 787, 727; Arteriosclerosis 10000, 2720, 2170, 1800, 1600, 1500, 880, 787, 776, 727, 20, 808.79; Arthritis, 10000, 1550, 880, 802, 787, 727, 28, 20, 7.7, 3.0, 1.2; Arthritis, rheumatoid of the muscles and tendons 250, 1.2, 1875, 787, 727, 808.79; Arthritis, arthrosis and parathyroid disturbances affecting calcium metabolism, 9.6; Arthritis, arthralgia due to gout, 9.4; Arthritis focal origin 9.4; Asthma, 2720, 2170, 1800, 1600, 1500, 880, 787, 727, 808.79; Astrocytoma 9.2, 8.25, 7.7, 2170, 2127, 880, 690, 666; Ataxia 2720, 2170, 1800, 1600, 1500, 880, 787, 776, 727, 650, 625, 600, 2309, 125, 95, 72, 20, 805.6, 811.91, spastic ataxia, 9.2, 8.25, 7.7; Athlete's Foot 465, 1550, 880, 802, 727, 20, Autoimmune Disorders 10000, 3352, 880, 787, 776, 727, 650, 625, 600, 465, 250, 28, 20, 9.6, 9.4, 7.7, 3.0, 522, 146, 1850; Auto-intoxication 0.5, 522, 146, 10000, 880, 787, 727, 20, 555, 333, 999; Back pain 26, 186, 3040, 522, 146, 555, 333, 999, 10000, 760, 2352, 880, 787, 727; Bad Breath 2352, 880, 787, 727, 20; Bedsores 880, 2352, 787, 727, 485 followed by 1.2 and 73; Bilioussness 2352, 10000, 880, 832, 787, 727, 465; Bladder and prostate complaints 880, 2352, 787, 727, 465, 20, 9.4 Boils 880, 2352, 787, 727, 465, 660, 589, 333; Bone trauma 380, 2352, 10000, 880, 787, 727; Bone/periodontal disease 47.5, 1800, 1600, 1875, 880, 787, 776, 727; Bone protuberance 1.2 and 250; Brachial Neuralgia 0.5; Breast 880, 2352, 787, 776, 3083; Bright's Syndrome 1500, 880, 787, 727; Bronchial asthma 0.5, 522, 146, 125, 95, 72, 966, 20; Bronchitis 880, 727, 9.4, 9.35; Bronchial Pneumonia 2352, 880, 787, 776, 727; Bubonic Plague 500, plus secondary infections: 880, 787, 727, 20; Burns 190, 10000, 880, 787, 727, 465, 200; Bursitis 880, 787, 727; Cancer 465, 2127, 880, 787, 3083, 589; Cardiac Edema 9.2; Carpal Tunnel Syndrome 1550 thru 465, 6.3, 148, 2.5, 727, 522, 146, 444; Cataract 10000, 880, 787, 727; Catarrh 1550, 802, 880, 787, 727, 444, 20; Cerebral Palsy 880, 787, 727, 522, 146; Chicken Pox 2352, 1800, 1600, 1500, 880, 787, 727, 20; Cholera 880, 802, 450, 880, 787, 727; Chronic Fatigue Syndrome 10000, 2720, 1800, 880, 802, 787, 776, 727, 465, 966, 330, 148, 125, 95, 73, 72, 35, 27.5, 20, 6.3, 4.9, 428, 660; Circulatory stasis 40; Circulation disturbance problems 40.9.4; Concentration 7.83; Cold in head, chest 10,000, 1550, 880, 802, 787, 776, 727, 444, 20; Colic 1550, 832, 802, 787, 727, 20; Colitis 1550, 880, 832, 802, 440; Constipation 1550, 880, 802, 787, 776, 727, 20; Convulsions 10,000, 880, 787, 727; Contusions 9.1, 110; Costalgia 160, 26, 3040, 880, 787, 727, 10,000, 2353; Cramps 10,000, 880, 787, 727, 26; Cystitis (of urinary bladder) 1550, 880, 802, 787, 727, 465, 20; Deafness 10,000, 1550, 880, 802, 787, 727, 20; Dental foci 3040, 95, 190, 47.5, 1600, 1550, 2353,

1500, 880, 832, 787, 776, 760, 727, 1316, 1916, 465, 2720, 2489, 1800, 522, 146, 555, 333, 999; Depression, anxiety, trembling, weakness 3.5; Depression 35, 220, 1.1, 73; Detoxification 6.3, 148, 2.5, 727, 522, 146, 444, 666, 428, 880, 787, 555, 333, 999; Diabetes 837, 260, 1885, 125, 95, 1957, 522, 428, 49.4 for secondary diabetes infections 10,000, 2720, 2170, 1800, 1550, 880, 802, 727, 465, 20; Diabetic loading 35, 700; Diarrhea/Dysentery 1550, 880, 802, 787, 727, 465; Distorsion 9.1, 110; Dizziness 5.8; Dupuytren's Contracture 1.2, 250; Dysmenorrhea 26, 4.9, 2352, 880, 787, 727, 465; Dyspepsia 1550, 880, 802, 787, 727; Ears 9.2, 880, 787, 727, 20; Eczema 9.2, 1550, 802, 787, 727; Eczema 9.4; Enuresis 10,000, 880, 787, 727, 465; Edema 522, 146, 154.3, 444, 440, 880, 787, 727, 465; Epicondylitis 251.2, 3040, 160, 26; Epididymitis 1500, 880, 787, 7827, 20; Epilepsy 10,000, 880, 802, 787, 727, 700, 650, 600, 125, 20; Epstein-Barr 428, 660, 776, 465, 880, 787, 727; Erysipelas 660, 2000, 10,000, 880, 787, 727, 465, 20; Erythema nosodum 9.4; Esophagus 880, 787, 727; Eustachian tube inflammation 1550, 880, 802, 787, 776, 727, 465; Eye inflammation 1.2, 80; Eye Disorders 1600, 10,000, 880, 787, 727, 20; Eyes 1600, 880, 787, 727; Facial paralysis 10,000, 880, 787, 727; Fatigue 2720, 1800, 428, 660, 465, 125, 95, 72, 27.5, 20, 2309; Fever 880, 787, 727, 20; Fibroma 2127, 2008, 727, 690, 666, 1550, 802, 465; Fibrosis of lung 457.5; Fistulal Ulcer 880, 832, 787, 727; Flashes, hot (complications) 10,000, 880, 787, 727; Flatulence 1550, 880, 802, 787, 760, 727, 465; Flu 1550, 880, 802, 787, 727, 20; Foot (blisters) 10,000, 880, 787, 727, 600; Fractures 450, 10,000, 880, 787, 727; Frostbite 880, 787, 727; Frozen shoulder 10,000, 880, 802, 787, 727; Fungal infection 465, 1550, 880, 802, 727, 20; Furunkulosis herpes, skin diseases 1200, 1550, 802, 2000, 787, 727; Gall bladder dystonia with osteitis 2.6, 3040, 880, 787, 727, 20; Gallstones 2.65, 3040, 880, 787, 727, 20; Gangrene 880, 787, 727, 20, 73; Gas (intestinal) 1550, 880, 802, 787, 760, 727, 465; Gastritis and flatus 880, 832, 787, 727, 20; Gout 9.4, 3040, 880, 787, 727, 20; Gravel in urine 2.6, 3040, 880, 787, 727, 20; Gums (inflammation, gingivitis, pyorrhea) 1550, 880, 802, 787, 727, 465, 20; Hair loss 10,000, 880, 787, 727, 465, 146; Hallucinations 10,000, 880, 787, 727, 20; Hang-over 10,000, 146, 610; Hay Fever 880, 787, 727, 20; Headaches 160, 10, 6.3, 522, 146, 4.9, 3040, 727, 148, 125, 95, 72, 20, 10, 5.8, 2.5, 444, 685, 880, 787, 650, 625, 600, 3040, 428, 660, 555, 333, 999; For headaches caused by parasites 160, 125, 95, 73, 20, 727, 3040; Headaches (urogenitaliy caused) 160, 9.4, 3040, 555, 333, 999, 160, 9.6, 3040; Head Injuries 160, 9.6, 3040, 880, 787, 727, 522, 72, 5.8, 4.9; Hearing problems 10,000, 2158, 880, 787, 760, 727, 20; Heart (lab animals only) 80, 160, 20, 73, 3.9, 10,000, 880, 787, 727, 465, 4125 95, 20; Heartburn 832, 2720, 2170, 2127, 1800, 1600, 2352, 1500, 880, 787, 727, 685, 465, 2309, 125, 95, 72, 20; Hemorrhage 1550, 802; Hemorrhoids 1550, 880, 802, 727; Hepatitis (generic) 29.2, 1550, 880, 802, 727; Hernia 10,000, 787, 727; Herpes (zoster) 2720, 2170, 1800, 1600, 1500, 1550, 802, 1865, 880, 787, 727, 20; High blood pressure 10,000, 880, 787, 727, 9.2; Hip Pain (as in coxarthrosis) 3040, 880, 787, 727, 160, 26; Hives 1800, 880, 787, 727, 522, 146, 4.9; Hoarseness 880, 760, 727; Hyperacidity of stomach 7.83 on solar plexus and 230 on stomach; Hydrocele 880, 787, 727; Hyperosmia 20, 10,000, 522, 1865, 810.35; Hypothyroid 3.0, 0.5; Hypoacidity of stomach 15, 20; Hypothyroid 12, 35; Hypertension 20, 95, 9.2, 6.0, 10,000, 880, 787, 727; Hypertension, spastic 95; Impotence 9.4, 2127, 2008, 465, 880, 802, 787, 2083, 125, 95, 73, 72, 20, 1875; Indigestion 10,000, 880, 787, 727, 465, 444, 20, 125, 95, 72, 4.9;

Infantile Paralysis 1500, 880, 787, 727; Infections 1600, 1550, 1500, 880, 832, 802, 787, 776, 760, 727, 700, 690, 666, 650, 625, 600, 465, 444, 125, 95, 72, 20, 1865, 500, 450, 440, 428, 660; Infertility 2127, 2008, 465, 880, 802, 787, 2083, 1875, 9.4; Inflammation 1.5; Influenza 1550, 802, 1500, 880, 787, 727, 20; Insomnia 3.6, 3.0, 1550, 1500, 880, 802; Intercostal neuralgia 3040, 160, 26, 1550, 880, 802, 787, 776, 727, 125, 20, 1865, 444; Intermittent claudication 45, 48; Itching 880, 787, 727, 465, 444, 125, 95, 72, 7.83, 2309, 2000; Jaundice 1600, 1550, 1500, 880, 802, 650, 625, 600, 2309, 146, 250, 125, 95, 72, 20; Kidney insufficiency 50, 440, 1600, 1550, 1500, 880, 802, 787, 727, 650, 625, 600, 465, 73, 40, 10, 2309, 522, 146, 250, 148, 125, 95, 72, 20, 555, 333, 999; Knee/Joint pain 3040, 160, 26, 1550, 880, 802, 787, 727, 28, 20, 7.7, 3.0, 1.2, 250, 1.2, 9.6, 9.4; Larynx 440, 465, 444, 1550, 880, 802, 787, 727, 28, 7.7, 3.0, 1.2, 250, 1.2, 9.6, 9.4; Leprosy 600, 1600, 1550, 1500, 880, 832, 802, 787, 776, 760, 727, 700, 690, 666, 650, 625, 600, 465, 444, 20, 500, 450, 440, 428, 660; Leukemia 2127, 2008, 880, 787, 727, 690, 666, 589; Leukoderma 880, 787, 727, 444, 20, 589, 2000; Locomotor Dysfunction 10,000, 880, 787, 776, 727, 650, 625, 600, 2309, 125, 95, 72, 20; convulsions, spasticity 28, 8.25, 7.7; Low Blood Pressure 880, 787, 727, 20; Lumbago 10,000, 880, 787, 727, 125, 95, 72, 2309, 9.2, 8.25, 7.7; Lungs 1550, 802, 880, 787, 776, 727, 125, 95, 72, 20, 2309, 450, 589; Lupus Erythematosus 880, 787, 727, 776, 1850; Luxation 119.1; Lymph stasis 154.3, 2.5, 727, 522, 146, 444, 428, 880, 787; Malaria 555; Measles 333; Mental Disorders 444, 148, 20, 6.3, 2.5, 522, 146, 10,000, 125, 95, 72, 4.9, 727, 428, 444, 465 to 1550; Menieres Disease 2352, 880, 787, 727, 465, 428, 589, 333; Meningitis 2352, 880, 832, 787, 727, 650, 625, 600, 465, 2309, 125, 95, 72, 20, 428, 660; Menstrual problems 26, 3040, 880, 787, 727, 465, 160, 20; Migraine 170; Motion Sickness 650, 625, 600, 465, 444, 1865, 522, 146, 125, 95, 72, 20; Mouth eruptions 465, 2127, 2008, 727, 690, 666; Herpes sores 2489, 1800, 465, 1550, 1500, 880, 787, 727, 1850, 428; Mucous membrane inflammation 380; Multiple Sclerosis use auto immune frequencies plus 880, 787, 727; Muscles, to relax 760, 6.8; Muscular Dystrophy 10,000, 2352, 880, 787, 776, 727, 650, 625, 600, 465, 250, 28, 20, 9.6, 9.4, 7.7, 3.0, 522, 610, 1850, 146, 555, 333, 999; Muscular pain 3040, 26, 160, 320, 250, 240, 80, 40, 10, 5.8, 2.5, 1.5, 1.2, 1.0, 0.5; Muscle pain from parasites 2309, 125, 811, 108, 72, 20; Mumps 10,000, 727, 2720, 2489, 2127, 2008, 14, 428, 880, 787, 727, 20; Nausea 880, 832, 787, 727, 20, 4.9; Nephritis/nephrosis 880, 787, 727, 10, 40, 73, 465; Nervousness 3.0; Nerve disorders 10,000, 2720, 2489, 2170, 1800, 1600, 2352, 880, 787, 148, 6.3, 2.5, 522, 610, 727, 650, 625, 600, 125, 95, 72, 20, 428, 660, 333; Neuralgia 3.9, 555, 333, 999; Neurosis 28; Nose infection, congestion 2352, 880, 787, 776, 727, 444, 440, 465, 20; Numbness 10,000, 2720, 2489, 2170, 1800, 1600, 2352, 880, 787, 727, 650, 625, 600, 440, 660, 333; Oral Lesions 2720, 2489, 2008, 1800, 1600, 2352, 880, 787, 776, 760, 727, 465, 444, 522, 146, 555, 333, 999; Orchitis 2720, 2489, 2170, 2127, 2008, 1800, 1600, 2352, 1500, 880, 832, 787, 776, 2083, 650, 625, 600, 125, 95, 72, 20, 14; Osteomyelitis 2.6, 2720, 2489, 2170, 2127, 2008, 1800, 1600, 2352, 1500, 880, 832, 787, 776, 2083; Otosclerosis 9.2; Ovarian disorders 1875, 465, 444, 26, 2720, 2489, 2170, 2127, 2008, 1800, 1600, 2352, 1500, 880, 832, 787, 776, 2083, 20; Pain 3040, 95, 160, 26; Pain of infection 3040, 160, 26, 95, 880, 2352, 787, 776, 727, 4.9; Pain of cancer 3040, 160, 26, 95, 2127, 2008, 2083; Pancreatic insufficiency 15, 10, 250, 1875, 465, 444, 26, 2720, 2489, 2170, 2127, 2008, 1800, 1600, 2352, 1500, 880, 832, 787, 776, 2083, 20; Paralysis, spastic 10,000, 880,

787, 776, 727, 650, 625, 600, 2309, 125, 95, 72, 20, 7.7; Paralysis, nonspastic 10,000, 880, 787, 776, 727, 650, 625, 600, 444, 1865, 125, 95, 72, 20, 9.2, 8.25; Parasites 125, 95, 72, 20, 2309; Pelvic Inflammatory Disease 2720, 2489, 2170, 2127, 2008, 1800, 1600, 2352, 787, 776, 2083, 650, 625, 600, 465, 966, 95, 72, 450, 428, 14, 589, 1570; Pericarditis 2720, 2170, 1600, 880, 2352, 787, 760, 727, 625, 465, 160, 125, 73, 26, 20, 3.9, 95, 80, 73, 72, 2309; Pharyngitis 440, 589, 380, 1600, 2352, 776, 787, 880, 727, 20, 522, 146; Pleurisy 1550, 802, 880, 787, 776, 727, 125, 95, 72, 20, 2309, 450, 589; Pneumonia 1550, 802, 880, 787, 776, 727, 20, 450, 589; Poliomyelitis 805.6 secondary complications 2158, 428, 1500, 880, 787, 727; Polyps 2720, 2489, 2170, 2127, 2008, 1800, 1600, 2083, 1875, 465, 444, 20, 522, 146; Pre-op and post-operation prevention and control of nosocomial and idiopathic infection 2170, 1800, 1600, 2352, 1500, 880, 832, 802, 787, 776, 727, 465, 444, 522, 146, 428, 589, 333; Prostatitis 510, 522, 146, 2720, 2489, 2170, 2127, 2008, 2352, 787, 776, 2083, 465, 125, 95, 72, 20, 1249; Prostate complaints 9.4, 2127, 2008, 727, 690, 666, 465, 880, 802, 787, 727, 125, 95, 73, 72, 20; Protritis 880, 787, 727, 465, 444, 125, 95, 72, 20, 2309, 2000; Prostate tumor (malignant) 2127, 2008, 2083; Psychosomatic pain 727, 160, 26, 148, 455, 511, 146, 555, 333, 999; Pyorrhea 2720, 2489, 2008, 1800, 1600, 2352, 880, 787, 776, 727, 465, 444, 522, 146; Rabies 807.23; Raynaud's Disease 727, 20; Rheumatoid arthritis of the muscles and tendons 251.2, 9.6, 9.4, 1875, 10,000, 2352, 880, 787, 727, 465, 20, 7.7, 3.0, 1.2; Rhinitis 10,000, 2352, 1500, 880, 787, 776, 727, 444, 7.83, 522, 146, 465; Sarcoma 2127, 2008, 880, 787, 727, 690, 666, 589; Scars 465 to 1568; Scarlet Fever 880, 787, 727, 1356; Sciatica or Ischia 2352, 880, 787, 727, 690, 666, 10; Sedative Effect 2.5; Sexual Dysfunction 9.4, 2127 2008, 465, 880, 802, 787, 2083, 125, 95, 73, 72, 20, 1875; Shingles 2352, 802, 1800, 1865, 1880, 787, 727, 26; Sinusitis 2352, 880, 787, 727, 125, 72, 522; Sleeping Sickness 124; Slipped discs 26, 125, 880, 787, 333, 787, 727, 95, 72, 522, 146, 555, 333, 999; Smallpox 811.91, 1550, 802, 880, 787, 727; Sneezing 880, 787, 727, 465, 146; Sore Throat 2720, 2489, 1800, 1600, 2352, 880, 787, 776, 727, 465, 589, 808.79; Spasms, muscle 6.8; Spastic paresis 650, 625, 30.87.48; Spleen, enlarged 35, 787, 10,000, 2720, 2170, 1800, 1550, 880, 802, 727, 465, 20; Spondylitis, acute 1.2, 10, 7.7, 28; Staphylococci infection 727; Stiff neck 4.9; Spastic stiff neck 9.2, 6.0; Stiff muscles in general 320, 250, 240, 160, 125, 80, 40, 20, 10, 5.8, 2.5, 1.5, 1.2, 1.0, 0.5, 1800, 2352, 880, 787, 776, 727; Stomach disorders 2127, 2008, 880, 787, 727, 2083, 125, 95, 72, 20, 3.9, 450, 1570, 1770; Stones 3.5; Streptococci infection 880; Streptothrix infection 787; Sty 10,000, 880, 787, 727, 20; Subluxation induced disorders 9.6; Sun Allergy 3.0, 330; Sunstroke 444, 428, 190, 3040, 95, 522, 146, 880, 20; Surgery—Prevention and control of nosocomial and idiopathic infection 2170, 1800 1600, 2352, 1500, 880, 832, 802, 787, 776, 727, 465, 444, 522, 146, 428, 589, 333, 3040, 95, 160, 26, 544; Swelling 522, 146, 6.3, 148, 444, 428, 880, 787, 727; Swollen glands 880, 787, 727, 20; Syphilis 1875; Tachycardia (also used for pain of arthritis, headache, Tui-na, and facial toning) 1.2; Teeth 3040, 95, 190, 47.5, 1600, 2352, 1500, 880, 10,000, 2489, 787, 776, 760, 727, 1875, 465, 2720, 2489, 1800; Tendo myopathy 3040, 320, 250, 160, 80, 40, 26, 20, 10, 5.8, 2.5, 1.5, 1.2, 1.0, 0.5, 522, 146; Tetanus 811.91, 600, 880, 787, 727; Thrush 465; Tonsillitis 74.2, 73, 2352, 1500, 880, 832, 787, 776, 727, 1875, 465, 808.79, 589; Toothache 3040, 95, 190, 47.5, 2720, 2489, 1800, 1600, 2352, 1500, 880, 832, 787, 776, 727, 666, 650, 600, 465, 522, 146, 555, 333, 999; Tooth extraction, afterward 3040,

95, 47.5, 7.83; Toxins 444, 160, 522, 146, 555, 333, 999, 10,000, 880, 787, 727, 125, 95, 72, 20; Trauma 3040, 95, 190, 760, 880, 787, 727, 465; Trench Mouth 1875, 880, 787, 776, 727, 650, 625, 465; Trigeminal neuralgia 2720, 2489, 2170, 1800, 1600, 2352, 880, 832, 787, 776, 760, 727, 650, 146, 7.83, 27.5, 428, 660, 589, 333, 999, 770, 2000, 805.6, 807.23, 811.91, 810.35; Tuberculosis 2352, 1875, 776, 2127, 2008, 465, 589; Typhoid 1570, 1770, 2352, 690, 1800; Ulcers 2489, 2170, 2127, 1800, 1600, 880, 832, 802, 787, 776, 727, 1.2 through 5.8, 1.2, 73; Urethritis 2720, 2170, 2127, 1800, 1600, 2352, 1500, 880, 832, 787, 776, 727, 1875, 455, 2309, 125, 95, 72, 1.2; Urticaria 1800, 444, 125, 95, 72, 522, 146, 787, 727, 6.3; Varicose 1.2, 28; Vegetative Dystonia 40; Vein thrombosis 685; Vertigo 60, 5.8; Warts 3040, 5, 2170, 2127, 2008, 1800, 1600, 1500, 2083; Worms 125, 95, 6.0, 20, 2309; Wound healing 880, 787, 727, 220, 190, 20, 40; Yeast infection 465; Yellow Fever 999; Scan brain for troubled areas 7.83; For relaxation 4.9; To stimulate mental clarity 12; Facial toning 9.6; Kidney meridian 9.2.

6. The system as defined in claim 2 wherein said frequency of said generated signal comprises a selected frequency in cycles per second from the group of specific frequencies for inactivating micro-organisms and viruses which may be present and enhancing specific tissue regeneration in a mammal, said group of frequencies and the purpose thereof consisting of Adynamia, geriatric 60, 27.5; Abdominal inflammation 380, 1.2, 2727, 2489, 2170, 1800, 1600, 2302, 880, 832, 787, 776, 727, 465, 444, 1865, 146, 125, 95, 72, 20, 450, 440, 428, 660, 589, 2000, 1570, 1770, 805.6, 811.91; Abdominal pain 3.0, 440, 95, 160, 26, 522, 146, 555, 333, 999, 10000, 3040, 1865, 125, 95, 72; Abscesses 2720, 2170, 880, 787, 727, 190, 500; Acidosis 10000, 880, 802, 787, 776, 727, 20, 146; Acne 2720, 2170, 1800, 1600, 1500, 880, 787, 727; Actinomycosis 465, 10000, 787, 727, 20; Acupuncture disturbance field 18; Acute Pain 160, 26, 3040, 95, 522, 146, 555, 333, 999, 10000, 880, 787, 727, 690, 666; Adenoids 1550, 802, 880, 787, 776, 727, 444, 20, 428, 660, 589, 333, 14, 2000, 770, 780, 805.6, 808.79, 2720, 2170; Adhesions 2720, 2170, 1550, 802, 880, 787, 776, 760, 727, 190; Adrenal stimulant 20; Allergy 10000, 3.0 and 330, 522, 146, 555, 333, 999, 125, 95, 72, 20, 444, 1865; Alopecia 10000, 880, 787, 727, 465, 146; Amenorrhea 444, 522, 146, 555, 333, 999, 125, 95, 72, 880, 10000, 3040, 787, 760, 727, 20; Anal itching 465, 444, 125, 95, 72, 444, 186, 20, 2000, 522, 146, 555, 333, 999, 10000, 880, 787, 760, 727; Angina 787, 776, 727, 465, 428, 660, 589, 333, 14, 2000, 780, 1570, 1770, 805.6; Ankylosing spondylitis 3040, 95, 3352, 880, 787, 776, 727, 1875, 28, 1.2, 10, 35, 28, 7.7, 1.2, 110, 100, 60, 808.79, 428, 660, 333, 14, 2000; Antiseptic effect 2352, 880, 787, 760, 727, 465, 444, 450, 428, 660, 589, 555, 333, 14, 2000, 1570, 1770; Apoplexy 40, 20, 1800, 1875, 125, 95, 72, 20, 2309, 522, 146, 555, 333, 999, 880, 787, 727; Appetite, lack of 465, 2309, 125, 95, 72, 20, 10000, 880, 787, 727; Arteriosclerosis 10000, 2720, 2170, 1800, 1600, 1500, 880, 787, 776, 727, 20, 808.79; Arthritis, 10000, 1550, 880, 802, 787, 727, 28, 20, 7.7, 3.0, 1.2; Arthritis, rheumatoid of the muscles and tendons 250, 1.2, 1875, 787, 727, 808.79; Arthritis, arthrosis and parathyroid disturbances affecting calcium metabolism, 9.6; Arthritis, arthralgia due to gout, 9.4; Arthritis focal origin 9.4; Asthma, 2720, 2170, 1800, 1600, 1500, 880, 787, 727, 808.79; Astrocytoma 9.2, 8.25, 7.7, 2170, 2127, 880, 690, 666; Ataxia 2720, 2170, 1800, 1600, 1500, 880, 787, 776, 727, 650, 625, 600, 2309, 125, 95, 72, 20, 805.6, 811.91, spastic ataxia, 9.2, 8.25, 7.7; Athlete's Foot 465, 1550, 880, 802, 727, 20, Autoimmune

Disorders 10000, 3352, 880, 787, 776, 727, 650, 625, 600, 465, 250, 28, 20, 9.6, 9.4, 7.7, 3.0, 522, 146, 1850; Auto-intoxication 0.5, 522, 146, 10000, 880, 787, 727, 20, 555, 333, 999; Back pain 26, 186, 3040, 522, 146, 555, 333, 999, 10000, 760, 2352, 880, 787, 727; Bad Breath 2352, 880, 787, 727, 20; Bedsores 880, 2352, 787, 727, 485 followed by 1.2 and 73; Biliousness 2352, 10000, 880, 832, 787, 727, 465; Bladder and prostate complaints 880, 2352, 787, 727, 465, 20, 9.4 Boils 880, 2352, 787, 727, 465, 660, 589, 333; Bone trauma 380, 2352, 10000, 880, 787, 727; Bone/periodontal disease 47.5, 1800, 1600, 1875, 880, 787, 776, 727; Bone protuberance 1.2 and 250; Brachial Neuralgia 0.5; Breast 880, 2352, 787, 776, 3083; Bright's Syndrome 1500, 880, 787, 727; Bronchial asthma 0.5, 522, 146, 125, 95, 72, 966, 20; Bronchitis 880, 727, 9.4, 9.35; Bronchial Pneumonia 2352, 880, 787, 776, 727; Bubonic Plague 500, plus secondary infections: 880, 787, 727, 20; Burns 190, 10000, 880, 787, 727, 465, 200; Bursitis 880, 787, 727; Cancer 465, 2127, 880, 787, 3083, 589; Cardiac Edema 9.2; Carpal Tunnel Syndrome 1550 thru 465, 6.3, 148, 2.5, 727, 522, 146, 444; Cataract 10000, 880, 787, 727; Catarrh 1550, 802, 880, 787, 727, 444, 20; Cerebral Palsy 880, 787, 727, 522, 146; Chicken Pox 2352, 1800, 1600, 1500, 880, 787, 727, 20; Cholera 880, 802, 450, 880, 787, 727; Chronic Fatigue Syndrome 10000, 2720, 1800, 880, 802, 787, 776, 727, 465, 966, 330, 148, 125, 95, 73, 72, 35, 27.5, 20, 6.3, 4.9, 428, 660; Circulatory stasis 40; Circulation disturbance problems 40, 9.4; Concentration 7.83; Cold in head, chest 10,000, 1550, 880, 802, 787, 776, 727, 444, 20; Colic 1550, 832, 802, 787, 727, 20; Colitis 1550, 880, 832, 802, 440; Constipation 1550, 880, 802, 787, 776, 727, 20; Convulsions 10,000, 880, 787, 727; Contusions 9.1, 110; Costalgia 160, 26, 3040, 880, 787, 727, 10,000, 2353; Cramps 10,000, 880, 787, 727, 26; Cystitis (of urinary bladder) 1550, 880, 802, 787, 727, 465, 20; Deafness 10,000, 1550, 880, 802, 787, 727, 20; Dental foci 3040, 95, 190, 47.5, 1600, 1550, 2353, 1500, 880, 832, 787, 776, 760, 727, 1316, 1916, 465, 2720, 2489, 1800, 522, 146, 555, 333, 999; Depression, anxiety, trembling, weakness 3.5; Depression 35, 220, 1.1, 73; Detoxification 6.3, 148, 2.5, 727, 522, 146, 444, 666, 428, 880, 787, 555, 333, 999; Diabetes 837, 260, 1885, 125, 95, 1957, 522, 428, 49.4 for secondary diabetes infections 10,000, 2720, 2170, 1800, 1550, 880, 802, 727, 465, 20; Diabetic loading 35, 700; Diarrhea/Dysentery 1550, 880, 802, 787, 727, 465; Distorsion 9.1, 110; Dizziness 5.8; Dupuytren's Contracture 1.2, 250; Dysmenorrhea 26, 4.9, 2352, 880, 787, 727, 465; Dyspepsia 1550, 880, 802, 787, 727; Ears 9.2, 880, 787, 727, 20; Eczema 9.2, 1550, 802, 787, 727; Eczema 9.4; Enuresis 10,000, 880, 787, 727, 465; Edema 522, 146, 154.3, 444, 440, 880, 787, 727, 465; Epicondylitis 251.2, 3040, 160, 26; Epididymitis 1500, 880, 787, 7827, 20; Epilepsy 10,000, 880, 802, 787, 727, 700, 650, 600, 125, 20; Epstein-Barr 428, 660, 776, 465, 880, 787, 727; Erysipelas 660, 2000, 10,000, 880, 787, 727, 465, 20; Erythema nosodum 9.4; Esophagus 880, 787, 727; Eustachian tube inflammation 1550, 880, 802, 787, 776, 727, 465; Eye inflammation 1.2, 80; Eye Disorders 1600, 10,000, 880, 787, 727, 20; Eyes 1600, 880, 787, 727; Facial paralysis 10,000, 880, 787, 727; Fatigue 2720, 1800, 428, 660, 465, 125, 95, 72, 27.5, 20, 2309; Fever 880, 787, 727, 20; Fibroma 2127, 2008, 727, 690, 666, 1550, 802, 465; Fibrosis of lung 457.5; Fistular Ulcer 880, 832, 787, 727; Flashes, hot (complications) 10,000, 880, 787, 727; Flatulence 1550, 880, 802, 787, 760, 727, 465; Flu 1550, 880, 802, 787, 727, 20; Foot (blisters) 10,000, 880, 787, 727, 600; Fractures 450, 10,000, 880, 787, 727; Frostbite 880, 787, 727; Frozen shoulder 10,000, 880, 802, 787, 727; Fungal

infection 465, 1550, 880, 802, 727, 20; Furunkulosis herpes, skin diseases 1200, 1550, 802, 2000, 787, 727; Gall bladder dystonia with osteitis 2.6, 3040, 880, 787, 727, 20; Gallstones 2.65, 3040, 880, 787, 727, 20; Gangrene 880, 787, 727, 20, 73; Gas (intestinal) 1550, 880, 802, 787, 760, 727, 465; Gastritis and flatus 880, 832, 787, 727, 20; Gout 9.4, 3040, 880, 787, 727, 20; Gravel in urine 2.6, 3040, 880, 787, 727, 20; Gums (inflammation, gingivitis, pyorrhea) 1550, 880, 802, 787, 727, 465, 20; Hair loss 10,000, 880, 787, 727, 465, 146; Hallucinations 10,000, 880, 787, 727, 20; Hang-over 10,000, 146, 610; Hay Fever 880, 787, 727, 20; Headaches 160, 10, 6.3, 522, 146, 4.9, 3040, 727, 148, 125, 95, 72, 20, 10, 5.8, 2.5, 444, 685, 880, 787, 650, 625, 600, 3040, 428, 660, 555, 333, 999; For headaches caused by parasites 160, 125, 95, 73, 20, 727, 3040; Headaches (urogenitali caused) 160, 9.4, 3040, 555, 333, 999, 160, 9.6, 3040; Head Injuries 160, 9.6, 3040, 880, 787, 727, 522, 72, 5.8, 4.9; Hearing problems 10,000, 2158, 880, 787, 760, 727, 20; Heart (lab animals only) 80, 160, 20, 73, 3.9, 10,000, 880, 787, 727, 465, 4125 95, 20; Heartburn 832, 2720, 2170, 2127, 1800, 1600, 2352, 1500, 880, 787, 727, 685, 465, 2309, 125, 95, 72, 20; Hemorrhage 1550, 802; Hemorrhoids 1550, 880, 802, 727; Hepatitis (generic) 29.2, 1550, 880, 802, 727; Hernia 10,000, 787, 727; Herpes (zoster) 2720, 2170, 1800, 1600, 1500, 1550, 802, 1865, 880, 787, 727, 20; High blood pressure 10,000, 880, 787, 727, 9.2; Hip Pain (as in coxarthrosis) 3040, 880, 787, 727, 160, 26; Hives 1800, 880, 787, 727, 522, 146, 4.9; Hoarseness 880, 760, 727; Hyperacidity of stomach 7.83 on solar plexus and 230 on stomach; Hydrocele 880, 787, 727; Hyperosmia 20, 10,000, 522, 1865, 810.35; Hyperthyroid 3.0, 0.5; Hypoacidity of stomach 15, 20; Hypothyroid 12, 35; Hypertension 20, 95, 9.2, 6.0, 10,000, 880, 787, 727; Hypertension, spastic 95; Impotence 9.4, 2127, 2008, 465, 880, 802, 787, 2083, 125, 95, 73, 72, 20, 1875; Indigestion 10,000, 880, 787, 727, 465, 444, 20, 125, 95, 72, 4.9; Infantile Paralysis 1500, 880, 787, 727; Infections 1600, 1550, 1500, 880, 832, 802, 787, 776, 760, 727, 700, 690, 666, 650, 625, 600, 465, 444, 125, 95, 72, 20, 1865, 500, 450, 440, 428, 660; Infertility 2127, 2008, 465, 880, 802, 787, 2083, 1875, 9.4; Inflammation 1.5; Influenza 1550, 802, 1500, 880, 787, 727, 20; Insomnia 3.6, 3.0, 1550, 1500, 880, 802; Intercostal neuralgia 3040, 160, 26, 1550, 880, 802, 787, 776, 727, 125, 20, 1865, 444; Intermittent claudication 45, 48; Itching 880, 787, 727, 465, 444, 125, 95, 72, 7.83, 2309, 2000; Jaundice 1600, 1550, 1500, 880, 802, 650, 625, 600, 2309, 146, 250, 125, 95, 72, 20; Kidney insufficiency 50, 440, 1600, 1550, 1500, 880, 802, 787, 727, 650, 625, 600, 465, 73, 40, 10, 2309, 522, 146, 250, 148, 125, 95, 72, 20, 555, 333, 999; Knee/Joint pain 3040, 160, 26, 1550, 880, 802, 787, 727, 28, 20, 7.7, 3.0, 1.2, 250, 1.2, 9.6, 9.4; Larynx 440, 465, 444, 1550, 880, 802, 787, 727, 28, 7.7, 3.0, 1.2, 250, 1.2, 9.6, 9.4; Leprosy 600, 1600, 1550, 1500, 880, 832, 802, 787, 776, 760, 727, 700, 690, 666, 650, 625, 600, 465, 444, 20, 500, 450, 440, 428, 660; Leukemia 2127, 2008, 880, 787, 727, 690, 666, 589; Leukoderma 880, 787, 727, 444, 20, 589, 2000; Locomotor Dysfunction 10,000, 880, 787, 776, 727, 650, 625, 600, 2309, 125, 95, 72, 20; convulsions, spasticity 28, 8.25, 7.7; Low Blood Pressure 880, 787, 727, 20; Lumbago 10,000, 880, 787, 727, 125, 95, 72, 2309, 9.2, 8.25, 7.7; Lungs 1550, 802, 880, 787, 776, 727, 125, 95, 72, 20, 2309, 450, 589; Lupus Erythematosus 880, 787, 727, 776, 1850; Luxation 119.1; Lymph stasis 154.3, 2.5, 727, 522, 146, 444, 428, 880, 787; Malaria 555; Measles 333; Mental Disorders 444, 148, 20, 6.3, 2.5, 522, 146, 10,000, 125, 95, 72, 4.9, 727, 428, 444, 465 to 1550; Menieres Disease 2352, 880, 787, 727, 465, 428, 589, 333;

Meningitis 2352, 880, 832, 787, 727, 650, 625, 600, 465, 2309, 125, 95, 72, 20, 428, 660; Menstrual problems 26, 3040, 880, 787, 727, 465, 160, 20; Migraine 170; Motion Sickness 650, 625, 600, 465, 444, 1865, 522, 146, 125, 95, 72, 20; Mouth eruptions 465, 2127, 2008, 727, 690, 666; Herpes sores 2489, 1800, 465, 1550, 1500, 880, 787, 727, 1850, 428; Mucous membrane inflammation 380; Multiple Sclerosis use auto immune frequencies plus 880, 787, 727; Muscles, to relax 760, 6.8; Muscular Dystrophy 10,000, 2352, 880, 787, 776, 727, 650, 625, 600, 465, 250, 28, 20, 9.6, 9.4, 7.7, 3.0, 522, 610, 1850, 146, 555, 333, 999; Muscular pain 3040, 26, 160, 320, 250, 240, 80, 40, 10, 5.8, 2.5, 1.5, 1.2, 1.0, 0.5; Muscle pain from parasites 2309, 125, 811, 108, 72, 20; Mumps 10,000, 727, 2720, 2489, 2127, 2008, 14, 428, 880, 787, 727, 20; Nausea 880, 832, 787, 727, 20, 4.9; Nephritis/nephrosis 880, 787, 727, 10, 40, 73, 465; Nervousness 3.0; Nerve disorders 10,000, 2720, 2489, 2170, 1800, 1600, 2352, 880, 787, 148, 6.3, 2.5, 522, 610, 727, 650, 625, 600, 125, 95, 72, 20, 428, 660, 333; Neuralgia 3.9, 20 555, 333, 999; Neurosis 28; Nose infection, congestion 2352, 880, 787, 776, 727, 444, 440, 465, 20; Numbness 10,000, 2720, 2489, 2170, 1800, 1600, 2352, 880, 787, 727, 650, 625, 600, 440, 660, 333; Oral Lesions 2720, 2489, 2008, 1800, 1600, 2352, 880, 787, 776, 760, 727, 465, 444, 522, 146, 555, 333, 999; Orchitis 2720, 2489, 2170, 2127, 2008, 1800, 1600, 2352, 1500, 880, 832, 787, 776, 2083, 650, 625, 600, 125, 95, 72, 20, 14; Osteomyelitis 2.6, 2720, 2489, 2170, 2127, 2008, 1800, 1600, 2352, 1500, 880, 832, 787, 776, 2083; Otosclerosis 9.2; Ovarian disorders 1875, 30 465, 444, 26, 2720, 2489, 2170, 2127, 2008, 1800, 1600, 2352, 1500, 880, 832, 787, 776, 2083, 20; Pain 3040, 95, 160, 26; Pain of infection 3040, 160, 26, 95, 880, 2352, 787, 776, 727, 4.9; Pain of cancer 3040, 160, 26, 95, 2127, 2008, 2083; Pancreatic insufficiency 15, 10, 250, 1875, 465, 444, 26, 2720, 2489, 2170, 2127, 2008, 1800, 1600, 2352, 1500, 880, 832, 787, 776, 2083, 20; Paralysis, spastic 10,000, 880, 787, 776, 727, 650, 625, 600, 2309, 125, 95, 72, 20, 7.7; Paralysis, nonspastic 10,000, 880, 787, 776, 727, 650, 625, 600, 444, 1865, 125, 95, 72, 20, 9.2, 8.25; Parasites 125, 95, 72, 20, 2309; Pelvic Inflammatory Disease 2720, 2489, 2170, 2127, 2008, 1800, 1600, 2352, 787, 776, 2083, 650, 625, 600, 465, 966, 95, 72, 450, 428, 14, 589, 1570; Pericarditis 2720, 2170, 1600, 880, 2352, 787, 760, 727, 625, 465, 160, 125, 73, 26, 20, 3.9, 95, 80, 73, 72, 2309; Pharyngitis 440, 589, 380, 1600, 2352, 776, 787, 880, 727, 20, 522, 146; Pleurisy 1550, 802, 880, 787, 776, 727, 125, 95, 72, 20, 2309, 450, 589; Pneumonia 1550, 802, 880, 787, 776, 727, 20, 450, 589; Poliomyelitis 805.6 secondary complications 2158, 428, 1500, 880, 787, 727; Polyps 2720, 50 2489, 2170, 2127, 2008, 1800, 1600, 2083, 1875, 465, 444, 20, 522, 146; Pre-op and post-operation prevention and control of nosocomial and idiopathic infection 2170, 1800, 1600, 2352, 1500, 880, 832, 802, 787, 776, 727, 465, 444, 522, 146, 428, 589, 333; Prostatitis 510, 522, 146, 2720, 2489, 2170, 2127, 2008, 2352, 787, 776, 2083, 465, 125, 95, 72, 20, 1249; Prostate complaints 9.4, 2127, 2008, 727, 690, 666, 465, 880, 802, 787, 727, 125, 95, 73, 72, 20; Protritis 880, 787, 727, 465, 444, 125, 95, 72, 20, 2309, 2000; Prostate tumor (malignant) 2127, 2008, 2083; Psychosomatic pain 727, 160, 26, 148, 455, 511, 146, 555, 333, 999; Pyorrhea 2720, 2489, 2008, 1800, 1600, 2352, 880, 787, 776, 727, 465, 444, 522, 146; Rabies 807.23; Raynaud's Disease 727, 20; Rheumatoid arthritis of the muscles and tendons 251.2, 9.6, 9.4, 1875, 10,000, 2352, 880, 787, 727, 465, 20, 7.7, 3.0, 1.2; Rhinitis 10,000, 2352, 1500, 880, 787, 776, 727, 444, 7.83, 522, 146, 465; Sarcoma 2127, 2008, 880, 787, 727, 690, 666, 589; Scars 465 to 1568; Scarlet

Fever 880, 787, 727, 1356; Sciatica or Ischia 2352, 880, 787, 727, 690, 666, 10; Sedative Effect 2.5; Sexual Dysfunction 9.4, 2127 2008, 465, 880, 802, 787, 2083, 125, 95, 73, 72, 20, 1875; Shingles 2352, 802, 1800, 1865, 1880, 787, 727, 26; Sinusitis 2352, 880, 787, 727, 125, 72, 522; Sleeping Sickness 124; Slipped discs 26, 125, 880, 787, 333, 787, 727, 95, 72, 522, 146, 555, 333, 999; Smallpox 811.91, 1550, 802, 880, 787, 727; Sneezing 880, 787, 727, 465, 146; Sore Throat 2720, 2489, 1800, 1600, 2352, 880, 787, 776, 727, 465, 589, 808.79; Spasms, muscle 6.8; Spastic paresis 650, 625, 30.87.48; Spleen, enlarged 35, 787, 10,000, 2720, 2170, 1800, 1550, 880, 802, 727, 465, 20; Spondylitis, acute 1.2, 10, 7.7, 28; Staphylococci infection 727; Stiff neck 4.9; Spastic stiff neck 9.2, 6.0; Stiff muscles in general 320, 250, 240, 160, 125, 80, 40, 20, 10, 5.8, 2.5, 1.5, 1.2, 1.0, 0.5, 1800, 2352, 880, 787, 776, 727; Stomach disorders 2127, 2008, 880, 787, 727, 2083, 125, 95, 72, 20, 3.9, 450, 1570, 1770; Stones 3.5; Streptococci infection 880; Streptothrix infection 787; Sty 10,000, 880, 787, 727, 20; Subluxation induced disorders 9.6; Sun Allergy 3.0, 330; Sunstroke 444, 428, 190, 3040, 95, 522, 146, 880, 20; Surgery—Prevention and control of nosocomial and idiopathic infection 2170, 1800 1600, 2352, 1500, 880, 832, 802, 787, 776, 727, 465, 444, 522, 146, 428, 589, 333, 3040, 95, 160, 26, 544; Swelling 522, 146, 6.3, 148, 444, 428, 880, 787, 727; Swollen glands 880, 787, 727, 20; Syphilis 1875; Tachycardia (also used for pain of arthritis, headache, Tui-na, and facial toning) 1.2; Teeth 3040, 95, 190, 47.5, 1600, 2352, 1500, 880, 10,000, 2489, 787, 776, 760, 727, 1875, 465, 2720, 2489, 1800; Tendo myopathy 3040, 320, 250, 160, 80, 40, 26, 20, 10, 5.8, 2.5, 1.5, 1.2, 1.0, 0.5, 522, 146; Tetanus 811.91, 600, 880, 787, 727; Thrush 465; Tonsillitis 74.2, 73, 2352, 1500, 880, 832, 787, 776, 727, 1875, 465, 808.79, 589; Toothache 3040, 95, 190, 47.5, 2720, 2489, 1800, 1600, 2352, 1500, 880, 832, 787, 776, 727, 666, 650, 600, 465, 522, 146, 555, 333, 999; Tooth extraction, afterward 3040, 95, 47.5, 7.83; Toxins 444, 160, 522, 146, 555, 333, 999, 10,000, 880, 787, 727, 125, 95, 72, 20; Trauma 3040, 95, 190, 760, 880, 787, 727, 465; Trench Mouth 1875, 880, 787, 776, 727, 650, 625, 465; Trigeminal neuralgia 2720, 2489, 2170, 1800, 1600, 2352, 880, 832, 787, 776, 760, 727, 650, 146, 7.83, 27.5, 428, 660, 589, 333, 999, 770, 2000, 805.6, 807.23, 811.91, 810.35; Tuberculosis 2352, 1875, 776, 2127, 2008, 465, 589; Typhoid 1570, 1770, 2352, 690, 1800; Ulcers 2489, 2170, 2127, 1800, 1600, 880, 832, 802, 787, 776, 727, 1.2 through 5.8, 1.2, 73; Urethritis 2720, 2170, 2127, 1800, 1600, 2352, 1500, 880, 832, 787, 776, 727, 1875, 455, 2309, 125, 95, 72, 1.2; Urticaria 1800, 444, 125, 95, 72, 522, 146, 787, 727, 6.3; Varicoses 1.2, 28; Vegetative Dystonia 40; Vein thrombosis 685; Vertigo 60, 5.8; Warts 3040, 5, 2170, 2127, 2008, 1800, 1600, 1500, 2083; Worms 125, 95, 6.0, 20, 2309; Wound healing 880, 787, 727, 220, 190, 20, 40; Yeast infection 465; Yellow Fever 999; Scan brain for troubled areas 7.83; For relaxation 4.9; To stimulate mental clarity 12; Facial toning 9.6; Kidney meridian 9.2.

7. The system as defined in claim 3 wherein said frequency of said generated signal comprises a selected frequency in cycles per second from the group of specific frequencies for inactivating micro-organisms and viruses which may be present and enhancing specific tissue regeneration in a mammal, said group of frequencies and the purpose therefor consisting of Adynamia, geriatric 60, 27.5; Abdominal inflammation 380, 1.2, 2727, 2489, 2170, 1800, 1600, 2302, 880, 832, 787, 776, 727, 465, 444, 1865, 146, 125, 95, 72, 20, 450, 440, 428, 660, 589, 2000, 1570, 1770, 805.6, 811.91; Abdominal pain 3.0, 440, 95, 160, 26, 522,

146, 555, 333, 999, 10000, 3040, 1865, 125, 95, 72; Abscesses 2720, 2170, 880, 787, 727, 190, 500; Acidosis 10000, 880, 802, 787, 776, 727, 20, 146; Acne 2720, 2170, 1800, 1600, 1500, 880, 787, 727; Actinomycosis 465, 10000, 787, 727, 20; Acupuncture disturbance field 18; Acute Pain 160, 26, 3040, 95, 522, 146, 555, 333, 999, 10000, 880, 787, 727, 690, 666; Adenoids 1550, 802, 880, 787, 776, 727, 444, 20, 428, 660, 589, 333, 14, 2000, 770, 780, 805.6, 808.79, 2720, 2170; Adhesions 2720, 2170, 1550, 802, 880, 787, 776, 760, 727, 190; Adrenal stimulant 20; Allergy 10000, 3.0 and 330, 522, 146, 555, 333, 999, 125, 95, 72, 20, 444, 1865; Alopecia 10000, 880, 787, 727, 465, 146; Amenorrhea 444, 522, 146, 555, 333, 999, 125, 95, 72, 880, 10000, 3040, 787, 760, 727, 20; Anal itching 465, 444, 125, 95, 72, 444, 186, 20, 2000, 522, 146, 555, 333, 999, 10000, 880, 787, 760, 727; Angina 787, 776, 727, 465, 428, 660, 589, 333, 14, 2000, 780, 1570, 1770, 805.6; Ankylosing spondylitis 3040, 95, 3352, 880, 787, 776, 727, 1875, 28, 1.2, 10, 35, 28, 7.7, 1.2, 110, 100, 60, 808.79, 428, 660, 333, 14, 2000; Antiseptic effect 2352, 880, 787, 760, 727, 465, 444, 450, 428, 660, 589, 555, 333, 14, 2000, 1570, 1770; Apoplexy 40, 20, 1800, 1875, 125, 95, 72, 20, 2309, 522, 146, 555, 333, 999, 880, 787, 727; Appetite, lack of 465, 2309, 125, 95, 72, 20, 10000, 880, 787, 727; Arteriosclerosis 10000, 2720, 2170, 1800, 1600, 1500, 880, 787, 776, 727, 20, 808.79; Arthritis, 10000, 1550, 880, 802, 787, 727, 28, 20, 7.7, 3.0, 1.2; Arthritis, rheumatoid of the muscles and tendons 250, 1.2, 1875, 787, 727, 808.79; Arthritis, arthrosis and parathyroid disturbances affecting calcium metabolism, 9.6; Arthritis, arthralgia due to gout, 9.4; Arthritis focal origin 9.4; Asthma, 2720, 2170, 1800, 1600, 1500, 880, 787, 727, 808.79; Astrocytoma 9.2, 8.25, 7.7, 2170, 2127, 880, 690, 666; Ataxia 2720, 2170, 1800, 1600, 1500, 880, 787, 776, 727, 650, 625, 600, 2309, 125, 95, 72, 20, 805.6, 811.91, spastic ataxia, 9.2, 8.25, 7.7; Athlete's Foot 465, 1550, 880, 802, 727, 20, Autoimmune Disorders 10000, 3352, 880, 787, 776, 727, 650, 625, 600, 465, 250, 28, 20, 9.6, 9.4, 7.7, 3.0, 522, 146, 1850; Auto-intoxication 0.5, 522, 146, 10000, 880, 787, 727, 20, 555, 333, 999; Back pain 26, 186, 3040, 522, 146, 555, 333, 999, 10000, 760, 2352, 880, 787, 727; Bad Breath 2352, 880, 787, 727, 20; Bedsores 880, 2352, 787, 727, 485 followed by 1.2 and 73; Biliousness 2352, 10000, 880, 832, 787, 727, 465; Bladder and prostate complaints 880, 2352, 787, 727, 465, 20, 9.4 Boils 880, 2352, 787, 727, 465, 660, 589, 333; Bone trauma 380, 2352, 10000, 880, 787, 727; Bone/periodontal disease 47.5, 1800, 1600, 1875, 880, 787, 776, 727; Bone protuberance 1.2 and 250; Brachial Neuralgia 0.5; Breast 880, 2352, 787, 776, 3083; Bright's Syndrome 1500, 880, 787, 727; Bronchial asthma 0.5, 522, 146, 125, 95, 72, 966, 20; Bronchitis 880, 727, 9.4, 9.35; Bronchial Pneumonia 2352, 880, 787, 776, 727; Bubonic Plague 500, plus secondary infections: 880, 787, 727, 20; Burns 190, 10000, 880, 787, 727, 465, 200; Bursitis 880, 787, 727; Cancer 465, 2127, 880, 787, 3083, 589; Cardiac Edema 9.2; Carpal Tunnel Syndrome 1550 thru 465, 6.3, 148, 2.5, 727, 522, 146, 444; Cataract 10000, 880, 787, 727; Catarrh 1550, 802, 880, 787, 727, 444, 20; Cerebral Palsy 880, 787, 727, 522, 146; Chicken Pox 2352, 1800, 1600, 1500, 880, 787, 727, 20; Cholera 880, 802, 450, 880, 787, 727; Chronic Fatigue Syndrome 10000, 2720, 1800, 880, 802, 787, 776, 727, 465, 966, 330, 148, 125, 95, 73, 72, 35, 27.5, 20, 6.3, 4.9, 428, 660; Circulatory stasis 40; Circulation disturbance problems 40.9.4; Concentration 7.83; Cold in head, chest 10,000, 1550, 880, 802, 787, 776, 727, 444, 20; Colic 1550, 832, 802, 787, 727, 20; Colitis 1550, 880, 832, 802, 440; Constipation 1550, 880, 802, 787, 776, 727, 20; Convulsions

10,000, 880, 787, 727; Contusions 9.1, 110; Costalgia 160, 26, 3040, 880, 787, 727, 10,000, 2353; Cramps 10,000, 880, 787, 727, 26; Cystitis (of urinary bladder) 1550, 880, 802, 787, 727, 465, 20; Deafness 10,000, 1550, 880, 802, 787, 727, 20; Dental foci 3040, 95, 190, 47.5, 1600, 1550, 2353, 1500, 880, 832, 787, 776, 760, 727, 1316, 1916, 465, 2720, 2489, 1800, 522, 146, 555, 333, 999; Depression, anxiety, trembling, weakness 3.5; Depression 35, 220, 1.1, 73; Detoxification 6.3, 148, 2.5, 727, 522, 146, 444, 666, 428, 880, 787, 555, 333, 999; Diabetes 837, 260, 1885, 125, 95, 1957, 522, 428, 49.4 for secondary diabetes infections 10,000, 2720, 2170, 1800, 1550, 880, 802, 727, 465, 20; Diabetic loading 35, 700; Diarrhea/Dysentery 1550, 880, 802, 787, 727, 465; Distorsion 9.1, 110; Dizziness 5.8; Dupuytren's Contracture 1.2, 250; Dysmenorrhea 26, 4.9, 2352, 880, 787, 727, 465; Dyspepsia 1550, 880, 802, 787, 727; Ears 9.2, 880, 787, 727, 20; Eczema 9.2, 1550, 802, 787, 727; Eczema 9.4; Enuresis 10,000, 880, 787, 727, 465; Edema 522, 146, 154.3, 444, 440, 880, 787, 727, 465; Epicondylitis 251.2, 3040, 160, 26; Epididymitis 1500, 880, 787, 727, 20; Epilepsy 10,000, 880, 802, 787, 727, 700, 650, 600, 125, 20; Epstein-Barr 428, 660, 776, 465, 880, 787, 727; Erysipelas 660, 2000, 10,000, 880, 787, 727, 465, 20; Erythema nosodum 9.4; Esophagus 880, 787, 727; Eustachian tube inflammation 1550, 880, 802, 787, 776, 727, 465; Eye inflammation 1.2, 80; Eye Disorders 1600, 10,000, 880, 787, 727, 20; Eyes 1600, 880, 787, 727; Facial paralysis 10,000, 880, 787, 727; 9.6, 3040, 880, 787, 727, 522, 72, 5.8, 4.9; Hearing problems 10,000, 2158, 880, 787, 760, 727, 20; Heart (lab animals only) 80, 160, 20, 73, 3.9, 10,000, 880, 787, 727, 465, 4125 95, 20; Heartburn 832, 2720, 2170, 2127, 1800, 1600, 2352, 1500, 880, 787, 727, 685, 465, 2309, 125, 95, 72, 20; Hemorrhage 1550, 802; Hemorrhoids 1550, 880, 802, 727; Hepatitis (generic) 29.2, 1550, 880, 802, 727; Hernia 10,000, 787, 727; Herpes (zoster) 2720, 2170, 1800, 1600, 1500, 1550, 802, 1865, 880, 787, 727, 20; High blood pressure 10,000, 880, 787, 727, 9.2; Hip Pain (as in coxarthrosis) 3040, 880, 787, 727, 160, 26; Hives 1800, 880, 787, 727, 522, 146, 4.9; Hoarseness 880, 760, 727; Hyperacidity of stomach 7.83 on solar plexus and 230 on stomach; Hydrocele 880, 787, 727; Hyperosmia 20, 10,000, 522, 1865, 810.35; Hyperthyroid 3.0, 0.5; Hypoacidity of stomach 15, 20; Hypothyroid 12, 35; Hypertension 20, 95, 9.2, 6.0, 10,000, 880, 787, 727; Hypertension, spastic 95; Impotence 9.4, 2127, 2008, 465, 880, 802, 787, 2083, 125, 95, 73, 72, 20, 1875; Indigestion 10,000, 880, 787, 727, 465, 444, 20, 125, 95, 72, 4.9; Infantile Paralysis 1500, 880, 787, 727; Infections 1600, 1550, 1500, 880, 832, 802, 787, 776, 760, 727, 700, 690, 666, 650, 625, 600, 465, 444, 125, 95, 72, 20, 1865, 500, 450, 440, 428, 660; Infertility 2127, 2008, 465, 880, 802, 787, 2083, 1875, 9.4; Inflammation 1.5; Influenza 1550, 802, 1500, 880, 787, 727, 20; Insomnia 3.6, 3.0, 1550, 1500, 880, 802; Intercoastal neuralgia 3040, 160, 26, 1550, 880, 802, 787, 776, 727, 125, 20, 1865, 444; Intermittent claudication 45, 48; Itching 880, 787, 727, 465, 444, 125, 95, 72, 7.83, 2309, 2000; Jaundice 1600, 1550, 1500, 880, 802, 650, 625, 600, 2309, 146, 250, 125, 95, 72, 20; Kidney insufficiency 50, 440, 1600, 1550, 1500, 880, 802, 787, 727, 650, 625, 600, 465, 73, 40, 10, 2309, 522, 146, 250, 148, 125, 95, 72, 20, 555, 333, 999; Knee/Joint pain 3040, 160, 26, 1550, 880, 802, 787, 727, 28, 20, 7.7, 3.0, 1.2, 250, 1.2, 9.6, 9.4; Larynx 440, 465, 444, 1550, 880, 802, 787, 727, 28, 7.7, 3.0, 1.2, 250, 1.2, 9.6, 9.4; Leprosy 600, 1600, 1550, 1500, 880, 832, 802, 787, 776, 760, 727, 700, 690, 666, 650, 625, 600, 465, 444, 20, 500, 450, 440, 428, 660; Leukemia 2127, 2008, 880, 787, 727, 690, 666, 589; Leukoderma 880, 787,

727, 444, 20, 589, 2000; Locomotor Dysfunction 10,000, 880, 787, 776, 727, 650, 625, 600, 2309, 125, 95, 72, 20; convulsions, spasticity 28, 8.25, 7.7; Low Blood Pressure 880, 787, 727, 20; Lumbago 10,000, 880, 787, 727, 125, 95, 72, 2309, 9.2, 8.25, 7.7; Lungs 1550, 802, 880, 787, 776, 727, 125, 95, 72, 20, 2309, 450, 589; Lupus Erythematosus 880, 787, 727, 776, 1850; Luxation 119.1; Lymph stasis 154.3, 2.5, 727, 522, 146, 444, 428, 880, 787; Malaria 555; Measles 333; Mental Disorders 444, 148, 20, 6.3, 2.5, 522, 146, 10,000, 125, 95, 72, 4.9, 727, 428, 444, 465 to 1550; Menieres Disease 2352, 880, 787, 727, 465, 428, 589, 333; Meningitis 2352, 880, 832, 787, 727, 650, 625, 600, 465, 2309, 125, 95, 72, 20, 428, 660; Menstrual problems 26, 3040, 880, 787, 727, 465, 160, 20; Migraine 170; Motion Sickness 650, 625, 600, 465, 444, 1865, 522, 146, 125, 95, 72, 20; Mouth eruptions 465, 2127, 2008, 727, 690, 666; Herpes sores 2489, 1800, 465, 1550, 1500, 880, 787, 727, 1850, 428; Mucous membrane inflammation 380; Multiple Sclerosis use auto immune frequencies plus 880, 787, 727; Muscles, to relax 760, 6.8; Muscular Dystrophy 10,000, 2352, 880, 787, 776, 727, 650, 625, 600, 465, 250, 28, 20, 9.6, 9.4, 7.7, 3.0, 522, 610, 1850, 146, 555, 333, 999; Muscular pain 3040, 26, 160, 320, 250, 240, 80, 40, 10, 5.8, 2.5, 1.5, 1.2, 1.0, 0.5; Muscle pain from parasites 2309, 125, 811, 108, 72, 20; Mumps 10,000, 727, 2720, 2489, 2127, 2008, 14, 428, 880, 787, 727, 20; Nausea 880, 832, 787, 727, 20, 4.9; Nephritis/nephrosis 880, 787, 727, 10, 40, 73, 465; Nervousness 3.0; Nerve disorders 10,000, 2720, 2489, 2170, 1800, 1600, 2352, 880, 787, 148, 6.3, 2.5, 522, 610, 727, 650, 625, 600, 125, 95, 72, 20, 428, 660, 333; Neuralgia 3.9, 555, 333, 999; Neurosis 28; Nose infection, congestion 2352, 880, 787, 776, 727, 444, 440, 465, 20; Numbness 10,000, 2720, 2489, 2170, 1800, 1600, 2352, 880, 787, 727, 650, 625, 600, 440, 660, 333; Oral Lesions 2720, 2489, 2008, 1800, 1600, 2352, 880, 787, 776, 760, 727, 465, 444, 522, 146, 555, 333, 999; Orchitis 2720, 2489, 2170, 2127, 2008, 1800, 1600, 2352, 1500, 880, 832, 787, 776, 2083, 650, 625, 600, 125, 95, 72, 20, 14; Osteomyelitis 2.6, 2720, 2489, 2170, 2127, 2008, 1800, 1600, 2352, 1500, 880, 832, 787, 776, 2083; Otosclerosis 9.2; Ovarian disorders 1875, 465, 444, 26, 2720, 2489, 2170, 2127, 2008, 1800, 1600, 2352, 1500, 880, 832, 787, 776, 2083, 20; Pain 3040, 95, 160, 26; Pain of infection 3040, 160, 26, 95, 880, 2352, 787, 776, 727, 4.9; Pain of cancer 3040, 160, 26, 95, 2127, 2008, 2083; Pancreatic insufficiency 15, 10, 250, 1875, 465, 444, 26, 2720, 2489, 2170, 2127, 2008, 1800, 1600, 2352, 1500, 880, 832, 787, 776, 2083, 20; Paralysis, spastic 10,000, 880, 787, 776, 727, 650, 625, 600, 2309, 125, 95, 72, 20, 7.7; Paralysis, nonspastic 10,000, 880, 787, 776, 727, 650, 625, 600, 444, 1865, 125, 95, 72, 20, 9.2, 8.25; Parasites 125, 95, 72, 20, 2309; Pelvic Inflammatory Disease 2720, 2489, 2170, 2127, 2008, 1800, 1600, 2352, 787, 776, 2083, 650, 625, 600, 465, 966, 95, 72, 450, 428, 14, 589, 1570; Pericarditis 2720, 2170, 1600, 880, 2352, 787, 760, 727, 625, 465, 160, 125, 73, 26, 20, 3.9, 95, 80, 73, 72, 2309; Pharyngitis 440, 589, 380, 1600, 2352, 776, 787, 880, 727, 20, 522, 146; Pleurisy 1550, 802, 880, 787, 776, 727, 125, 95, 72, 20, 2309, 450, 589; Pneumonia 1550, 802, 880, 787, 776, 727, 20, 450, 589; Poliomyelitis 805.6 secondary complications 2158, 428, 1500, 880, 787, 727; Polyps 2720, 2489, 2170, 2127, 2008, 1800, 1600, 2083, 1875, 465, 444, 20, 522, 146; Pre-op and post-operation prevention and control of nosocomial and idiopathic infection 2170, 1800, 1600, 2352, 1500, 880, 832, 802, 787, 776, 727, 465, 444, 522, 146, 428, 589, 333; Prostatitis 510, 522, 146, 2720, 2489, 2170, 2127, 2008, 2352, 787, 776, 2083, 465, 125, 95, 72, 20, 1249; Prostate complaints 9.4, 2127, 2008, 727, 690,

666, 465, 880, 802, 787, 727, 125, 95, 73, 72, 20; Proctitis 880, 787, 727, 465, 444, 125, 95, 72, 20, 2309, 2000; Prostate tumor (malignant) 2127, 2008, 2083; Psychosomatic pain 727, 160, 26, 148, 455, 511, 146, 555, 333, 999; Pyorrhea 2720, 2489, 2008, 1800, 1600, 2352, 880, 787, 776, 727, 465, 444, 522, 146; Rabies 807.23; Raynaud's Disease 727, 20; Rheumatoid arthritis of the muscles and tendons 251.2, 9.6, 9.4, 1875, 10,000, 2352, 880, 787, 727, 465, 20, 7.7, 3.0, 1.2; Rhinitis 10,000, 2352, 1500, 880, 787, 776, 727, 444, 7.83, 522, 146, 465; Sarcoma 2127, 2008, 880, 787, 727, 690, 666, 589; Scars 465 to 1568; Scarlet Fever 880, 787, 727, 1356; Sciatica or Ischia 2352, 880, 787, 727, 690, 666, 10; Sedative Effect 2.5; Sexual Dysfunction 9.4, 2127 2008, 465, 880, 802, 787, 2083, 125, 95, 73, 72, 20, 1875; Shingles 2352, 802, 1800, 1865, 1880, 787, 727, 26; Sinusitis 2352, 880, 787, 727, 125, 72, 522; Sleeping Sickness 124; Slipped discs 26, 125, 880, 787, 333, 787, 727, 95, 72, 522, 146, 555, 333, 999; Smallpox 811.91, 1550, 802, 880, 787, 727; Sneezing 880, 787, 727, 465, 146; Sore Throat 2720, 2489, 1800, 1600, 2352, 880, 787, 776, 727, 465, 589, 808.79; Spasms, muscle 6.8; Spastic paresis 650, 625, 30.87.48; Spleen, enlarged 35, 787, 10,000, 2720, 2170, 1800, 1550, 880, 802, 727, 465, 20; Spondylitis, acute 1.2, 10, 7.7, 28; Staphylococci infection 727; Stiff neck 4.9; Spastic stiff neck 9.2, 6.0; Stiff muscles in general 320, 250, 240, 160, 125, 80, 40, 20, 10, 5.8, 2.5, 1.5, 1.2, 1.0, 0.5, 1800, 2352, 880, 787, 776, 727; Stomach disorders 2127, 2008, 880, 787, 727, 2083, 125, 95, 72, 20, 3.9, 450, 1570, 1770; Stones 3.5; Streptococci infection 880; Streptothrix infection 787; Sty 10,000, 880, 787, 727, 20; Subluxation induced disorders 9.6; Sun Allergy 3.0, 330; Sunstroke 444, 428, 190, 3040, 95, 522, 146, 880, 20; Surgery—Prevention and control of nosocomial and idiopathic infection 2170, 1800 1600, 2352, 1500, 880, 832, 802, 787, 776, 727, 465, 444, 522, 146, 428, 589, 333, 3040, 95, 160, 26, 544;

Swelling 522, 146, 6.3, 148, 444, 428, 880, 787, 727; Swollen glands 880, 787, 727, 20; Syphilis 1875; Tachycardia (also used for pain of arthritis, headache, Tui-na, and facial toning) 1.2; Teeth 3040, 95, 190, 47.5, 1600, 2352, 1500, 880, 10,000, 2489, 787, 776, 760, 727, 1875, 465, 2720, 2489, 1800; Tendo myopathy 3040, 320, 250, 160, 80, 40, 26, 20, 10, 5.8, 2.5, 1.5, 1.2, 1.0, 0.5, 522, 146; Tetanus 811.91, 600, 880, 787, 727; Thrush 465; Tonsillitis 74.2, 73, 2352, 1500, 880, 832, 787, 776, 727, 1875, 465, 808.79, 589; Toothache 3040, 95, 190, 47.5, 2720, 2489, 1800, 1600, 2352, 1500, 880, 832, 787, 776, 727, 666, 650, 600, 465, 522, 146, 555, 333, 999; Tooth extraction, afterward 3040, 95, 47.5, 7.83; Toxins 444, 160, 522, 146, 555, 333, 999, 10,000, 880, 787, 727, 125, 95, 72, 20; Trauma 3040, 95, 190, 760, 880, 787, 727, 465; Trench Mouth 1875, 880, 787, 776, 727, 650, 625, 465; Trigeminal neuralgia 2720, 2489, 2170, 1800, 1600, 2352, 880, 832, 787, 776, 760, 727, 650, 146, 7.83, 27.5, 428, 660, 589, 333, 999, 770, 2000, 805.6, 807.23, 811.91, 810.35; Tuberculosis 2352, 1875, 776, 2127, 2008, 465, 589; Typhoid 1570, 1770, 2352, 690, 1800; Ulcers 2489, 2170, 2127, 1800, 1600, 880, 832, 802, 787, 776, 727, 1.2 through 5.8, 1.2, 73; Urethritis 2720, 2170, 2127, 1800, 1600, 2352, 1500, 880, 832, 787, 776, 727, 1875, 455, 2309, 125, 95, 72, 1.2; Urticaria 1800, 444, 125, 95, 72, 522, 146, 787, 727, 6.3; Varicoses 1.2, 28; Vegetative Dystonia 40; Vein thrombosis 685; Vertigo 60, 5.8; Warts 3040, 5, 2170, 2127, 2008, 1800, 1600, 1500, 2083; Worms 125, 95, 6.0, 20, 2309; Wound healing 880, 787, 727, 220, 190, 20, 40; Yeast infection 465; Yellow Fever 999; Scan brain for troubled areas 7.83; For relaxation 4.9; To stimulate mental clarity 12; Facial toning 9.6; Kidney meridian 9.2.

* * * * *